

# Corning® Cryogenic Storage Solutions

CORNING

A new and improved way to freeze your cells

official distributor

SZABO-SCANDIC HandelsgmbH  
Quellenstraße 110, A-1100 Wien  
T. +43(0)1 489 3961-0  
F. +43(0)1 489 3961-7  
mail@szabo-scandic.com  
www.szabo-scandic.com



**SZABO  
SCANDIC**

Part of Europa Biosite




# A New Standard in Cell Cryopreservation

At Corning, we continuously look for ways to help our customers improve or streamline steps in the cell culture workflow. One such area is cell cryopreservation. While current methods exist, they require chemicals and maintenance.

Now there is a new and improved way to freeze cells for cryogenic storage – we call it Corning® CoolCell®.

Corning CoolCell is an alcohol-free cell freezing container, which controls the rate of freezing to  $-1^{\circ}\text{C}/\text{minute}$  when placed in a  $-80^{\circ}\text{C}$  freezer. CoolCell has been performance tested with a variety of cell types including stem cells, primary cells, PBMC cell lines, insect cells, and yeast. The patent-pending CoolCell technology utilizes a thermo-conductive alloy core and highly insulative outer material to control the rate of heat removal and provide reproducible cell cryopreservation. CoolCell units are easy to use and deliver comparable results.

## Corning CoolCell Reproducibility

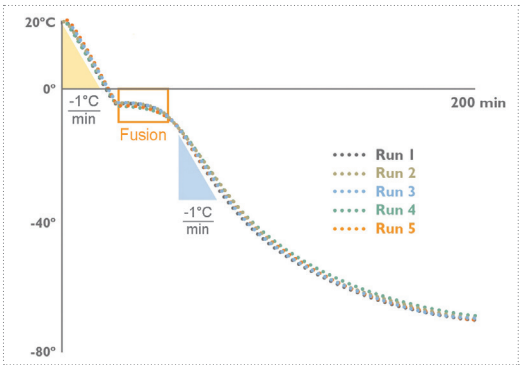


**Performance test:** A temperature probe was placed into a 2.0 mL cryogenic vial containing 1.0 mL of water and the tube was inserted into a room temperature Corning CoolCell. The CoolCell was placed into a  $-80^{\circ}\text{C}$  freezer and the temperature rate and profile was recorded over a 3-hour period. The test was repeated 5 consecutive times.

**Conclusion:** Corning CoolCell generated identical fusion time and cooling profiles over five consecutive freeze cycles.


### Alcohol-free with No Ongoing Costs or Maintenance

Isopropanol (IPA) containers used for cryogenic freezing require costly alcohol replacements every 5 uses, can be cumbersome to handle, and may have inconsistent freezing rates. Corning CoolCell is different, because it's a reusable, alcohol-free way to uniformly freeze your cells at a lower cost of use. With CoolCell, you can depend on high reproducibility and high cell viability, to ensure you preserve the most cells possible for your research.

## Corning CoolCell Performance Versus IPA Container

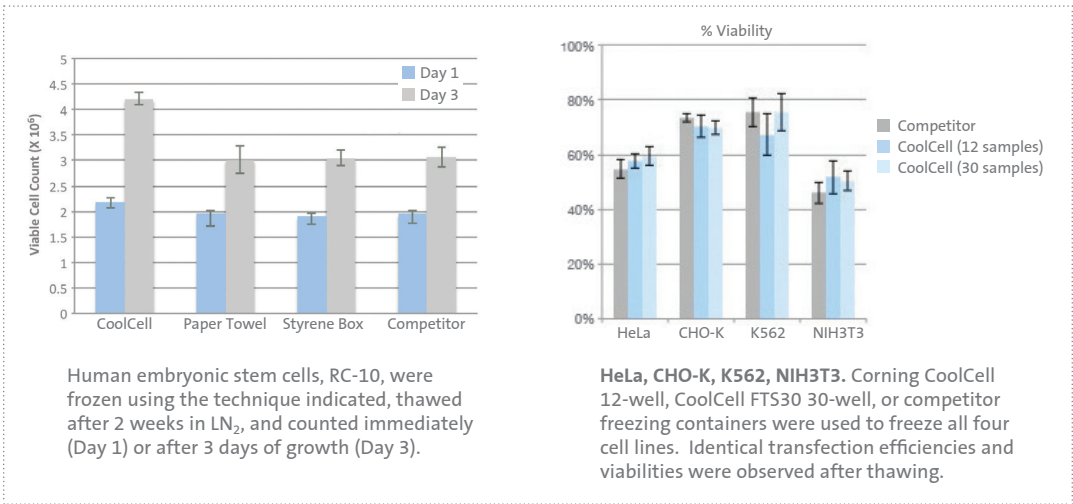

## How Corning CoolCell LX Works

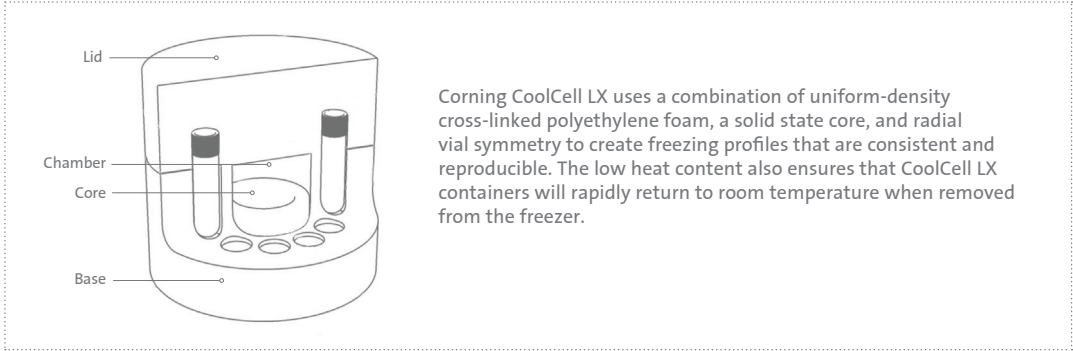



New



Corning 2D bar coded cryogenic tube, 2 mL, self-standing (Cat. Nos. 8670, 8671)



DMSO Media



Cryopreservation bags

## Features

Unique features of controlled-rate freezing with Corning® CoolCell® include:

- ▶ Ease of use
- ▶ Alcohol and fluid-free freezing
- ▶ Lower cost of use than alcohol-based devices
- ▶ High cell recovery and cell viability
- ▶ Reproducibility
- ▶ Simple, consistent way to standardize controlled-rate freezing

## Keep Your Samples Safe

Use Corning cryogenic vials and DMSO media with Corning CoolCell to further protect your valuable cell lines, biological, and aqueous solutions in ultra-low temperature storage. Choose from external or internal cryogenic vial thread caps or assorted color cap options to suit your needs. For added convenience, Corning also offers reusable cryogenic racks.

## New Improve Your Sample Management

Manage and manipulate multiple storage tubes more efficiently with Corning 1D/2D bar coded cryogenic vials. Corning synchronized permanent, laser-etched 2D bar codes on the bottom of the vial as well as linear bar codes on the side of the vial can be used for permanent sample identification and ensure a reliable chain of custody. Corning 1D/2D bar coded cryogenic vials are temperature-resistant polypropylene vials that can withstand temperatures down to -196°C and are compatible with most scanning and capper/decapper systems.

## Sort and Move With Ease

Corning cryogenic vial grippers feature a unique design to grasp internal- or external-threaded cryogenic vials. Easily sort or move vials while maintaining sterility and protecting fingers from frozen vials, dry ice, and liquid nitrogen.

## Also Available for Cell Freezing

### Corning Cryopreservation Bags

Corning also offers cryogenic storage containers designed for the storage, preservation, and transfer of cells. Features include: a unique bag film material that remains flexible at low temperatures and proprietary port designs that allow for increased flexibility. Learn more at [www.corning.com/lifesciences](http://www.corning.com/lifesciences).



## Ordering Information

### Corning® CoolCell® Containers

Cat. No.	Description	Capacity (Vials)	Exposed Vial Tops	Qty/Pk	Qty/Cs
432000	CoolCell, purple	12	No	1	1
432001	CoolCell LX, purple	12	Yes	1	1
432002	CoolCell LX, green	12	Yes	1	1
432003	CoolCell LX, orange	12	Yes	1	1
432004	CoolCell LX, pink	12	Yes	1	1
432005	CoolCell 5 mL LX, purple	12	Yes	1	1
432006	CoolCell FTS30, purple	30	Yes	1	1
432007	CoolCell FTS30, orange	30	Yes	1	1
432008	CoolCell FTS30, green	30	Yes	1	1
432009	CoolCell FTS30, pink	30	Yes	1	1
432010	CoolCell SV2	12	Yes	1	1
432011	CoolCell SV10	6	Yes	1	1
432012	CoolCell SV2 stem cell cryopreservation system Kit includes Corning CoolCell SV2, Corning XT Starter, and Corning CoolRack® SV2 module	—	—	1	1
432013	CoolCell SV10 stem cell cryopreservation system Kit includes Corning CoolCell SV10, Corning XT Starter, and Corning CoolRack SV10 module	—	—	1	1

### Corning CoolCell Accessories

432076	CoolCell filler vials, 2.0 mL	—	—	6	6
432077	CoolCell filler vials, 5.0 mL	—	—	6	6
432078	CoolCell FTS30 vial module	30	—	10	10
432136	Cryogenic vial grippers, multi-colored	—	—	5	5

### Corning Cryogenic Vials and Accessories

#### External Thread Cryogenic Vials

Cat. No.	Capacity (mL)	Style	Self-standing	Qty/Pk	Qty/Cs
New 8671	2.0	1D and 2D bar coded, round bottom	Yes	50	500
New 8676	2.0	1D bar coded, round bottom	Yes	50	500
430658	1.2	Conical bottom	Yes	50	500
430659	2.0	Round bottom	Yes	50	500
430662	4.0	Round bottom	Yes	50	500
430663	5.0	Round bottom	Yes	50	500

#### Internal Thread Orange Cap Cryogenic Vials

New 8670	2.0	1D and 2D bar coded, round bottom	Yes	50	500
New 8672	2.0	1D bar coded, round bottom	Yes	50	500
430487	1.2	Conical bottom	Yes	50	500
430488	2.0	Round bottom	Yes	50	500
430489	2.0	Round bottom	No	50	500
430490	4.0	Round bottom	No	50	500
430491	4.0	Round bottom	Yes	50	500
430656	5.0	Round bottom	Yes	50	500

**Warning!** Do not use cryogenic vials for storage in the liquid phase of liquid nitrogen. Only store vials in the vapor phase above the liquified gas. Always use appropriate safety equipment when removing vials from cryogenic storage.

#### Cryogenic Vial Racks

Cat. No.	Description	Qty/Pk	Qty/Cs
8673	Cryogenic storage box, Polycarbonate, holds 81 vials, designed to fit Corning 2D bar coded cryogenic vials	5	10
8674	Cryogenic storage box, Polycarbonate, holds 100 vials, designed to fit Corning 2D bar coded cryogenic vials	5	10

#### Cryopreservation Bags

Cat. No.	Size (mL)	Fill Volume	Qty/Cs
91-200-88	50	10 - 20 mL	1
91-200-89	250	30 - 70 mL	1
91-200-90	500	55 - 100 mL	1
91-200-91	750	80 - 190 mL	1

## official distributor

**SZABO-SCANDIC HandelsgmbH**  
 Quellenstraße 110, A-1100 Wien  
 T. +43(0)1 489 3961-0  
 F. +43(0)1 489 3961-7  
 mail@szabo-scandic.com  
 www.szabo-scandic.com



**SZABO SCANDIC**  
 Part of Europa Biosite

For more specific information on claims, visit the Certificates page at [www.corning.com/lifesciences](http://www.corning.com/lifesciences).

**Warranty/Disclaimer:** Unless otherwise specified, all products are for research use only. Not intended for use in diagnostic or therapeutic procedures. Corning Life Sciences makes no claims regarding the performance of these products for clinical or diagnostic applications.

For additional product or technical information, visit [www.corning.com/lifesciences](http://www.corning.com/lifesciences), or contact our Scientific Support Team at [ScientificSupport@corning.com](mailto:ScientificSupport@corning.com).

#### Corning Incorporated Life Sciences Europe

Corning BV  
 Fogstraat 12  
 1060 LJ Amsterdam  
 The Netherlands  
 Phone: +31 (0) 20 659 60 51  
 Fax: +31 (0) 20 659 76 73  
[CSEurope@corning.com](mailto:CSEurope@corning.com)  
[www.corning.com/lifesciences](http://www.corning.com/lifesciences)

#### Support Offices

##### EUROPE

**France**  
 t 0800 916 882  
 f 0800 918 636  
**Germany**  
 t 0800 101 1153  
 f 0800 101 2427

**The Netherlands**  
 t 31 20 655 79 28  
 f 31 20 659 76 73  
**United Kingdom**

**All Other European Countries**  
 t 31 (0) 20 659 60 51  
 f 31 (0) 20 659 76 73

**CORNING** | **FALCON** | **AXYGEN** | **GOSSELIN**

For a listing of trademarks, visit [www.corning.com/clstrademarks](http://www.corning.com/clstrademarks).  
 All other trademarks are the property of their respective owners.