



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.202.1401  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

**Data Sheet**  
***Mouse RIPK1 Kinase Assay Kit***  
**Size : 96-well**  
**Catalog #79897**

**Background:** Mouse RIPK1 is a member of the receptor-interacting protein kinase (RIP) family. It plays an important role in cell survival as well as cell death, including apoptosis and necroptosis. RIPK1 is highly expressed in human Alzheimer's disease (AD) brains and in AD mouse models. Importantly, it was shown that inhibition of RIPK1 results in decreased amyloid and inflammation levels as well as accelerated amyloid beta peptide (A $\beta$ ) degradation by microglia *in vitro*.

**Description:** The Mouse *RIPK1 Kinase Assay Kit* is designed to measure Mouse RIPK1 kinase activity for screening and profiling applications using ADP-Glo<sup>®</sup> Kinase Assay as a detection reagent. The Mouse *RIPK1 Kinase Assay Kit* comes in a convenient 96-well format, with enough purified recombinant Mouse RIPK1 enzyme, RIPK1 substrate, ATP and kinase assay buffer for 100 enzyme reactions.

**COMPONENTS:**

Catalog #	Reagent	Amount	Storage	
100340	Mouse RIPK1	10 $\mu$ g	-80°C	<b>Avoid multiple freeze/ thaw cycles!</b>
79334	5x Kinase assay buffer	1.5 ml	-20°C	
79686	ATP (500 $\mu$ M)	100 $\mu$ l	-20°C	
	RIPK1 substrate (MBP)	100 $\mu$ l	-20°C	
79696	96-well plate, white	1	Room Temp.	

**MATERIALS OR INSTRUMENTS REQUIRED BUT NOT SUPPLIED:**

ADP-Glo<sup>®</sup> Kinase Assay (Promega #V6930)  
Dithiothreitol (DTT, 1 M; optional)  
Microplate reader capable of reading luminescence  
Adjustable micropipettor and sterile tips  
30°C incubator

**APPLICATIONS:** Useful for studying enzyme kinetics and screening small molecular inhibitors for drug discovery and HTS applications.

**STABILITY:** Up to 6 months when stored as recommended.

**CONTRAINDICATION:** Avoid >1% DMSO

**REFERENCE:**

Ofengeim D., *et al. Proc. Natl. Acad. Sci. USA* **114(41)**: E8788-E8797 (2017)

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.202.1401  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

#### ASSAY PROTOCOL:

***All samples and controls should be tested in duplicate.***

- 1) Thaw **5x Kinase assay buffer**, **ATP** and **RIPK1 substrate (MBP)**. (Optional: If desired, add DTT to **5x Kinase assay buffer** to make a 10 mM concentration; e.g. add 10  $\mu$ l of 1 M DTT to 1 ml **5x Kinase assay buffer**. Prepare only enough 5x Kinase assay buffer with DTT as required for the assay, as any excess 5x kinase buffer/DTT cannot be stored and should be discarded.)
- 2) Prepare the master mixture (12.5  $\mu$ l per well): N wells x (3  $\mu$ l **5x Kinase assay buffer** + 0.5  $\mu$ l **ATP (500  $\mu$ M)** + 1  $\mu$ l **RIPK1 substrate (MBP)** + 8  $\mu$ l distilled water). Add 12.5  $\mu$ l to every well.

	Positive Control	Test Inhibitor	Blank
5x Kinase assay buffer	3 $\mu$ l	3 $\mu$ l	3 $\mu$ l
ATP (500 $\mu$ M)	0.5 $\mu$ l	0.5 $\mu$ l	0.5 $\mu$ l
RIPK1 substrate (MBP)	1 $\mu$ l	1 $\mu$ l	1 $\mu$ l
Water	8 $\mu$ l	8 $\mu$ l	8 $\mu$ l
Test Inhibitor	–	2.5 $\mu$ l	–
Inhibitor Buffer (e.g. 10% DMSO(aq))	2.5 $\mu$ l	–	2.5 $\mu$ l
1x Kinase buffer	–	–	10 $\mu$ l
Mouse RIPK1 (10 ng/ $\mu$ l)	10 $\mu$ l	10 $\mu$ l	–
Total	25 $\mu$ l	25 $\mu$ l	25 $\mu$ l

- 3) Prepare 10X concentrated inhibitor in an aqueous-based solution. *Note: Final DMSO concentration must be  $\leq$ 1%. Higher DMSO levels can significantly decrease the enzyme activity. For example, to test an inhibitor dissolved in 100% DMSO at 10  $\mu$ M, dilute 1 mM inhibitor with water to make a 100  $\mu$ M inhibitor in 10% DMSO(aq). Then, add 2.5  $\mu$ l of the 100  $\mu$ M solution to the assay to make a 1% DMSO concentration in the final reaction mixture.*
- 4) Add 2.5  $\mu$ l of Inhibitor solution of each well labeled as “Test Inhibitor.” For the “Positive Control” and “Blank,” add 2.5  $\mu$ l of the same solution without inhibitor (Inhibitor buffer, e.g. 10% DMSO(aq)).
- 5) Prepare 3 ml of **1x Kinase assay buffer** by mixing 600  $\mu$ l of **5x Kinase assay buffer** with 2400  $\mu$ l water. 3 ml of **1x Kinase assay buffer** is sufficient for 100 reactions.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.202.1401  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)


- 6) To the wells designated as "Blank," add 10  $\mu$ l of **1x Kinase assay buffer**.
- 7) Thaw **Mouse RIPK1 enzyme** on ice. Upon first thaw, briefly spin tube containing enzyme to recover full content of the tube. Calculate the amount of **Mouse RIPK1** required for the assay and dilute enzyme to 10 ng/ $\mu$ l with **1x Kinase assay buffer**. Store remaining undiluted enzyme in aliquots at -80°C. *Note: Mouse RIPK1 enzyme is sensitive to freeze/thaw cycles. Avoid multiple freeze/thaw cycles. Do not re-use thawed aliquots or diluted enzyme.*
- 8) Initiate reaction by adding 10  $\mu$ l of **diluted Mouse RIPK1 enzyme** to the wells designated "Positive Control" and "Test Inhibitor Control." Carefully shake the plate well and incubate it at 30°C for 50 minutes.
- 9) Thaw ADP-Glo reagent.
- 10) After the 50 minutes reaction, add 25  $\mu$ l of ADP-Glo reagent to each well. Cover plate with aluminum foil and incubate the plate at room temperature for 45 minutes.
- 11) Thaw Kinase Detection reagent.
- 12) After the 45 minutes incubation, add 50  $\mu$ l of Kinase Detection reagent to each well. Cover plate with aluminum foil and incubate the plate at room temperature for another 45 minutes.
- 13) Measure luminescence using the microplate reader. "Blank" value should be subtracted from all wells.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)

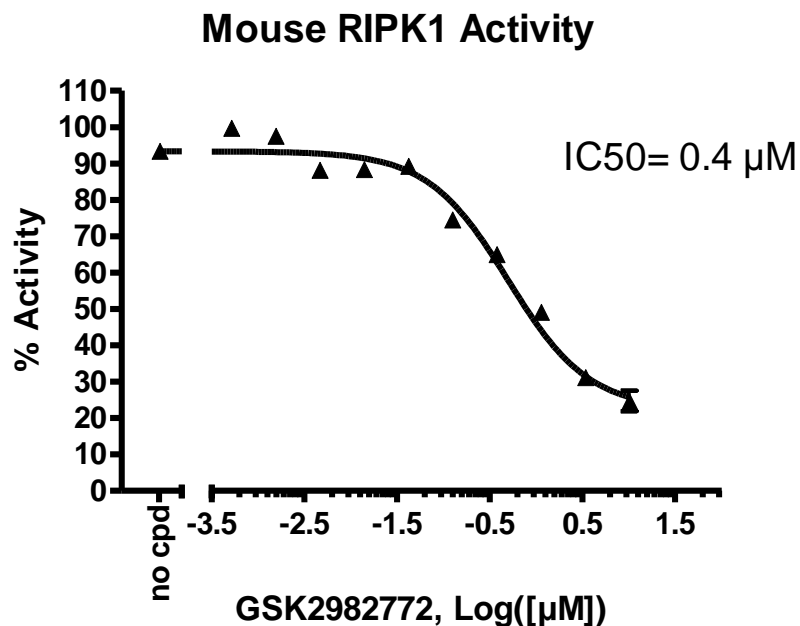

**Example of Assay Results:**

Inhibition of Mouse RIPK1 enzyme by GSK2982772, measured using the *Mouse RIPK1 kinase assay kit* (BPS Bioscience #79897). *Data shown is lot-specific. For lot-specific information, please contact BPS Bioscience, Inc. at [info@bpsbioscience.com](mailto:info@bpsbioscience.com)*

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.202.1401  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

**RELATED PRODUCTS:**

<b><u>Product Name</u></b>	<b><u>Catalog #</u></b>	<b><u>Size</u></b>
RIPK1, GST-Tag (Mouse)	100340	10 µg
RIPK1, GST-Th-Tag	40371	10 µg
RIPK2, GST-tag	40173	10 µg
RIPK5, GST-tag	40174	10 µg
RIPK1 Kinase Assay Kit	79560	96 rxns.
RIPK2 Kinase Assay Kit	79737	96 rxns.
LRRK2 (RIPK7), GST-Tag	40842	10 µg
LRRK2 (Y1699C), GST-Tag	40838	10 µg
LRRK2 (Y1699G), GST-Tag	40839	10 µg
LRRK2 (R1441H), GST-Tag	40837	10 µg
LRRK2 (R1441G), GST-Tag	40836	10 µg
LRRK2 (R1441C), GST-Tag	40835	10 µg
LRRK2 (G2385R), GST-Tag	40833	10 µg
LRRK2 (G2019S), GST-Tag	40832	10 µg
LRRK2 (D1994A), GST-Tag	40831	10 µg
LRRK2 (I2020T), GST-Tag	40834	10 µg

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [\*\*info@bpsbioscience.com\*\*](mailto:info@bpsbioscience.com)

Please visit our website at: [\*\*www.bpsbioscience.com\*\*](http://www.bpsbioscience.com)