



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



LipidLaunch™ LNP-102 Uptake Kit (Near-Infrared Fluorescence)

Item No. 39065

www.caymanchem.com

Customer Service 800.364.9897

Technical Support 888.526.5351

1180 E. Ellsworth Rd · Ann Arbor, MI · USA



TABLE OF CONTENTS

| | | |
|----------------------|----|-----------------------------------|
| GENERAL INFORMATION | 3 | Materials Supplied |
| | 3 | Safety Data |
| | 4 | Precautions |
| | 4 | If You Have Problems |
| | 5 | Storage and Stability |
| | 5 | Materials Needed But Not Supplied |
| INTRODUCTION | 6 | Background |
| PROTOCOL PREPARATION | 8 | Protocol |
| | 9 | Reagent Preparation |
| PROTOCOL | 10 | Performing the Protocol |
| RESOURCES | 14 | References |
| | 15 | Notes |
| | 15 | Warranty and Limitation of Remedy |

GENERAL INFORMATION

Materials Supplied

| Item Number | Item Name | Quantity/Size | Storage Temperature |
|-------------|---|---------------|---------------------|
| 33474 | SM-102 | 1 vial/25 mg | -20°C |
| 15100 | 1,2-Distearoyl- <i>sn</i> -glycero-3-PC | 1 vial/10 mg | -20°C |
| 9003100 | Cholesterol | 1 vial/25 mg | -20°C |
| 33945 | DMG-PEG(2000) | 1 vial/5 mg | -20°C |
| 34954 | 1,1-Dioctadecyl-3,3,3,3-tetramethylindotricarbocyanine (iodide) | 1 vial/5 mg | -20°C |

If any of the items listed above are damaged or missing, please contact our Customer Service department at (800) 364-9897 or (734) 971-3335. We cannot accept any returns without prior authorization.



WARNING: THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

Safety Data

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent *via* email to your institution.

Precautions

Please read these instructions carefully before beginning this assay.

If You Have Problems

Technical Service Contact Information

Phone: 888-526-5351 (USA and Canada only) or 734-975-3888

Email: techserv@caymanchem.com

In order for our staff to assist you quickly and efficiently, please be ready to supply the lot number of the kit (found on the outside of the box).

Storage and Stability

This kit will perform as specified if stored as directed in the Materials Supplied section (see page 3) and used before the expiration date indicated on the outside of the box.

Materials Needed But Not Supplied

1. A device capable of measuring fluorescence with excitation and emission wavelengths of 750 and 780 nm, respectively
2. Absolute ethanol
3. 50 mM sodium acetate, pH 5.0
4. Nucleic acid payload
5. Commercial microfluidic device, off-the-shelf microfluidic device, extrusion kit, or sonicator
6. PBS, pH 7.4

Background

Lipid nanoparticles (LNPs) are a subset of lipid-based drug delivery (LBDD) systems that utilize ionizable cationic lipids, such as SM-102, for the delivery of nucleic acid (e.g. siRNA, mRNA, cyclic dinucleotides) payloads to cells.¹ They consist of a lipid shell composed of structural phospholipids, cholesterol, and PEGylated lipids that surround an internal aqueous core, where the ionizable cationic lipids organize into inverted micelles around the encapsulated nucleic acids.² Ionizable cationic lipids are near-neutral at physiological pH and cationic in acidic environments (theoretical pK_a for SM-102 = 6.75), a property that promotes the encapsulation of negatively charged nucleic acids during LNP preparation and facilitates intracellular delivery after cellular uptake.^{1,3} LNPs are internalized into cells *via* endocytosis.⁴ SM-102 becomes protonated and positively charged in the acidic environment of the endosomal compartment, promoting LNP endosomal escape and intracellular delivery.⁵

Cayman's LipidLaunch™ LNP-102 Uptake Kit (Near-Infrared Fluorescence) is intended to serve as a starting point for laboratories to explore the uptake of LNPs *in vivo* without the need for specialized equipment. Optimal preparation conditions for the encapsulation of nucleic acids with LNPs must be determined by the end user. Adjustment of the following parameters may facilitate this process:

- Lipid molar ratio
- Lipid:nucleic acid (w:w) ratio
- Ionizable cationic lipid nitrogen:nucleotide phosphate (N:P) molar ratio
- Aqueous buffer: identity and ionic strength
- Particle size: extrusion size or microfluidic operating parameters, as applicable
- LNP preparation method

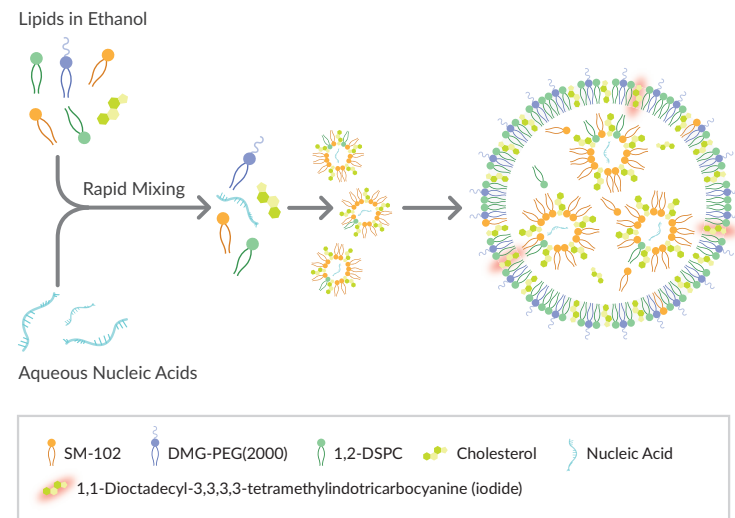


Figure 1. Schematic of nucleic acid-containing LNP formation with the LipidLaunch™ LNP-102 Uptake Kit (Near-Infrared Fluorescence)

PROTOCOL PREPARATION

Protocol

An example for preparing nucleic acid-containing LNPs with an ethanolic lipid mixture containing SM-102 (Item No. 33474), 1,2-distearoyl-*sn*-glycero-3-PC (1,2-DSPC; Item No. 15100), cholesterol (Item No. 9003100), 1,1-Dioctadecyl-3,3,3,3-tetramethylindotricarbocyanine (iodide) (DiR; Item No. 34954), and DMG-PEG(2000) (Item No. 33945) at lipid molar ratios of 50:10:37.5:1:1.5, respectively, is shown below. These molar ratios are derived from those used in mRNA-based vaccines. This example is shown with a lipid:nucleic acid (w:w) ratio of 10:1 and an ethanol:aqueous ratio of 1:3. The end user may scale volumes and adjust lipid molar and lipid:nucleic acid ratios as desired. It is possible to produce multiple small batches of LNPs using the parameters in this example and the reagents provided in the kit.

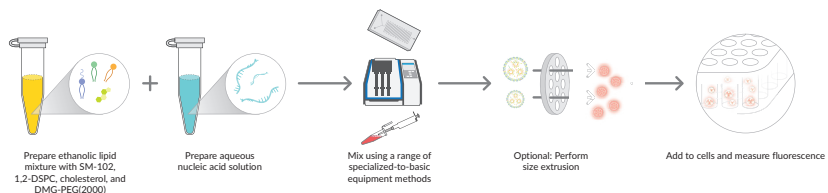


Figure 2. LipidLaunch™ LNP-102 Uptake Kit (Near-Infrared Fluorescence) workflow

Reagent Preparation

1. Ethanolic Lipid Mixture

Prepare individual lipid stock solutions of the four lipids supplied as crystalline solids in absolute ethanol. SM-102 is ready to use as supplied. Bring all stock solutions to room temperature prior to use and ensure they are well-dissolved. Transfer the appropriate volume of each lipid mixture component to a single tube as listed in the table below to prepare the ethanolic lipid mixture. Mix by pipetting several times.

| Lipid Mixture Component | Stock Solutions | | Working Mixture | | |
|-------------------------|-----------------|---------|-----------------|------|-----------------|
| | mg/ml | MW | Molar Ratio | mg | Required Volume |
| SM-102 | 100 | 710.2 | 50.0 | 3.40 | 34 μ l |
| 1,2-DSPC | 25 | 790.2 | 10.0 | 0.76 | 30 μ l |
| Cholesterol | 5 | 386.7 | 37.5 | 1.39 | 278 μ l |
| DMG-PEG(2000) | 1 | 2,509.2 | 1.5 | 0.36 | 360 μ l |
| DiR | 1 | 1,013.4 | 1.0 | 0.10 | 97 μ l |
| Absolute ethanol | | | | | 201 μ l |
| Total | | | | 6.00 | 1 ml |

Table 1. Preparation of ethanolic lipid mixture

2. Aqueous Nucleic Acid Solution

In a separate tube, bring nucleic acid payload to a concentration of 50-200 μ g/ml in 3 ml of 50 mM sodium acetate, pH 5.0. Optimal nucleic acid:lipid ratio should be determined empirically for each payload.

Performing the Protocol

Several methods are suitable for laboratory-scale, small-volume LNP production. These are described briefly below, and they may be adapted for use with a range of specialized-to-basic equipment. The procedures are performed at room temperature unless otherwise indicated.

1. Mixing

Commercial Microfluidic Device Mixing: Mix the ethanolic lipid mixture with the aqueous nucleic acid solution using a microfluidic device or chip with a staggered herringbone-, T-, or Y-channel design. The inlet flow rates can be controlled with syringe pumps and should be optimized by the end user.

Off-the-Shelf Microfluidic Device Mixing: Mix the ethanolic lipid mixture with the aqueous nucleic acid solution using an off-the-shelf microfluidic mixing device. These mixers can be assembled with common and inexpensive materials. Two inlets, composed of individual syringes containing the ethanolic lipid mixture and the aqueous nucleic acid solution, can be connected to opposite ends of a T- or Y-connector (2 mm I.D.) with appropriate tubing (1.5 mm I.D.) and fittings. A single outlet fitted with the appropriate tubing will direct the LNPs into a collection tube. The inlet flow rates can be controlled with syringe pumps and should be optimized by the end user.

Solvent-injection Mixing: Rapidly inject the ethanolic lipid mixture into the aqueous nucleic acid solution using a syringe with the needle placed in the center of the solution. Stir at 400 rpm for 30 minutes above the 1,2-DSPC transition temperature ($T_m = 55^\circ\text{C}$).

Hand Mixing: Hand mix the ethanolic lipid mixture with the aqueous nucleic acid solution *via* pipette by rapidly transferring the ethanolic lipid mixture into the aqueous nucleic acid solution. Mix by repeated pipetting for 15 seconds. Leave undisturbed for 10 minutes.

2. Final Preparation

- Perform size extrusion if necessary. Size extrusion is typically required following mixing techniques that produce large and heterogenous LNPs (e.g. off-the-shelf microfluidic device mixing, solvent-injection, and hand mixing) to yield a narrower LNP size distribution.
- Dialyze LNPs in neutral buffer (e.g. PBS, pH 7.4) against 1,000 volumes of buffer using the appropriate molecular weight cut-off (MWCO) membrane overnight.
- If desired, LNP solutions may be concentrated by centrifugation using the appropriate MWCO filter.
- Filter-sterilize LNP solutions with a 0.22 μm filter and store at 4°C until use. The LNP solutions will be stable at 4°C for one week. If longer storage is required, stabilizing reagents such as sucrose may be added for storage at -20°C . Stability testing should be conducted for each LNP formulation method to optimize storage conditions.

3. Characterization and Validation

A variety of techniques are available to characterize LNPs prior to *in vitro* or *in vivo* use. Contact Cayman Services for *in vitro* testing of your LNPs.

| Attribute | Assay(s) |
|------------------------------------|---|
| Particle size and distribution | Dynamic light scattering (DLS) |
| Zeta potential | Laser doppler electrophoresis |
| Lipid quantification and integrity | RP-HPLC, SE-HPLC, IP-HPLC |
| LNP morphology | Microscopy (cryo TEM, ESEM, AFM) |
| LNP uptake | Cell-based imaging or fluorescence detection |
| Translation or knockdown analyses | Cell-based reporter assays, Western blotting |
| Encapsulation efficiency | RNA-compatible nucleic acid dye with and without detergent, such as 0.5% Triton X-100 |

Table 2. LNP attributes and corresponding assays Adapted from Schoenmaker, L., *et al.*⁶

4. Uptake analysis

DiR incorporated into the LNPs will fluoresce at excitation/emission wavelengths of 750/780 nm, respectively. For *in vitro* studies, dilute prepared LNPs in medium containing serum about 1:100 to 1:500, then add to cells. After 24-48 hours, exchange medium for PBS and image or capture total fluorescence using a plate reader. Formulated nanoparticles can also be used *in vivo* to determine organ targeting with near-infrared-compatible *in vivo* imaging systems.

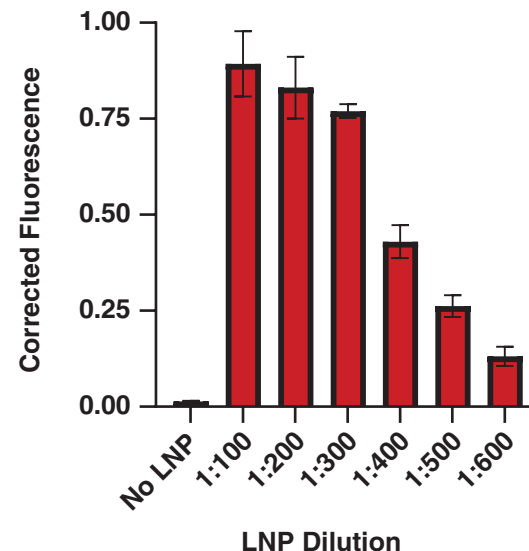


Figure 3. LNP-102 uptake in hepatocytes. LNP-102 was formulated in a Nunchuck device from Unchained Labs using a flow rate ratio of 3 and a total flow rate of 15 ml/min. After dialysis against PBS, LNPs were diluted in cell culture media and added to Huh7 hepatocytes for 24 hours. Hoescht 33342 (Item No. 15547) and DiR fluorescence were measured using a Tecan Infinite M Plex microplate reader. Hoescht 33342 was used to calculate the corrected fluorescence of LNP-102 uptake.

References

1. Mitchell, M.J., Billingsley, M.M., Haley, R.M., *et al.* Engineering precision nanoparticles for drug delivery. *Nat. Rev. Drug Discov.* **20**(2), 101-124 (2021).
2. Guevara, M.L., Persano, F., and Persano, S. Advances in lipid nanoparticles for mRNA-based cancer immunotherapy. *Front. Chem.* **8**, 589959 (2020).
3. Buschmann, M.D., Carrasco, M.J., Alishetty, S., *et al.* Nanomaterial delivery systems for mRNA vaccines. *Vaccines (Basel)* **9**(1), 65 (2021).
4. Degors, I.M.S., Wang, C., Rehman, Z.U., *et al.* Carriers break barriers in drug delivery: Endocytosis and endosomal escape of gene delivery vectors. *Acc. Chem. Res.* **52**(7), 1750-1760 (2019).
5. Carrasco, M.J., Alishetty, S., Alameh, M.-G., *et al.* Ionization and structural properties of mRNA lipid nanoparticles influence expression in intramuscular and intravascular administration. *Commun. Biol.* **4**(1), 956 (2021).
6. Schoenmaker, L., Witzigmann, D., Kulkarni, J.A., *et al.* mRNA-lipid nanoparticle COVID-19 vaccines: Structure and stability. *Int. J. Pharm.* **601**, 120586 (2021).

Warranty and Limitation of Remedy

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

This document is copyrighted. All rights are reserved. This document may not, in whole or part, be copied, photocopied, reproduced, translated, or reduced to any electronic medium or machine-readable form without prior consent, in writing, from Cayman Chemical Company.

©02/15/2024, Cayman Chemical Company, Ann Arbor, MI, All rights reserved.
Printed in U.S.A.

