



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## EVs SEC

### Size Exclusion Chromatography Columns for EVs isolation from biological fluids

#### INTRODUCTION

Extracellular vesicles (EVs) are thought to provide a means of intercellular communication (2,3) and of transmission of macromolecules between cells allowing the spread of proteins, lipids, mRNA, miRNA and DNA and as contributing factors in the development of several diseases. EVs can also modulate cancer microenvironment (4) and the immune response (5,6) and are recognized as potential markers of human diseases, including cancer. In particular, due to their presence and stability in most body fluids and the similarity of their content with tumor cells, EVs have great potential as non-invasive biomarkers for liquid biopsy. Nevertheless, the use of EVs for diagnostic purposes has been limited by the lack of reproducible isolation methods (7). About isolation methods, ultracentrifugation remains by far like the most widely used primary method despite is a time consuming method, that provide aggregated proteins and nucleic acid may be pelleted and is not very suitable for small sample volumes. In this sense, size exclusion chromatography (SEC) has been described as most efficient method for isolating EVs from complex biological fluids by single-step (8), with a good recovery and with almost complete removal of contaminants, such as proteins and lipoproteins.

#### PRODUCT DESCRIPTION

Immunostep has developed SEC columns for EVs isolation from complex biological fluids such as: plasma, serum, urine and cerebrospinal fluid.

#### SPECIFICATIONS

Reference	SEC7012	SEC3512
Separation size	70-1000nm	35-350nm
Volumetric Flow Rate at RT	>0.75 ml/min	>0,75 ml/min
Sample volume	<1 ml (500 µl optimal volume)	
Column volume	10 ml	
Void volume	3 ml	
pH stability working range	3-13	
pH stability cleaning in place (CIP)	2-14	
Shelf life	12 months at 2-8°C	

#### REAGENTS PROVIDED

DESCRIPTION	REFERENCE	COLUMNS (UNITS)
12 ml SEC columns, which a diameter of 1.56 cm and 7.9 cm high, packed with a resin for a separation size of <b>35-350 nm</b> .	SEC3512-4	4
	SEC3512-8	8
12 ml SEC columns, which a diameter of 1.56 cm and 7.9 cm high, packed with a resin for a separation size of <b>70-1000 nm</b> .	SEC7012-4	4
	SEC7012-8	8

#### APPROPRIATE STORAGE AND HANDLING CONDITIONS

Store the columns at 2-8 °C in the presence of a bacteriostatic agent (e.g. filtered PBS containing Sodium Azide 0.09% or 20% ethanol). This one can be introduced during the rinsing step. Keep columns upright.

#### REAGENTS NOT PROVIDED

Running buffer is not supplied.

Some examples of buffers that can be used are:

- PBS 20 mM, 0.32% trisodium citrate (pH 7.4, 0.22 mm filtered and degassed).
- PBS 20 mM, 150 mM NaCl ( pH 7.2, 0.22 mm filtered and degassed)

It is recommended to use buffers with an ionic strength of 0.15 M or higher, because in low ionic strength buffers, the concentration of competing buffer ions is low, and contaminant proteins spend most of their time adsorbed to binding sites on the stationary phase.

For buffer choice it is also important consider the downstream protocol for isolated vesicles. For example citrate buffers could interfere with molecular biology protocols.

#### EVIDENCE OF DETERIORATION

Columns are reusable up to 5 times and should not be used if any evidence of deterioration is observed. For more information, please contact our technical service: [tech@immunostep.com](mailto:tech@immunostep.com)



#### RECOMMENDATIONS AND WARNINGS

- The column anti-bacterial solution contains sodium azide. In acid conditions, it is transformed into hydrazoic acid, a highly toxic compound. Azide compounds must be diluted in running water before being discarded. These conditions are recommended so as to avoid deposits in plumbing, where explosive conditions could develop. The safety data sheet (SDS) is available online at [www.immunostep.com](http://www.immunostep.com)*
- The biological samples should be handled in the same way as those capable of transmitting infection. Appropriate handling procedures should be guaranteed.*
- Do not use after the expiry date indicated on the column.*
- Temper running buffer and columns before sample purification for an adequate Volumetric Flow Rate.*
- For professional use only.*


#### WARRANTY

Warranted only to conform to the quantity and contents stated on the label or in the product labelling at the time of delivery to the customer. Immunostep disclaims hereby other warranties. Immunostep's sole liability is limited to either the replacement of the products or refund of the purchase price.

## PROTOCOL

### Preparation for use

1. Place the column in a rack and level it, make sure the column is vertical. Column box can be used as a

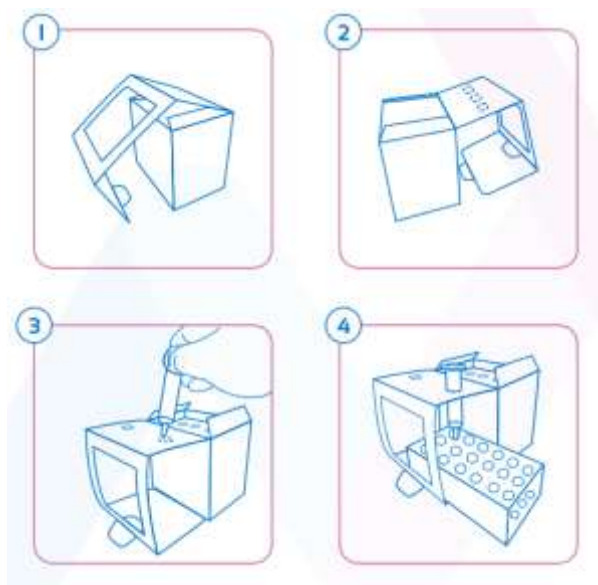

rack, just following next scheme (Fig.1).

Figure 1. Graphic instructions to use column box as well as rack, avoiding to get a specific column rack in the laboratory and facilitating the protocol.

2. Remove the bottom cap.
3. Leave the top-cap carefully.

### Column equilibration

4. Rinse the column with 21 mL of running buffer.
5. Make sure enough time is given for the column to be in the operational temperature range.
6. Do not allow the column to run dry. The top filter must stay wet.
7. Use only fresh filtered (0.2 µL) buffer to avoid particulate contamination.

### Sample fraction collection

8. When the column is used according to the protocol, for 500µl fractions, the first four fractions (2.0 mL), is the void volume which does not contain vesicles. These vesicles elute predominantly in fractions 5, 6 and 7 with removal of protein contamination. Fractions beyond 8 usually contain higher protein and lower vesicle levels (Fig 2).

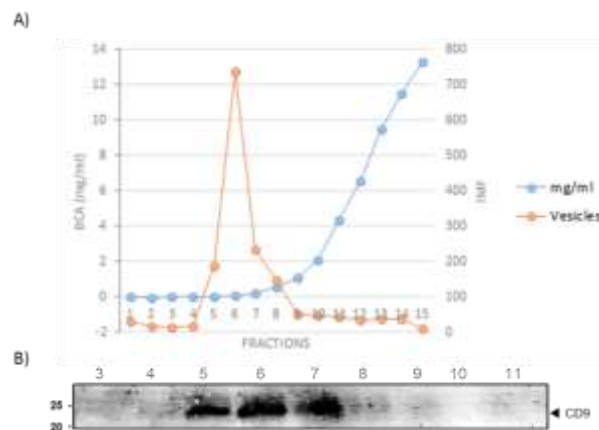


Figure 2. Amount of EVs and protein in each fraction from the column. A) Comparative of protein (BCA) vs vesicles (FACS CD63+/CD9+) content. B) Western Blot. SEC fractions were loaded on SDS-PAGE and immunoblotted for CD9 tetraspanin with anti-CD9 (VJ1/20), under non-reducing conditions

### Post collection of vesicles

9. After the collection of the vesicle fractions, rinse the column with 21 mL of running buffer to rinse out all the protein and small molecules before the next sample application. After that, rinse the column with 10 mL of a bacteriostatic agent and store as indicated in the storage section.

### Column cleaning

Cleaning is performed in order to remove any precipitated proteins or other contaminants that can build up on the column. Cleaning is performed using 0.5 ml of sodium hydroxide (500 mM NaOH). Note that the column should never be stored in sodium hydroxide. Equilibrate the columns immediately after the cleaning with 10 ml of water followed by 21 ml of running buffer.

### Fraction collection concentration

If the volume of the starting sample is small and / or the concentration and EVs is low in the sample, it is possible that the EVs will be much diluted in the collected fractions. In many cases, for downstream processing it is necessary to concentrate the fractions, for that, we recommended to use an Amicon Ultra-4, membrane PLHK Ultracel-PL, 100 kDa (Merck).

identifying in which fractions there is only protein and therefore there are no EVs. C) Flow cytometric analysis of elution fractions. Stain index: (MFI positive- MFI background)/ 2 $\sigma$  background)

## PERFORMANCE DATA

Below is a representative example of the isolation of EVs from a plasma by Immunostep SEC columns. Sample obtained can be subsequently analyzed by Abs (280) (Fig. 3A) and Abs (260/280) (Fig. 3B), which allows monitoring the elution of the sample, by conventional cytometry with the use of kits based on beads (Fig. 3C), such as Exostep (Ref. ExoS-25-P81), NTA (Fig. 4) or by PCR or RT-PCR (Fig. 5).

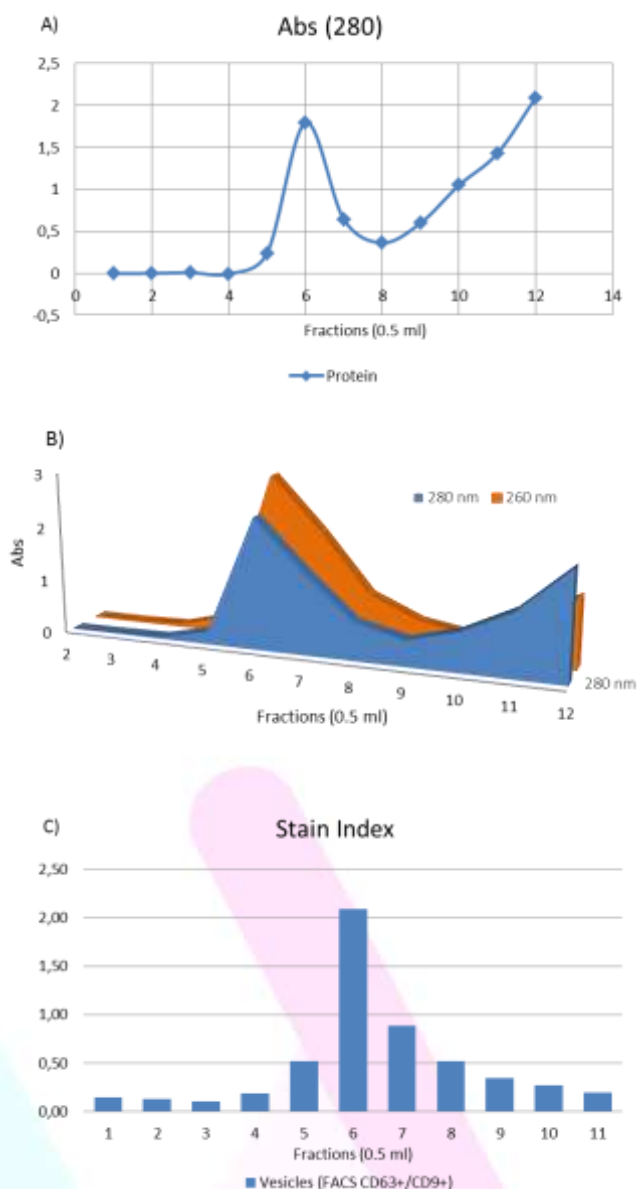


Figure 3. EVs isolation from plasma sample. A) Elution profile monitored by Abs (280). B) 260/280 ratio makes it possible to determine more clearly the fractions in which there is a presence of EVs, also


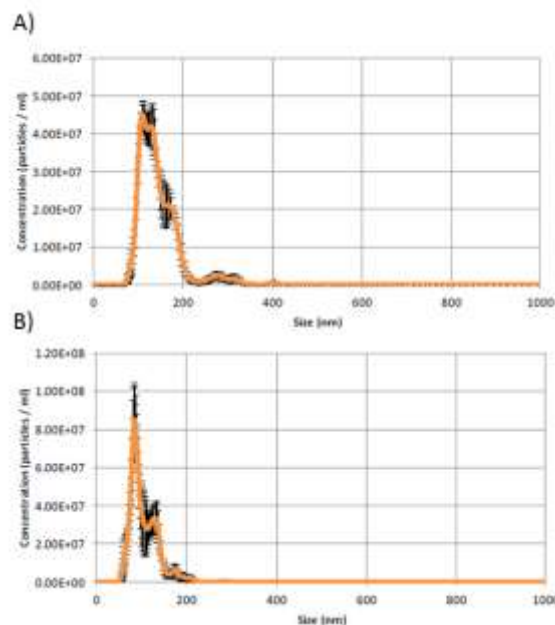


Figure 4. NanoSight analysis of the EVs recovered from Plasma by SEC columns. A) Column reference SEC1270 (70-1000 nm) B) Column reference SEC1235 (35-350 nm). EVs isolated by SEC resulting in a high level of purity.

A comparison was made using columns competitor, obtaining Immunostep columns a better performance in terms of recovery. A) Comparison of the diameter of the isolated particles by column type. The D50 value means 50% of the particles have that diameter or less. B) Concentration (particles / ml) obtained from purifying plasma (500  $\mu$ l).

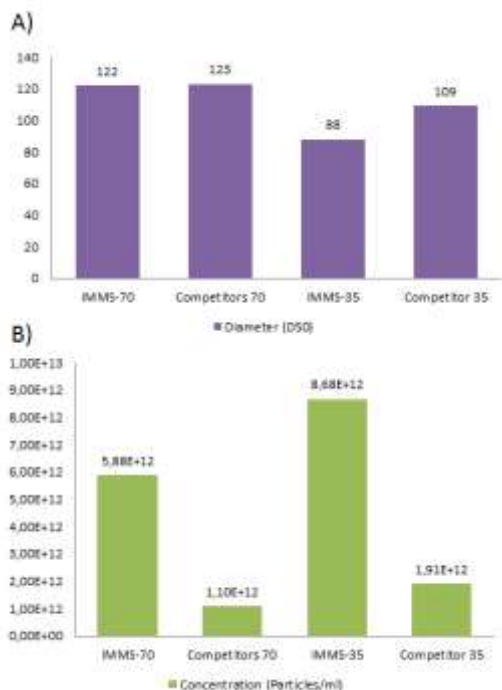

- Becker A, Thakur BK, Weiss JM1, Kim HS, Peinado H, Lyden D. Extracellular Vesicles in Cancer: Cell-to-Cell Mediators of Metastasis. *Cancer Cell* 2016 Dec 12;30(6):836-848.
- López-Cobo S, Campos-Silva C, Valés-Gómez M. Glycosyl-Phosphatidyl-Inositol (GPI)-Anchors and Metalloproteases: Their Roles in the Regulation of Exosome Composition and NKGD2-Mediated Immune Recognition. *Front Cell Dev Biol*. 2016 Sep 12;4:97.
- Jonathan M. Pitt, Guido Kroemer, Laurence Zitvogel Extracellular vesicles: masters of intercellular communication and potential clinical interventions. 2016 *J Clin Invest*. 2016;126(4):1139-1143
- Jara-Acevedo R, Campos-Silva C, Valés-Gómez M, Yáñez-Mó M, Suárez H, Fuentes M. Exosome beads array for multiplexed phenotyping in cancer. *J Proteomics*. 2019; Apr 30;198:87-97.
- Böing, A., van der Pol, E., Grootemaat, A., Coumans, F., Sturk, A. and Nieuwland, R. (2014). Single-step isolation of extracellular vesicles by size-exclusion chromatography. *Journal of Extracellular Vesicles*, 3(1), p.23430.

## MANUFACTURED BY



**Immunostep S.L.**  
Avda. Universidad de Coimbra, s/n  
Cancer Research Center (CIC)  
Campus Miguel de Unamuno  
37007 Salamanca (Spain)  
Tel. (+34) 923 294 827  
www.immunostep.com

Despite the columns are reusable up to 5 times, a single use is advisable when the downstream analysis is by PCR or RT-PCR.

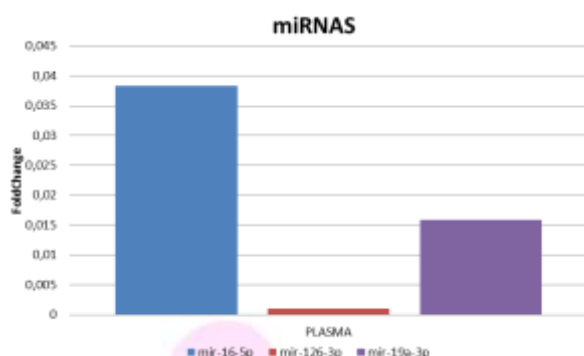


Figure 5. Analysis of the exosomal miRNA levels by quantitative RT-PCR. We followed the instructions in the QIAGEN Handbook to extract miRNAs, do the RT-PCR to obtain the cDNA and then the qPCR or quantitative PCR. Column protocol does not interfere with the extraction of RNAs, nor with RT-PCR and allows the detection of endogenous miRNAs miR-126-3p, miR-19a-3p and miR16-5p present in the sample.

## REFERENCES

- Yáñez-Mó M, Siljander P, Andreu Z, Bedina Zavec A, Borràs F, Buzas E et al. Biological properties of extracellular vesicles and their physiological functions. *Journal of Extracellular Vesicles*. 2015;4(1):27066.
- Pitt JM, André F, Amigorena S, Soria JC, Eggermont A, Kroemer G, Zitvogel L. Dendritic cell-derived exosomes for cancer therapy. *J Clin Invest*. 2016.
- Tkach M, Théry C. Communication by Extracellular Vesicles: Where We Are and Where We Need to Go. 2016 *Cell* 10;164(6):1226-32.