SZABO SCANDIC

Part of Europa Biosite

Produktinformation

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mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Anti-native DNA (nDNA) Antibody Test

Crithidia luciliae Substrate

IVD For in vitro Diagnostic Use

PRODUCT INSERT

REF 1106 48 Determination

REF 1106-2 96 Determination

REF 1106-6 120 Determination

INTENDED USE

An indirect immunofluorescence antibody test for the detection and quantitation of antibodies to native (double stranded) deoxyribonucleic acid (nDNA) in human serum.

SUMMARY AND EXPLANATION

Antibodies to nDNA are specific for systemic lupus erythematosus (SLE) and rarely occur in patients with rheumatoid arthritis, scleroderma or other autoimmune disorders¹. The frequency and titer of these antibodies fluctuate with disease activity and tend to disappear upon immunosuppressive treatment and during remission. There is a good correlation between the disease activity and anti-nDNA antibody levels²⁻⁸.

The two most commonly employed methods for detecting anti-nDNA antibodies are radioimmunoassay and immunofluorescence. The specificity and sensitivity of the *Crithidia luciliae* immunofluorescent method are comparable or even better than radioimmunoassay^{9,10}. The indirect immunofluorescent test using *Crithidia luciliae* as the antigenic substrate is a simple and specific method for detecting anti-nDNA antibodies.¹¹ *Crithidia luciliae* contains a kinetoplast, an organelle consisting of compact circular DNA which reacts bright apple-green when the test specimen is positive.¹²

PRINCIPLES OF PROCEDURE

In the indirect immunofluorescence method used in this kit, patients' sera are incubated on smears of *Crithidia luciliae* to allow binding of antibodies to the substrate. Any antibodies not bound are removed by rinsing. Bound antibodies of the IgG class are detected by incubation of the substrate with fluorescein-labeled, anti-human IgG conjugate. When observed under a fluorescence microscope equipped with appropriate filters, positive reactions appear as apple green fluorescence of the kinetoplast with or without associated nuclear staining^{12,13}. The titer, which is the reciprocal of the highest dilution giving a positive reaction, is determined by testing serial dilutions¹³.

PRODUCT INFORMATION

Storage and preparation

Store all reagents at 2-8°C. Reagents are ready for use after equilibration to room temperature.


Materials provided

ImmuGlo™ nDNA Kit	REF 1106
ImmuGlo™ nDNA Kit	REF 1106-2
ImmuGlo™ nDNA Kit	REF 1106-6

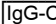
Kits contain sufficient reagents to perform 48 determinations each.

8 x	SORB SLD 6	6 well Substrate Slides, <i>Crithidia luciliae</i> (REF 1106)
16 x	SORB SLD 6	6 well Substrate Slides, <i>Crithidia luciliae</i> (REF 1106-2)
20 x	SORB SLD 6	6 well Substrate Slides, <i>Crithidia luciliae</i> (REF 1106-6)
1 x 0.5 ml	CONTROL + nDNA	nDNA Positive Control. Contains human serum.
1 x 0.5 ml	CONTROL -	Negative Control. Contains human serum.
1 x 5 ml	IgG-CONJ FITC EB	Anti-human IgG FITC Conjugate containing Evan's Blue. Protect from light. (2x REF 1106-2, 1106-6)
1 x 60 ml	BUF	Buffered Diluent.
2 vials	BUF WASH	Phosphate Buffered Saline (PBS). Dissolve each vial to 1 liter.
1 x 5.0 ml	MOUNTING MEDIUM	Mounting Medium. Do not freeze. (2x REF 1106-2, 1106-6)
1 x 12	COVER SLD	Coverslips. (2x REF 1106-2, 1106-6)

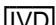EN
Optional Components

5 ml	 IgG-CONJ FITC	Anti-human IgG FITC conjugate. Protect from light.  2100X
1 ml	 EVANS	Evan's blue counterstain*.  2510.

Symbols used on labels:

 LOT	Lot number
 REF	Catalog number
 IVD	In vitro diagnostic use
	Use by
	Storage temperature
	Consult instructions for use
	Number of tests
	Manufacturer
	Date of Manufacture

 * Danger. May cause cancer. Obtain special instructions before use. Do not handle until all safety precautions have been read and understood. IF exposed or concerned: Get medical advice/attention. Wear protective gloves, protective clothing, eye protection, face protection. Store locked up. Dispose of contents/container to comply with local, state and federal regulations.

Material required but not provided

- Fluorescence microscope
- Micropipette or Pasteur pipette
- Serological pipettes
- Staining dish (e.g. Coplin jar)
- Small test tubes (e.g. 13 x 75 mm) and test tube rack
- Distilled or deionized water
- 1 liter container
- Wash bottle
- Paper towels
- Incubation chamber

WARNINGS AND PRECAUTIONS

For *in vitro* Diagnostic Use. AH human derived components used have been tested for HbsAg, HCV, HIV-1 and 2 and HTLV-I and found negative by FDA required tests. All human serum specimens and human derived products should be treated as potentially hazardous, regard-less of their origin. Follow good laboratory practices in storing, dispensing and disposing of these materials¹⁴.

WARNING – Sodium azide (NaN₃) may react with lead and copper plumbing to form highly explosive metal azides. Upon disposal of liquids, flush with large volumes of water to prevent azide buildup. Sodium azide may be toxic if ingested. If ingested, report incident immediately to laboratory director or poison control center.

Instructions should be followed exactly as they appear in this insert to ensure valid results. Do not interchange kit components with those from other sources other than the same catalog number from Immco Diagnostics Inc. Do not use beyond expiration date.

SPECIMEN COLLECTION AND PREPARATION

Only serum specimens should be used for this procedure. Grossly hemolyzed, lipemic or microbially contaminated specimens may interfere with the performance of this test and should not be used. Store specimens at 2-8°C for no longer than one week. For longer Storage, serum should be frozen at -20°C. Avoid repeated freezing and thawing of samples.

PROCEDURE

Test Method

A. Screening

1. Dilute each patient serum 1:10 with the Buffered Diluent provided (10 µl serum + 90 µl Diluent). Do not dilute Positive or Negative Controls. Save the undiluted sera to determine antibody titers if screening tests are found to be positive.

EN

2. Allow pouches containing substrate slides to equilibrate to room temperature for 10-15 minutes. Carefully remove the slides without touching the substrate.
3. Label the slides and place them in an incubation chamber lined with paper towels moistened with water to prevent drying.
4. Invert dropper vial and gently squeeze to apply 1 drop (approximately 50 μ l) of the Negative Control to well #1. Similarly apply 1 drop of Positive Control to well #2. Avoid overfilling the wells.
5. Using a micropipette or Pasteur Pipette, apply 1 drop of patient's diluted serum (approximately 50 μ l) to the other wells. Avoid overfilling the wells.
6. Place the lid on the incubation chamber and incubate slides 30 minutes at room temperature.
7. Remove a slide from the incubation chamber. Hold slide at tab end and rinse gently with approximately 10 ml PBS using a pipette, or rinse slide in beaker filled with PBS. Do not use wash bottle. Transfer slide immediately into Coplin jar and wash 10 minutes. Repeat process with all remaining slides.
8. Remove slide(s) from Coplin jar. Blot the edge of the slide on a paper towel to remove excess PBS. Place the slide in the incubation chamber. Immediately invert the Conjugate dropper vial and gently squeeze to apply 1 drop (approximately 50 μ l) to each well.
9. Repeat steps **7 and 8** for each slide.
10. Replace the lid on the incubation chamber. Incubate 30 minutes at room temperature.
11. Remove a slide from incubator. Hold the slide at the tab end and dip the slide in a beaker containing PBS to remove excess conjugate. Place slide(s) in a staining dish filled with PBS for 10 minutes. If optional conjugate without counterstain is used (see optional components in Materials Provided Section), 2-3 drops of Evan's Blue counterstain may be added to the final wash. Repeat for the remaining slides. NOTE: Improper washing may lead to increased background fluorescence.
12. Remove a slide from the staining dish. Blot the edge of the slide on a paper towel to remove excess PBS. **To prevent slide from drying, proceed immediately with next step while slide is still wet.**
13. Mount the coverslip by applying **3 drops** of Mounting Medium evenly on the coverslip and place coverslip over slide. Avoid applying undue pressure and prevent lateral movement of the coverslip.
14. Repeat steps 12 and 13 for each slide.
15. Examine for specific fluorescence under a fluorescence microscope at a magnification of 200x or greater.

Slides may be read as soon as prepared. However, because of the presence of antifading agent in the mounting medium, no significant loss of staining intensity occurs if reading is delayed for up to 48 hours. Slides should be stored in the dark at 2-8°C.

B. Endpoint Determination (titration)

A serum positive in the screening test may be further tested following steps 5 through 13 to determine the titer. Each test run should include the Positive and Negative Controls. Make serial two-fold dilutions starting at 1:10. The reciprocal of the highest dilution producing a positive reaction is the titer.

Preparation of Serial Dilutions

Number six tubes 1 through 6. Add 0.9 ml of Sample Diluent to tube 1 and 0.2 ml to tubes 2 through 6. Pipette 0.1 ml of undiluted serum to tube 1 and mix thoroughly. Transfer 0.2 ml from tube 1 to tube 2 and mix thoroughly. Continue transferring 0.2 ml from one tube to the next after mixing to yield the dilutions depicted in the following table:

Tubes	1	2	3	4	5	6
Serum	0.1 ml					
	+					
Buffered Diluent	0.9 ml	0.2 ml	0.2 ml	0.2 ml	0.2 ml	0.2 ml
		↻	↻	↻	↻	↻
Transfer	0.2 ml	0.2 ml	0.2 ml	0.2 ml	0.2 ml	
Final dilution	1:10	1:20	1:40	1:80	1:160	1:320 etc.

QUALITY CONTROL

Both a Positive and Negative Control should be included with each test run. The Negative Control should show no specific fluorescence of the kinetoplast. The Positive Control should have 2+ or greater staining intensity of this structure. If expected results are not obtained, the run should be repeated, continue to occur with the controls, these may be due to:

- Turbidity. Discard and use another control
- Problems with the optical system of the fluorescence microscope, improper alignment, bulb beyond useful life expectancy, etc.
- Allowing the slide to dry during the procedure.

INTERPRETATION OF RESULTS

The results of the tests for anti-*n*DNA antibodies should be reported as negative (<10), positive, (greater or equal to 320) or alternatively positive with specific endpoint titre.

EN

Read only fields which contain well separated *C. luciliae*. Observe for specific staining of the kinetoplast (see figure 1 at the end of this document). Staining of the nucleus or the polar body should not be interpreted as a positive DNA antibody test. Absence of specific staining of the kinetoplast is considered negative for antibodies to *n*DNA.

LIMITATION OF THE PROCEDURE

In some cases, sera positive for kinetoplast staining may either be very weak or negative at the initial screening dilution prozone phenomenon. In such doubtful cases the sera should be screened at higher dilutions and, if positive, antibody titers determined.

In rare instances false positive reactions may be observed. These may be due to the presence of high levels of lipoproteins or other proteins which bind to DNA^{12,13}.

The goat anti-human IgG FITC Conjugate supplied in this kit is primarily heavy chain specific but has some light chain activity. It reacts primarily with IgG class autoantibodies, but may, to a lesser degree, react with light chains of other classes such as IgM.

The clinician should consider the results of all positive indirect immunofluorescence tests along with the results of other laboratory tests and the clinical condition of the patient when making a diagnosis.

EXPECTED VALUES

As seen in Table 1 at the end of this document, anti-*n*DNA antibodies are not detected (titer <10) in sera from normal subjects or patients with scleroderma or rheumatoid arthritis. Anti-*n*DNA antibodies occur (titer >10) in over one half of patients with SLE.

PERFORMANCE CHARACTERISTICS

The ImmuGlo™ *n*DNA Antibody test was compared with another commercially available fluorescent antibody test using *C. luciliae* as a substrate. The comparison included 106 serum samples from normal subjects as well as from patients with the diagnosis of SLE, scleroderma, or rheumatoid arthritis. Sera were tested according to the procedure and screening dilution recommended by the manufacturer. These yielded comparable results as summarized in Table 2.

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REFERENCES • ΒΙΒΛΙΟΓΡΑΦΙΑ • LITERATUR • BIBLIOGRAPHIE • BIBLIOGRAFIA

1. Feltkamp TEW (Ed). The significance of the determination of anti-DNA and DNA/anti-DNA complexes. Scand J Rheumatol Suppl. 1975;11:7-64.
2. Tan EM. Antinuclear antibodies: diagnostic markers for autoimmune diseases and probes for cell biology. Adv Immunol. 1989; 44:93-151.
3. Tzioufas AG, Manoussakis MN, Drosos AA et al. Enzyme immunoassays for the detection of IgG and IgM anti-dsDNA antibodies: clinical significance and specificity. Clin Exp Rheumatol. 1987; 5:247-253.
4. Ruffati A, Calligaro A, Ross D et al. Anti-double stranded DNA antibodies in the healthy elderly: prevalence and characteristics. J Clin Immunol. 1990; 10:300-303.
5. Kadlubowski M, Jackson M, Yap PL and Neill G. Lack of specificity for antibodies to double stranded DNA found in four commercial kits. J Clin Path. 1991; 44:246-250.
6. Brinkman K, Termaat R, Van den Brink H et al. The specificity of the anti-dsDNA ELISA: a closer look. J Immunol Methods. 1991; 13:91-100.
7. Lange A. Evaluation of the simultaneous estimation of anti-dsDNA and anti-ssDNA antibodies for clinical purposes. Clin Exp Immunol. 1978; 31:472-481.
8. Borg EJ, Horst G, Hum EJ et al. Measurement of increases in antidouble stranded DNA antibody level as a predictor of disease exacerbation in systematic lupus erythematosus: a long term, prospective study. Arth Rheum. 1990; 33:634-643.
9. Aarden LA, Lakmaker F, de Groot ER, et al.: Detection of antibodies to DNA by radioimmunoassay and immunofluorescence. Scand. J. Rheu. 1975; Suppl. 11:12-19.
10. Fish F, Ziff M: A sensitive solid phase microradioimmunoassay for anti double-stranded DNA antibodies. Arthritis and Rheumatism 1981; 24:534-543.
11. Aarden LA, de Groot ER, Feltkamp TE: Immunology of DNA. III. *Crithidia luciliae*, a simple substrate for the detection of anti-dsDNA with immunofluorescence techniques. Ann. NY Acad. Sei. 1975; 254: 505-515.
12. Crowe W, Kusher I: An immunofluorescent method using *Crithidia luciliae* to detect antibodies to double-stranded DNA. Arthritis and Rheumatism 1977; 20:811-814.
13. Feltkamp TE, Kirkwood TB et al. The first international standard for antibodies to double stranded DNA. Ann Rheum Dis. 1988; 47:740-746.
14. Biosafety in Microbiological and Biomedical Laboratories. Centers for Disease Control/National Institute of Health, Fourth Edition, 1999, (HHS Pub. # (CDC) 93-8395).

Figure 1.C. *Luciliae*

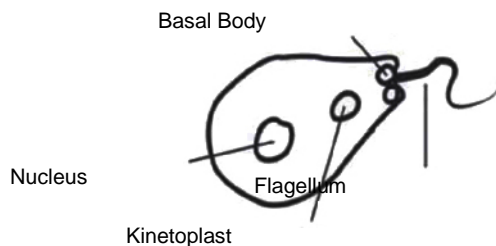
Table 1. Incidence of Antibodies to *n*DNA in Collagen-Vascular Disorders as detected by Indirect Immunofluorescence on *C. luciliae*

Clinical Condition	No. Tested	No. Positive	% Positive
SLE	28	19	68
Scleroderma	23	0	0
Rheumatoid Arthritis	8	0	0
Normal Controls	106	0	0

Table 2. Comparison of Kits using *C. luciliae* Substrate for the Detection of Antibodies to *n*DNA

Clinical Condition	n	Immco™		Other	
		Positive	% Positive	Positive	% Positive
SLE	28	19	68	13	46
Scleroderma	23	0	0	0	0
Rheumatoid Arthritis	8	0	0	0	0
Normal Controls	106	0	0	0	0

For technical assistance please contact:

 **IMMCO Diagnostics, Inc.**
60 Pineview Drive
Buffalo, NY 14228-2120 USA
Telephone: (716) 691-0091
Fax: (716) 691-0466
Toll Free USA/Canada: 1-800-537-TEST
E-Mail: info@immco.com
or your local product distributor

EMERGO Europe.
Prinsessegracht 20
2514 AP The Hague
The Netherlands