

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in





CALR & B2M Protein Protein Interaction Antibody Pair

Catalog #: DI0066 規格:[1 Set]

List All

Specification

Product Description:


This protein protein interaction antibody pair set comes with two antibodies to detect the protein-protein interaction, one against the CALR protein, and the other against the B2M protein for use in *in situ* Proximity Ligation Assay. See Publication Reference below.

Application Image


In situ Proximity Ligation Assay (Cell)




1. Incubate with target primary antibodies

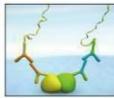

2. Add PLA probes PLUS and MINUS

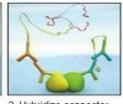

3. Hybridize connector oligos

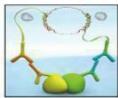

4. Ligation to form a complete DNA circle

Human



5. Rolling circle amplification




6. Add fluorescent probes to reveal interaction

Reactivity:

Quality Control Protein protein interaction immunofluorescence result.

Testing:

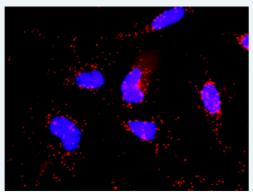

Representative image of Proximity Ligation Assay of protein-protein interactions between CALR and B2M. HeLa cells were stained with anti-CALR rabbit purified polyclonal antibody 1:1200 and anti-B2M mouse monoclonal antibody 1:50. Each red dot represents the detection of protein-protein interaction complex. The images were analyzed using an optimized freeware (BlobFinder) download from The Centre for Image Analysis at Uppsala University.

Supplied **Product:**

Antibody pair set content:

- 1. CALR rabbit purified polyclonal antibody (20 ug)
- 2. B2M mouse monoclonal antibody (40 ug)

*Reagents are sufficient for at least 30-50 assays using recommended protocols.

Storage Instruction:

Store reagents of the antibody pair set at -20°C or lower. Please aliquot to avoid repeated freeze thaw cycle. Reagents should be returned to -

20°C storage immediately after use.

MSDS:

<u>m</u>Download

Publication Reference

1. An analysis of protein-protein interactions in cross-talk pathways reveals CRKL as a novel prognostic marker in hepatocellular carcinoma. Liu CH, Chen TC, Chau GY, Jan YH, Chen CH, Hsu CN, Lin KT, Juang YL, Lu PJ, Cheng HC, Chen MH, Chang CF, Ting YS, Kao CY, Hsiao M, Huang CY. Mol Cell Proteomics. 2013 Feb 8. [Epub ahead of print]

Applications

In situ Proximity Ligation Assay (Cell)

B2M CALR

Gene Information

Entrez GeneID: 811

Gene Name: CALR

Gene Alias: CRT,FLJ26680,RO,SSA,cC1qR

calreticulin Gene

Description:

Omim ID: 109091

Gene Ontology: Hyperlink

Gene Summary: Calreticulin is a multifunctional protein that acts as a major Ca(2+)binding (storage) protein in the lumen of the endoplasmic reticulum. It is also found in the nucleus, suggesting that it may have a role in transcription regulation. Calreticulin binds to the synthetic peptide KLGFFKR, which is almost identical to an amino acid sequence in the DNA-binding domain of the superfamily of nuclear receptors. Calreticulin binds to antibodies in certain sera of systemic lupus and Sjogren patients which contain anti-Ro/SSA antibodies, it is highly conserved among species, and it is located in the endoplasmic and sarcoplasmic reticulum where it may bind calcium. The amino terminus of calreticulin interacts with the DNA-binding domain of the glucocorticoid receptor and prevents the receptor from binding to its specific glucocorticoid response element. Calreticulin can inhibit the binding of androgen receptor to its hormone-responsive DNA element and can inhibit androgen receptor and retinoic acid receptor transcriptional activities in vivo, as well as retinoic acid-induced neuronal differentiation. Thus, calreticulin can act as an important modulator of the regulation of gene transcription by nuclear hormone receptors. Systemic lupus erythematosus is associated with increased autoantibody titers against calreticulin but calreticulin is not a Ro/SS-A antigen. Earlier papers referred to calreticulin as an Ro/SS-A antigen but this was later disproven. Increased autoantibody titer against human calreticulin is found in infants with complete congenital heart block of both the IgG and IgM classes. [provided by RefSeq

Other **Designations:** Sicca syndrome antigen A (autoantigen Ro; calreticulin), autoantigen Ro

Gene Information

Entrez GeneID: 567

Gene Name: B₂M

> Page 2 of 4 2016/5/19

Gene Alias:

Gene

beta-2-microglobulin

Description:


Omim ID:

<u>109700</u>, <u>241600</u>

Gene Ontology: Hyperlink


Gene Summary: This gene encodes a serum protein found in association with the major histocompatibility complex (MHC) class I heavy chain on the surface of nearly all nucleated cells. The protein has a predominantly beta-pleated sheet structure that can form amyloid fibrils in some pathological conditions. A mutation in this gene has been shown to result in hypercatabolic hypoproteinemia

Other **Designations:** beta chain of MHC class I molecules, beta-2-microglobin

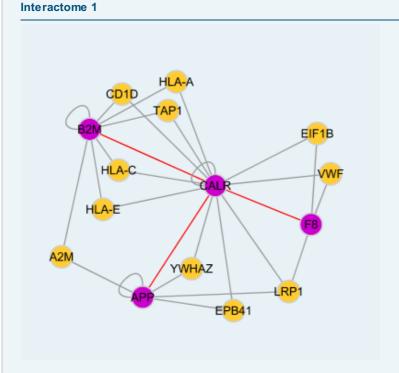

Interactome 2

Page 3 of 4 2016/5/19

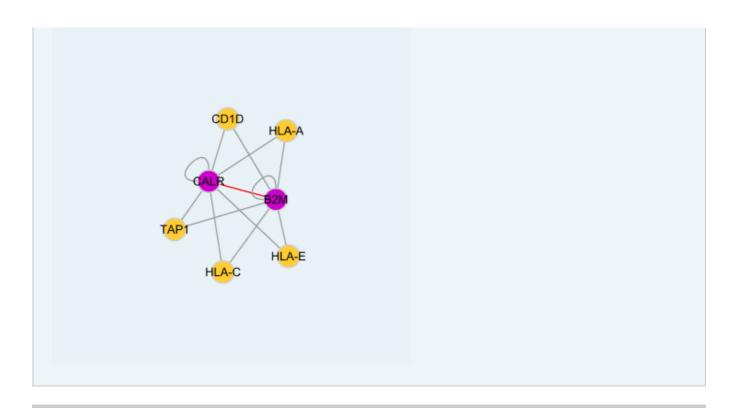

服務條款 | 隱私權政策 | 著作及商標 | 網站地圖

©2016 亞諾法生技股份有限公司 Abnova Corporation. 版權所有.

Page 4 of 4 2016/5/19