

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in





OSM & COL4A6 Protein Protein Interaction Antibody Pair

Catalog #: DI0105 規格:[1 Set]

List All

Specification

Product Description:

This protein protein interaction antibody pair set comes with two antibodies to detect the protein-protein interaction, one against the OSM protein, and the other against the COL4A6 protein for use in *in situ* Proximity Ligation Assay. See Publication Reference below.

Application Image

In situ Proximity Ligation Assay (Cell)




 Incubate with target primary antibodies



2. Add PLA probes PLUS and MINUS



Hybridize connector oligos

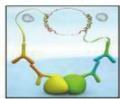

 Ligation to form a complete DNA circle



Rolling circle amplification




Add fluorescent probes
to reveal interaction

Reactivity:

Human

Quality Contro Testing:

Quality Control Protein protein interaction immunofluorescence result.

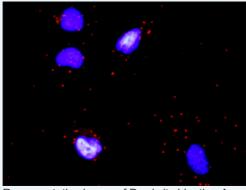

Representative image of Proximity Ligation Assay of protein-protein interactions between OSM and COL4A6. HeLa cells were stained with anti-OSM rabbit purified polyclonal antibody 1:1200 and anti-COL4A6 mouse monoclonal antibody 1:50. Each red dot represents the detection of protein-protein interaction complex. The images were analyzed using an optimized freeware (BlobFinder) download from The Centre for Image Analysis at Uppsala University.

Supplied Product:

Antibody pair set content:

- 1. OSM rabbit purified polyclonal antibody (20 ug)
- 2. COL4A6 mouse monoclonal antibody (40 ug)

*Reagents are sufficient for at least 30-50 assays using recommended protocols.

Storage Instruction:

Store reagents of the antibody pair set at -20°C or lower. Please aliquot to avoid repeated freeze thaw cycle. Reagents should be returned to -



Publication Reference

 An analysis of protein-protein interactions in cross-talk pathways reveals CRKL as a novel prognostic marker in hepatocellular carcinoma.
Liu CH, Chen TC, Chau GY, Jan YH, Chen CH, Hsu CN, Lin KT, Juang YL, Lu PJ, Cheng HC, Chen MH, Chang CF, Ting YS, Kao CY, Hsiao M, Huang CY. Mol Cell Proteomics. 2013 Feb 8. [Epub ahead of print]

Applications

In situ Proximity Ligation Assay (Cell)

COL4A6 OSM

Gene Information

Entrez GeneID: 5008

Gene Name: OSM

Gene Alias: MGC20461

Gene oncostatin M

Description:

Omim ID: <u>165095</u>

Gene Ontology: Hyperlink

Gene Summary: Oncostatin M is a member of a cytokine family that includes leukemia-

inhibitory factor, granulocyte colony-stimulating factor, and interleukin 6. This gene encodes a growth regulator which inhibits the proliferation of a number of tumor cell lines. It regulates cytokine production, including IL-6, G-CSF and GM-CSF from endothelial cells. [provided by RefSeq

Other Designations:

Gene Information

Entrez GeneID: 1288

Gene Name: COL4A6

Gene Alias: MGC88184

Gene collagen, type IV, alpha 6

Description:

Omim ID: <u>303631</u>, <u>308940</u>


Gene Ontology: Hyperlink

Gene Summary: This gene encodes one of the six subunits of type IV collagen, the major

structural component of basement membranes. Like the other members of the type IV collagen gene family, this gene is organized in a head-to-head conformation with another type IV collagen gene, alpha 5 type IV collagen, so that the gene pair shares a common promoter. Deletions in the alpha 5 gene that extend into the alpha 6 gene result in diffuse leiomyomatosis accompanying the X-linked Alport syndrome caused by the deletion in the alpha 5 gene. Two splice variants have been

identified for this gene. [provided by RefSeq

Page 2 of 3 2016/5/19

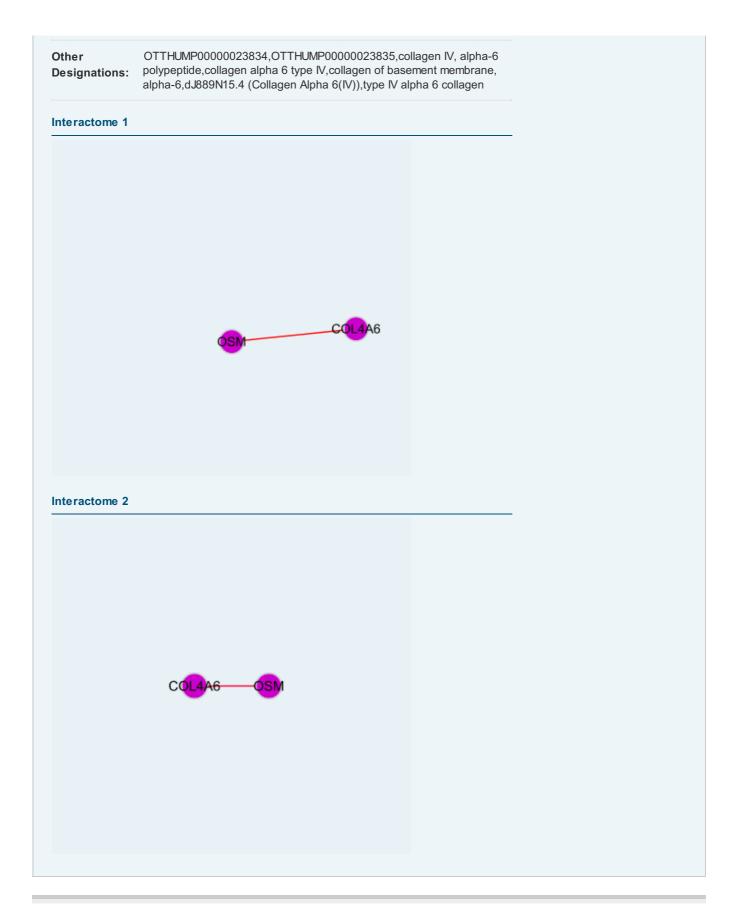

服務條款 | 隱私權政策 | 著作及商標 | 網站地圖

©2016 亞諾法生技股份有限公司 Abnova Corporation. 版權所有.

Page 3 of 3 2016/5/19