



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

CELLBANKER® I

Cryopreservation Medium (Serum Containing)

Cat # I1884

Qty: 4 x 20ml

Expiry Date: 3 years from manufacturing date (see label)




Cell-Freezing:

For optimum results, cells for cryopreservation should be in log phase of growth. Similar or standard freezing protocols may be substituted.

1. Examine and make sure the cell culture is free of contamination, in healthy situation and proper confluency, etc.
2. Perform a cell count to determine the viability of cells
3. Gently pellet the cells by centrifugation (3 - 5 minutes at 1,000~2,000rpm, 4°C). Remove the supernatant by using an aspirator.
4. Gently suspend the cells with CELLBANKER® I cryopreservation medium (1 ml for 5×10^5 - 5×10^6 cells).
5. Dispense the cell suspension in 1ml aliquots to cryopreservation vials that have been labeled with the cell line name, cell concentration, passage date and other essential information.
6. Place the vials directly in a -80°C for storage. If necessary, transfer the frozen vials to a liquid nitrogen storage tank after the vials have been frozen for at least 24 hours.
7. Optimum protocol may change with the cell types.

IMPORTANT: Optimum protocol may change with the cell types.

Procedure for Use:

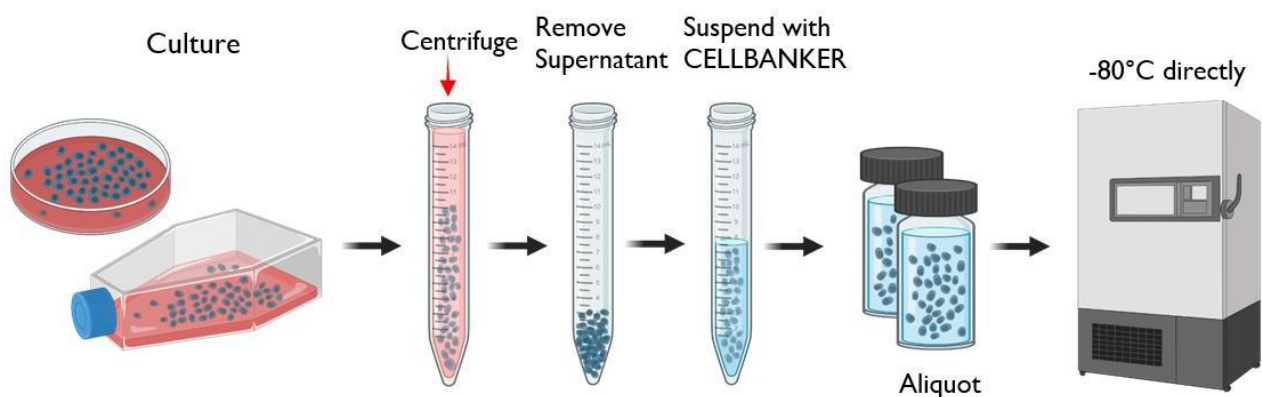

Thawing:

1. Remove the frozen cell from storage and quickly thaw in a 37°C shaking water bath.
2. Immediately dilute and gently mix each 1ml of cells with 10ml of complete cell culture medium.
3. Gently pellet the cells by centrifugation (3-5 minutes at 1,000 - 2,000rpm, 4°C). Remove the supernatant by aspirator.
4. Gently suspend the cells with appropriate volume of complete cell culture medium and plate in a culture flask.
5. Continue the further culture procedures according to standard protocols.

Guarantee of Quality:

1. Bacterial contamination free – Product has been tested and confirmed to be free of bacteria, fungi and mycoplasma.
2. Chemical Analysis: pH (7.0 to 8.5 at room temperature) Endotoxin (<5 EU/mL)
3. Performance test – Cell viability above 80% (JM404, SK-007) is guaranteed.

Storage of CELLBANKER® I:

1. CELLBANKER® I should be stored at 4°C or below.
2. For long-term storage CELLBANKER® I can be frozen. Repeated freezing and thawing may impair the quality of the product; it is recommended that CELLBANKER® I is aliquoted before freezing.

Precautions:

1. For research use only
2. Not for clinical or diagnostic use.
3. Performance of trial tests using cells of intended use before experiments is recommended.

Product Range:

Description	Pack Size
CELLBANKER® I – Serum Containing	20 ml
CELLBANKER® I – Serum Containing	4 x 20 ml
CELLBANKER® I – Serum Containing	100 ml
CELLBANKER® 2 – Serum Free	20 ml
CELLBANKER® 2 – Serum Free	4 x 20 ml
CELLBANKER® 2 – Serum Free	100 ml
STEM-CELLBANKER® - GMP	20 ml
STEM-CELLBANKER® - GMP	4 x 20 ml
STEM-CELLBANKER® - GMP	100 ml
STEM-CELLBANKER® - GMP - DMSO Free	20 ml
STEM-CELLBANKER® - GMP - DMSO Free	4 x 20 ml
STEM-CELLBANKER® - GMP - DMSO Free	100 ml
STEM-CELLBANKER® EX - GMP	100 ml
CELLOTION cell wash solution	100 ml

Citations:

Hata, M., Omi, M., Kobayashi, Y. et al. Transplantation of human dental pulp stem cells ameliorates diabetic polyneuropathy in streptozotocin-induced diabetic nude mice: the role of angiogenic and neurotrophic factors. Stem Cell Res Ther 11, 236 (2020). <https://doi.org/10.1186/s13287-020-01758-9>

Typical Experimental Results:

Cell Type	Preservation period (year)	Viability of cells (%)	
		-80°C	-196°C
Mouse			
Hybridoma	10	95	95
Myeloma	10	90	90
L929	10	90	90
FM3A	5	90	90
BALB/3T3	5	90	90
MI	5	90	90
YAC-I	5	90	-
Rat			
RLC-16	5	90	90
NRK	5	90	90
PC-12	5	90	-
Hamster			
CHO	5	90	90
V79	5	90	90
Monkey			
COS-I	5	90	90
Vero	5	90	90
Human			
Kidney-derived tumor cell	5	90	90
EBV transformed cell	5	90	90
HEL-derived fibroblast	5	90	90
Melanoma	5	90	90
Caco-2	3	90	-
C-5	5	90	90
CEM	5	90	90
K562	10	90	90
Jurkat	10	90	90
BALL-I	5	90	90
HUC-Fm	5	80	80



AMS BIOTECHNOLOGY (EUROPE) LTD
UK & REST OF THE WORLD
 184 Park Drive, Milton Park
 Abingdon, OX14 4SE
T: +44 (0) 1235 828 200
F: +44 (0) 1235 820 482



AMSBIO LLC
USA & CANADA
 1035 Cambridge Street,
 Cambridge, MA 02141
T: +1 (617) 945-5033 or
 +1 (800) 987-0985
F: +1 (617) 945-8218



AMS BIOTECHNOLOGY (EUROPE) LTD
SWITZERLAND
 Via Lisano 3, (Cp.683)
 CH-6900 Messagno
T: +41 (0) 91 604 55 22
F: +41 (0) 91 605 17 85



AMSBIO EUROPE BV
EU
 Berenkoog 41,
 1822 BH Alkmaar,
 Netherlands
T: +31 (0) 72 8080244
F: +31 (0) 72 8080142