



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

CELLOTION®

Cell wash and recovery solution

Cat # 11918 (previously 11898)

Qty: 100ml


Expiry Date: 3 years from manufacturing date (see label)



CELLOTION® is a substitution for general-use cell wash solution such as PBS buffer, saline or culture medium with standards cell wash procedures.

Example of use: Monocytes isolation with CELLOTION® Cell Wash and Recovery Solution.

1. Gently dilute the peripheral blood or umbilical cord blood (UCB) with CELLOTION® Cell Wash and Recovery Solution 1:1.
2. Gently add Ficoll-hypaque solution to the mixture 1:2. Do not mix and the blood should remain on top.)
3. Centrifuge the tube for 30minutes at 1,800rpm, at room temperature.
4. Collect the buffy coat (white blood cells ring) carefully with a sterile pipette tip without touching the Ficoll-hypaque solution.
5. Gently dilute the cells with CELLOTION®, then wash for 2 - 3 times in CELLOTION®.
6. Count the number of cells.

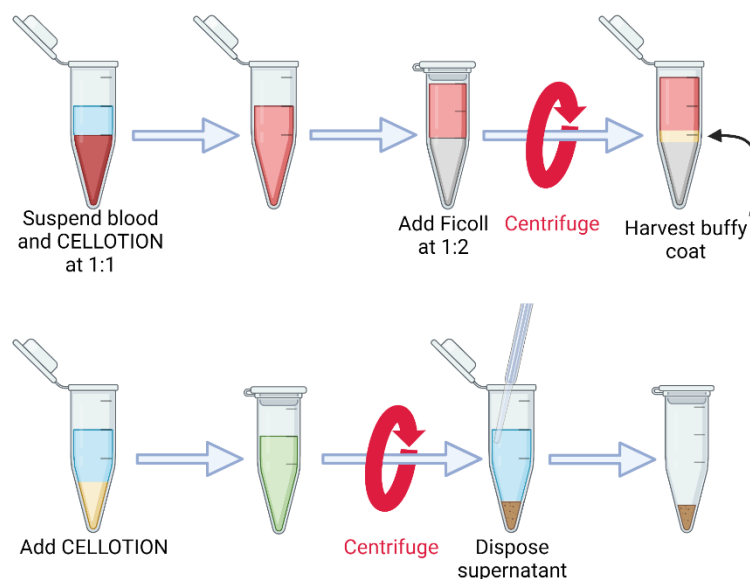

AMSBIO | www.amsbio.com | info@amsbio.com

AMSBIO LLC
USA & Canada 
1035 Cambridge Street,
Cambridge, MA 02141
T: +1 (617) 945-5033 or
T: +1 (800) 987-0985
F: +1 (617) 945-8218

AMSBIO Europe BV
EU 
Berenkoog 41,
1822 BH Alkmaar,
Netherlands
T: +31 (0) 72 8080244
F: +31 (0) 72 8080142

AMS Biotechnology (Europe) Ltd
UK & Rest of the World 
184 Park Drive, Milton Park
Abingdon OX14 4SE
T: +44 (0) 1235 828 200
F: +44 (0) 1235 820 482

AMS Biotechnology (Europe) Ltd
Switzerland 
Via Lisano 3,
(CP.683)
CH-6900 Massagno
T: +41 (0) 91 604 55 22
F: +41 (0) 91 605 17 85

AMSBIO Europe BV
Deutschland 
T: +49 (0) 69 779099
F: +49 (0) 69 13376880

Product Range:

Description	Pack Size
CELLBANKER® I – Serum Containing	20 ml
CELLBANKER® I – Serum Containing	4 x 20 ml
CELLBANKER® I – Serum Containing	100 ml
CELLBANKER® 2 – Serum Free	20 ml
CELLBANKER® 2 – Serum Free	4 x 20 ml
CELLBANKER® 2 – Serum Free	100 ml
STEM-CELLBANKER® - GMP	20 ml
STEM-CELLBANKER® - GMP	4 x 20 ml
STEM-CELLBANKER® - GMP	100 ml
STEM-CELLBANKER® - GMP - DMSO Free	20 ml
STEM-CELLBANKER® - GMP - DMSO Free	4 x 20 ml
STEM-CELLBANKER® - GMP - DMSO Free	100 ml
STEM-CELLBANKER® EX - GMP	100 ml
CELLOTION cell wash solution	100 ml

For research use only; not for use as a
diagnostic.

AMSBIO | www.amsbio.com | info@amsbio.com

AMSBIO LLC
USA & Canada 
1035 Cambridge Street,
Cambridge, MA 02141
T: +1 (617) 945-5033 or
T: +1 (800) 987-0985
F: +1 (617) 945-8218

AMSBIO Europe BV
EU 
Berenkoog 41,
1822 BH Alkmaar,
Netherlands
T: +31 (0) 72 8080244
F: +31 (0) 72 8080142

AMS Biotechnology (Europe) Ltd
UK & Rest of the World 
184 Park Drive, Milton Park
Abingdon OX14 4SE
T: +44 (0) 1235 828 200
F: +44 (0) 1235 820 482

AMS Biotechnology (Europe) Ltd
Switzerland 
Via Lisano 3,
(CP.683)
CH-6900 Massagno
T: +41 (0) 91 604 55 22
F: +41 (0) 91 605 17 85

AMSBIO Europe BV
Deutschland 
T: +49 (0) 69 779099
F: +49 (0) 69 13376880