



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

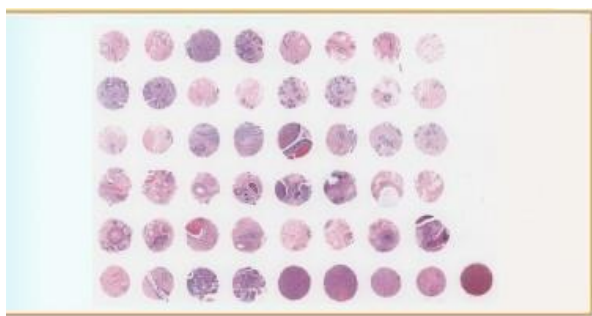
[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Datasheet

Product Code: BB08015

Product Description: Breast ductal intraepithelial neoplasia and cancer tissue array, including TNM, clinical stage and pathology grade, 24 cases/48 cores

Breast ductal intraepithelial neoplasia and cancer tissue array, including TNM, clinical stage and pathology grade, 24 cases/48 cores

Microarray Panel	Breast ductal intraepithelial neoplasia and carcinoma tissue microarray, containing 20 cases of breast ductal intraepithelial neoplasia and 4 breast invasive ductal carcinoma, duplicate cores per case												
Cores	48												
Cases	24												
Row number	6												
Column number	8												
Core Diameter (mm)	1.5												
Thickness (µm)	5		1	2	3	4	5	6	7	8			
Tissue Array Type	FFPE	A	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre		
Species	Human	B	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre		
		C	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre		
		D	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre		
		E	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre		
		F	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Bre	Adr		
Applications	Routine histology procedures including Immunohistochemistry (IHC) and In Situ Hybridization (ISH), protocols which can be found at our support page.												
Notes	1. TMA slides were sectioned and stored at 4°C and may not be fresh cut, but still suitable for IHC. Please request fresh cut if experiment involves phospho-specific antibodies, RNA studies, FISH or ISH, etc. A minimum of 3 slides per TMA must be purchased to cover the cost of trimming for fresh sectioning. 2. Most TMA slides were not coated with an extra layer of paraffin (tissue cores can be easily seen on the glass). To prevent tissue detachment during antigen retrieval, unbaked slides must be baked for at least 30 to 120 minutes at 60°C. before putting into xylene for de-paraffinization. Baked slides were sent out baked for 2 hours. In the following specsheet, “*” means invalid core; “-” means no applicable or negative in IHC markers.												

AMSBIO | www.amsbio.com | info@amsbio.com



UK & Rest of the World
184 Park Drive, Milton Park
Abingdon OX14 4SE, UK
T: +44 (0)1235 828 200
F: +44 (0) 1235 820 482



North America
1035 Cambridge Street,
Cambridge, MA 02141
T: +1 (617) 945-5033 or
T: +1 (800) 987-0985
F: +1 (617) 945-8218



Germany
Bockenheimer Landstr. 17/19
60325 Frankfurt/Main
T: +49 (0) 69 779099
F: +49 (0) 69 13376880



Switzerland
Centro Nord-Sud 2E
CH-6934 Bioggio-Lugano
T: +41(0) 91 604 55 22
F: +41(0) 91 605 17 85

Pos.	No.	Age	Sex	Organ/ Anatomic Site	Pathology diagnosis	TNM	Grade	Stage	Type
A1	1	30	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
A2	2	30	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
A3	3	50	F	Breast	Ductal intraepithelial neoplasia (DIN 1b)	-	-	-	DIN
A4	4	50	F	Breast	Ductal intraepithelial neoplasia (DIN 1b)	-	-	-	DIN
A5	5	50	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
A6	6	50	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
A7	7	56	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
A8	8	56	F	Breast	Ductal intraepithelial neoplasia	-	-	-	DIN
B1	9	54	F	Breast	Ductal intraepithelial neoplasia (DIN 1b)	-	-	-	DIN
B2	10	54	F	Breast	Ductal intraepithelial neoplasia (DIN 1b)	-	-	-	DIN
B3	11	48	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
B4	12	48	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
B5	13	22	F	Breast	Ductal intraepithelial neoplasia (DIN 1b)	-	-	-	DIN
B6	14	22	F	Breast	Ductal intraepithelial neoplasia (DIN 1b)	-	-	-	DIN
B7	15	37	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
B8	16	37	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
C1	17	49	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
C2	18	49	F	Breast	Ductal intraepithelial neoplasia (DIN 1a)	-	-	-	DIN
C3	19	49	F	Breast	Invasive ductal carcinoma	T2N1 M0	1	IIB	Malignant
C4	20	49	F	Breast	Invasive ductal carcinoma	T2N1 M0	1	IIB	Malignant
C5	21	42	F	Breast	Ductal intraepithelial neoplasia (DIN 2)	-	-	-	DIN
C6	22	42	F	Breast	Ductal intraepithelial neoplasia (DIN 2)	-	-	-	DIN
C7	23	38	F	Breast	Ductal intraepithelial neoplasia (DIN 1b)	-	-	-	DIN
C8	24	38	F	Breast	Ductal intraepithelial neoplasia (DIN 1b)	-	-	-	DIN
D1	25	36	F	Breast	Ductal intraepithelial neoplasia (DIN 2)	-	-	-	DIN
D2	26	36	F	Breast	Ductal intraepithelial neoplasia (DIN 2)	-	-	-	DIN

D3	27	52	F	Breast	Ductal intraepithelial neoplasia (DIN 3)	-	-	-	DIN
D4	28	52	F	Breast	Ductal intraepithelial neoplasia (DIN 3)	-	-	-	DIN
D5	29	44	F	Breast	Ductal intraepithelial neoplasia (DIN 2)	-	-	-	DIN
D6	30	44	F	Breast	Ductal intraepithelial neoplasia (DIN 2)	-	-	-	DIN
D7	31	44	F	Breast	Ductal intraepithelial neoplasia (DIN 1b)	-	-	-	DIN
D8	32	44	F	Breast	Ductal intraepithelial neoplasia (DIN 1b)	-	-	-	DIN
E1	33	46	F	Breast	Ductal intraepithelial neoplasia (DIN 2)	-	-	-	DIN
E2	34	46	F	Breast	Ductal intraepithelial neoplasia (DIN 2)	-	-	-	DIN
E3	35	67	F	Breast	Ductal intraepithelial neoplasia (DIN 2)	-	-	-	DIN
E4	36	67	F	Breast	Ductal intraepithelial neoplasia (DIN 2)	-	-	-	DIN
E5	37	40	F	Breast	Ductal intraepithelial neoplasia (DIN 3)	-	-	-	DIN
E6	38	40	F	Breast	Ductal intraepithelial neoplasia (DIN 3)	-	-	-	DIN
E7	39	43	F	Breast	Ductal intraepithelial neoplasia (DIN 3)	-	-	-	DIN
E8	40	43	F	Breast	Ductal intraepithelial neoplasia (DIN 3)	-	-	-	DIN
F1	41	30	F	Breast	Ductal intraepithelial neoplasia (DIN 3)	-	-	-	DIN
F2	42	30	F	Breast	Ductal intraepithelial neoplasia (DIN 3)	-	-	-	DIN
F3	43	63	F	Breast	Invasive ductal carcinoma	T2N0 M0	1	IIA	Malignant
F4	44	63	F	Breast	Invasive ductal carcinoma	T2N0 M0	1	IIA	Malignant
F5	45	54	F	Breast	Invasive ductal carcinoma	T1N0 M0	2	I	Malignant
F6	46	54	F	Breast	Invasive ductal carcinoma	T1N0 M0	2	I	Malignant
F7	47	42	F	Breast	Invasive ductal carcinoma	T2N1 M0	3	IIB	Malignant
F8	48	42	F	Breast	Invasive ductal carcinoma	T2N1 M0	3	IIB	Malignant

AMSBIO | www.amsbio.com | info@amsbio.com

 **UK & Rest of the World**
184 Park Drive, Milton Park
Abingdon OX14 4SE, UK
T: +44 (0)1235 828 200
F: +44 (0) 1235 820 482

 **North America**
1035 Cambridge Street,
Cambridge, MA 02141
T: +1 (617) 945-5033 or
T: +1 (800) 987-0985
F: +1 (617) 945-8218

 **Germany**
Bockenheimer Landstr. 17/19
60325 Frankfurt/Main
T: +49 (0) 69 779099
F: +49 (0) 69 13376880

 **Switzerland**
Centro Nord-Sud 2E
CH-6934 Bioggio-Lugano
T: +41(0) 91 604 55 22
F: +41(0) 91 605 17 85