



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

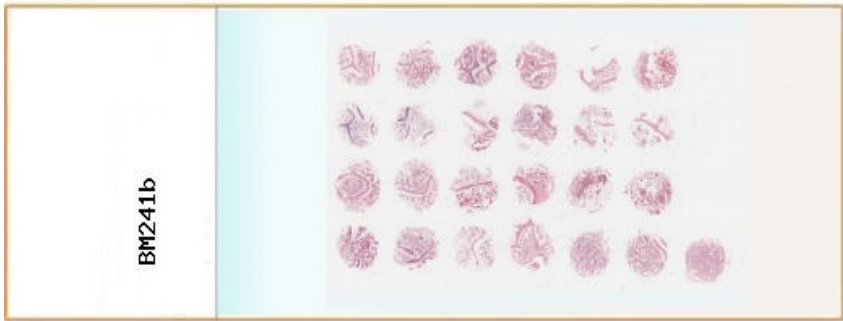
F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

BM241b : Bone marrow tissue array, 12 cases/24 cores, replacing BM241a

Microarray Panel	Bone marrow tissue microarray, containing 12 cases of bone marrow tissue, duplicated cores per case				
Cores	24	BM241b H&E 			
Cases	12				
Row number	4				
Column number	6				
Core Diameter (mm)	1.5				
Thickness (µm)	5				
Tissue Array Type	FFPE				
Species	Human				
Applications	Routine histology procedures including Immunohistochemistry (IHC) and In Situ Hybridization (ISH), protocols which can be found at our support page.				
Notes	1. TMA slides were sectioned and stored at 4°C and may not be fresh cut, but still suitable for IHC. Please request fresh cut if experiment involves phospho-specific antibodies, RNA studies, FISH or ISH, etc. A minimum of 3 slides per TMA must be purchased to cover the cost of trimming for fresh sectioning. 2. Most TMA slides were not coated with an extra layer of paraffin (tissue cores can be easily seen on the glass). To prevent tissue detachment during antigen retrieval, unbaked slides must be baked for at least 30 to 120 minutes at 60°C. before putting into xylene for de-paraffinization. Baked slides were sent out baked for 2 hours. In the following specsheet, "*" means invalid core; "-" means no applicable or negative in IHC markers.				

	1	2	3	4	5	6
A	Rib	Rib	Rib	Rib	Rib	Rib
B	Rib	Rib	Rib	Rib	Rib	Rib
C	Rib	Rib	Rib	Rib	Rib	Rib
D	Rib	Rib	Rib	Rib	Rib	Adr

Legend:

Rib - Rib

 - Malignant tumor,  - NAT

, tissue IDs are available in exported Excel files.

Pos.	No.	Age	Sex	Organ/Anatomic Site	Pathology diagnosis	TNM	Grade	Stage	Type	Image
A1	1	45	M	Rib	Bone marrow tissue				NAT	
A2	2	45	M	Rib	Bone marrow tissue				NAT	
A3	3	62	M	Rib	Bone marrow tissue				NAT	
A4	4	62	M	Rib	Bone marrow tissue				NAT	
A5	5	71	M	Rib	Bone marrow tissue				NAT	
A6	6	71	M	Rib	Bone marrow tissue				NAT	
B1	7	36	M	Rib	Bone marrow tissue				NAT	
B2	8	36	M	Rib	Bone marrow tissue				NAT	
B3	9	41	M	Rib	Bone marrow tissue				NAT	
B4	10	41	M	Rib	Bone marrow tissue				NAT	
B5	11	75	M	Rib	Bone marrow tissue				NAT	
B6	12	75	M	Rib	Bone marrow tissue				NAT	
C1	13	50	F	Rib	Bone marrow tissue				NAT	
C2	14	50	F	Rib	Bone marrow tissue				NAT	
C3	15	62	F	Rib	Bone marrow tissue				NAT	
C4	16	62	F	Rib	Bone marrow tissue				NAT	
C5	17	72	M	Rib	Bone marrow tissue				NAT	
C6	18	72	M	Rib	Bone marrow tissue				NAT	
D1	19	58	M	Rib	Bone marrow tissue				NAT	

D2	20	58	M	Rib	Bone marrow tissue				NAT	
D3	21	60	M	Rib	Bone marrow tissue				NAT	
D4	22	60	M	Rib	Bone marrow tissue				NAT	
D5	23	56	M	Rib	Bone marrow tissue				NAT	
D6	24	56	M	Rib	Bone marrow tissue				NAT	
-	0	42	M	Adrenal gland	Pheochromocytoma (tissue marker)		-		Malignant	



AMSBIO | www.amsbio.com | info@amsbio.com

UK & Rest of the World 
 184 Park Drive, Milton Park
 Abingdon OX14 4SE.
 T: +44 (0) 1235 828 200
 F: +44 (0) 1235 820 482

North America 
 1035 Cambridge Street,
 Cambridge, MA 02141.
 T: +1 (617) 945-5033 or
 T: +1 (800) 987-0985
 F: +1 (617) 945-8218

Europe 
 Berenkoog 41,
 1822 BH Alkmaar,
 Netherlands
 T: +31 (0) 72 8080244
 F: +31 (0) 72 8080142

Switzerland 
 Via Lisano 3,
 (CP.683)
 CH-6900
 T: +41 (0) 91 604 55 22
 F: +41 (0) 91 605 17 85

Germany 
 Bockenheimer Landstr. 17/19
 60325 Frankfurt/Main
 T: +49 (0) 69 779099
 F: +49 (0) 69 13376880