



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)



*Data Sheet***Ready-to-Hybridize Minipig Multiple
Tissue Northern Blot****Catalog #:** NN-MT-1**Quantity:** 1 blot**Storage Conditions:** Store in a sealed bag at room temperature or 4°C. It is good for several months.**Applications:** Minipig multiple tissues Northern blot can be used to:

- Analyze gene expression pattern in 15 multiple tissues simultaneously.
- Determine size and relative abundance of specific messages in the mixed tissues.
- Identify alternative spliced forms.

Quality Control: Every step of preparation of the blot, from harvesting tissues and extraction of RNA to the blotting, is carefully monitored to ensure the superior quality and performance.

- Blot is made from total RNA that is treated with RNase-free DNase to remove residual DNA. The purity and integrity of RNA are tested by formaldehyde-denatured agarose gel electrophoresis.
- The efficiency of transfer is checked by staining gel with ethidium bromide.
- The integrity of the blotted RNA is tested by beta-actin specific probe.

Description: Total RNA (20 µg each) was extracted from freshly dissected multiple Minipig tissues (Gottingen), fractionated through formaldehyde-denatured agarose gel, transferred to a positively charged nylon membrane, and cross-lined by UV light.

The band sizes of the RNA marker (major bands: 9kb, 6 kb, 5 kb, 4 kb, 3kb, 2.5 kb, 2kb, 1.5 kb, 1.0 kb and 0.5 kb) are marked on the membrane.

Tissues on the blot: RNAs were loaded in the following left to right order: (1) Brain (whole), (2) stomach (whole), (3) intestine (whole), (4) colon, (5) liver, (6) lung (7) kidney, (8) Heart, (9) ovary, (10) skeletal muscles, (11) spleen, (12) testis, (13) thymus, (14) prostate, and (15) Pancreas.**Hybridization:** Prior to use, wet blot in water for few minutes. Zyagen blots can be hybridized using isotopic and non-isotopic DNA or RNA probes. You can follow any general protocol of Northern blot analysis for hybridization, post-washing, and stripping (reusing blot) probe. However, we recommend protocols published in the following books: "Short protocols in molecular biology" by Frederick Ausubel et al. and "RNA Methodologies" by Robert Farrell.**Stripping Probe:** Blot is reusable for several times, just strip probe and re-hybridize with a new one. Keep blot wet until the previously hybridized probe has been removed. It is extremely difficult to completely remove hybridized probe from a membrane that has been allowed to dry.**Intended use of the products:** All products sold by AMSBIO are intended for research use only and must not be used for medical diagnostic or drug purposes. AMSBIO products may not resell, modified for resale or used to manufacture commercial products without written approval from AMSBIO.**Limited product warranty:** AMSBIO has a total quality assurance policy. If any product does not meet the specification stated in the product data sheet will be replaced at no charge. This warranty limits AMSBIO liability to the replacement of the product only. No other warranties of any kind, express or implied, including without limitation, implied warranties of merchantability or fitness for a particular purpose, are provided by AMSBIO. AMSBIO shall not liable for any incidental, in direct, special or consequential damages, even if it is aware of the possibility of such damages.**UK & Rest of World**

184 Milton Park, Abingdon
OX14 4SE, Oxon, UK
Tel: +44 (0) 1235 828 200
Fax: +44 (0) 1235 820 482

Switzerland

Centro Nord-Sud 2E
CH-6934 Bioggio-Lugano
Tel: +41 (0) 91 604 55 22
Fax: +41 (0) 91 605 17 85

Deutschland


Bockenheimer Landstr. 17/19
60325 Frankfurt/Main
Tel: +49 (0) 69 779099
Fax: +49 (0) 69 13376880

United States

23591 El Toro Rd, Suite #167
Lake Forest, CA 92630
Tel: +1 800 987 0985
Fax: +1 949 265 7703

amsbio

info@amsbio.com

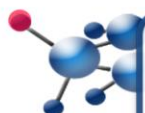

UK & Rest of World

184 Milton Park, Abingdon
OX14 4SE, Oxon, UK
Tel: +44 (0) 1235 828 200
Fax: +44 (0) 1235 820 482

Switzerland

Centro Nord-Sud 2E
CH-6934 Bioggio-Lugano
Tel: +41 (0) 91 604 55 22
Fax: +41 (0) 91 605 17 85

Deutschland

Bockenheimer Landstr. 17/19
60325 Frankfurt/Main
Tel: +49 (0) 69 779099
Fax: +49 (0) 69 13376880

North America

23591 El Toro Rd, Suite #147
Lake Forest, CA 92630
Tel: + 1 800 987 0985
Fax: + 1 949 265 7703

amsbio

info@amsbio.com

www.amsbio.com
AMS Biotechnology