

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



Rabbit IgG-F(ab')2 Fragment cross-adsorbed Antibody

F(ab')2 Goat Polyclonal Conjugate DyLight® 488

Antigen Affinity Purified

Catalog No. A120-212D2 Lot No. A120-212D2-10

APPLICATIONS IHC, ICC, F, IF

SPECIES REACTIVITY Rabbit. Minimum reactivity to human, mouse and rat

AMOUNT 1 ml

CONCENTRATION 0.5 mg/ml

STORAGE/SHELF LIFE 2 – 8° C / 1 year from date of receipt

PHYSICAL STATE Liquid

BUFFER Phosphate Buffered Saline (PBS) containing 0.2% BSA and 0.09% Sodium Azide

FLUOROPHORE/PROTEIN 5.6
ISOTYPE IgG
ORIGIN USA

PRODUCTIONAntiserum was cross adsorbed using human, mouse and rat immunosorbents to remove cross reactive antibodies. The antibody to rabbit IgG–F(ab')2 was isolated by affinity

chromatography using antigen coupled to agarose beads. F(ab')2 fragments were generated using a pepsin digestion. Fc fragments and whole IgG molecules have been removed.

F(ab')2 fragments were conjugated to DyLight® 488.

Antibody concentration was determined by extinction coefficient: absorbance at 280 nm of

1.4 equals 1.0 mg of IgG.

By immunoelectrophoresis and ELISA this antibody reacts specifically with F(ab')2 fragments of rabbit IgG. No antibody was detected against non-immunoglobulin serum proteins. Cross reactivity with IgA and IgM is negligible. Less than 1% cross reactivity to human, mouse and rat F(ab')2 was detected. This antibody may cross react with F(ab')2 fragments of

IgG from other species.

APPLICATIONS Centrifuge tube to remove product from lid. Optimal working dilutions should be determined

experimentally by the investigator. Prepare working dilution immediately before use.

Immunohistochemistry 1:25 – 1:200

Immunocytochemistry 1:25 – 1:200

Flow Cytometry 1:50 - 1:200

Immunofluorescence 1:25 – 1:500

APPLICATION NOTES Not all listed applications have been specifically tested by our laboratory.

DyLight® 488 is excited at 493 (in PBS) and emits at 518 (in PBS).

DyLight[®] is a trademark of Thermo Fisher Scientific Inc. and its subsidiaries.


ADDITIONAL INFO https://www.bethyl.com/product/A120-212D2

Use the link above to view SDS, a current list of citations, and other product specific information.

This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc. Eric McIntush, PhD | Chief Scientific Officer

Date: June 24, 2020



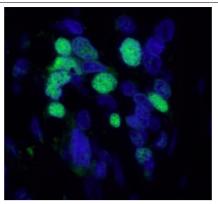

Detection of human PCNA by immunofluorescence. Sample: FFPE section of human ovarian carcinoma. Primary Antibody: Affinity purified rabbit anti-PCNA (Cat. No. IHC-00012 Lot1) used at a dilution of 1:100. Secondary Antibody: Green-fluorescent F(ab')2 goat anti-rabbit IgG-F(ab')2 Fragment cross-adsorbed Antibody DyLight® 488 Conjugated (A120-212D2 Lot 2) used at a dilution of 1:50 (10µg/ml). Counterstain: DAPI (blue)