

# Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



## Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

## Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

#### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



#### Pds5B Antibody



Antigen Affinity Purified Protein ID CAB69911.1

Catalog No. A300-537A-T GeneID 23047

Lot No. A300-537A-T-2

**APPLICATIONS** WB, IP, IHC, IHC-IF

**PRESUMED REACTIVITY** Based on 100% sequence identity, this antibody is predicted to react with Human

**AMOUNT** 10 μl

CONCENTRATION 200 μg/ml

STORAGE/SHELF LIFE 2 - 8°C / 1 year from date of receipt

PHYSICAL STATE Liquid

**BUFFER** Tris-buffered Saline containing 0.1% BSA and 0.09% Sodium Azide

ISOTYPE IgG
ORIGIN USA

**PRODUCTION** Antibody was affinity purified using an epitope specific to Pds5B immobilized on solid

**PROCEDURES** support.

The epitope recognized by A300-537A-T maps to a region between residues 1200 and 1250 of human Precocious Dissociation of Sisters-5B using the numbering given in entry

CAB69911.1 (GeneID 23047).

**APPLICATIONS** Centrifuge tube to remove product from lid. Optimal working dilutions should be determined

experimentally by the investigator. Prepare working dilution immediately before use.

Western Blot 1:2,000 - 1:10,000Immunoprecipitation  $2 - 10 \mu g/mg$  lysate

Immunohistochemistry 1:100 – 1:500. Epitope retrieval with citrate buffer pH 6.0 is

recommended for FFPE tissue sections.

Immunofluorescence 1:100 - 1:500. Epitope retrieval with citrate buffer pH 6.0 is

(IHC) recommended for FFPE tissue sections.

ADDITIONAL INFO https://www.fortislife.com/p/A300-537A-T


Use the link above to view SDS, a current list of citations, and other product specific information.


This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc.

Michael Spencer, PhD

Date: August 8, 2022

Phone: 800.338.9579 • Fax: 866.597.6105 • Web: www.fortislife.com Orders: orders@fortislife.com • Support: technical@fortislife.com Pds5B Antibody A300-537A-T

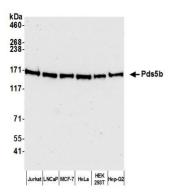

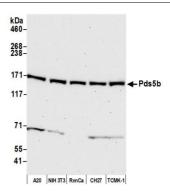

Detection of human Pds5b by western blot. Samples: Whole cell lysate (50 μg) from Jurkat, LNCaP, MCF-7, HeLa, HEK293T, and Hep-G2 cells prepared using NETN lysis buffer. Antibody: Affinity purified rabbit anti-Pds5b antibody A300-537A-T (lot A300-537A-T-2) used for WB at 1:1000. Detection: Chemiluminescence with an exposure time of 30 seconds.

Detection of mouse Pds5b by western blot. Samples: Whole cell lysate (50 μg) from A20, NIH 3T3, RenCa, CH27, and TCMK-1 cells prepared using NETN lysis buffer. Antibody: Affinity purified rabbit anti-Pds5b antibody A300-537A-T (lot A300-537A-T-2) used for WB at 1:1000. Detection: Chemiluminescence with an exposure time of 30 seconds.