

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



BTF Antibody

Antigen Affinity Purified



Protein ID NP 055554.1

Catalog No. A300-609A-T GeneID 9774

Lot No. A300-609A-T-1

APPLICATIONS WB, IP, IHC
SPECIES REACTIVITY Human

AMOUNT 10 μl

CONCENTRATION 200 μg/ml

STORAGE/SHELF LIFE 2 - 8°C / 1 year from date of receipt

PHYSICAL STATE Liquid

BUFFER Tris-buffered Saline containing 0.1% BSA and 0.09% Sodium Azide

ISOTYPE IgG
ORIGIN USA

PRODUCTION PROCEDURES

Antibody was affinity purified using an epitope specific to BTF immobilized on solid support.

The epitope recognized by A300-609A-T maps to a region between residues 300 and 350 of human BCL2-Associated Transcription Factor 1 using the numbering given in entry

NP_055554.1 (GeneID 9774).

APPLICATIONS Centrifuge tube to remove product from lid. Optimal working dilutions should be determined

experimentally by the investigator. Prepare working dilution immediately before use.

Western Blot 1:2,000 - 1:10,000Immunoprecipitation $1 - 4 \mu g/mg$ lysate

Immunohistochemistry 1:100 - 1:500. Epitope retrieval with citrate buffer pH 6.0 is

recommended for FFPE tissue sections.

ADDITIONAL INFO https://www.bethyl.com/product/A300-609A-T

Use the link above to view SDS, a current list of citations, and other product specific information.


IP-western blot protocol: https://www.bethyl.com/content/protocol_IP_WB

This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc.


Michael Spencer, PhD

Date: June 6, 2022

Phone: 800.338.9579 • Fax: 866.597.6105 • Web: www.bethyl.com Orders: orders@fortislife.com • Support: technical@fortislife.com BTF Antibody A300-609A-T

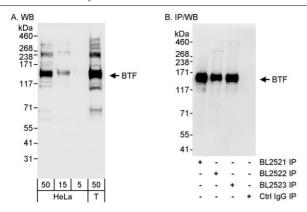

Detection of human BTF by western blot and immunoprecipitation. Samples: Whole cell lysate from HeLa (5, 15 and 50 μ g for WB; 1 mg for IP, 20% of IP loaded) and HEK293T (T; 50 μ g) cells. Antibodies: Affinity purified rabbit anti-BTF antibody BL2522 (Cat. No. A300–609A) used for WB at 0.04 μ g/ml (A) and 1 μ g/ml (B) and used for IP at 3 μ g/mg lysate (B). BTF was also immunoprecipitated using rabbit anti-BTF antibodies BL2521 (Cat. No. A300–608A) and BL2523 (Cat. No. A300–610A) at 3 μ g/mg lysate. Detection: Chemiluminescence with exposure times of 3 seconds (A and B).

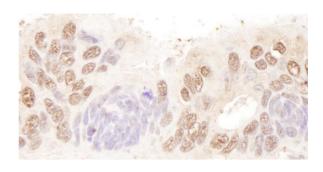

Detection of human BTF by immunohistochemistry. Sample: FFPE section of human ovarian carcinoma. Antibody: Affinity purified rabbit anti-BTF (Cat. No. A300-609A Lot1) used at a dilution of 1:200 (1µg/ml). Detection: DAB