



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)



NF-kappaB1 Antibody

Rabbit Polyclonal

Antigen Affinity Purified

Protein ID AAF35232.1

Catalog No. A301-819A-T

GeneID 4790

Lot No. A301-819A-T-2

APPLICATIONS	IP
SPECIES REACTIVITY	Human
AMOUNT	10 µl
CONCENTRATION	1000 µg/ml
STORAGE/SHELF LIFE	2 – 8°C / 1 year from date of receipt
PHYSICAL STATE	Liquid
BUFFER	Tris-citrate/phosphate buffer, pH 7 to 8 containing 0.09% Sodium Azide
ISOTYPE	IgG
ORIGIN	USA
PRODUCTION PROCEDURES	Antibody was affinity purified using an epitope specific to NF-kappaB1 immobilized on solid support.

The epitope recognized by A301-819A-T maps to a region between residue 1 and 50 of human nuclear factor of kappa light polypeptide gene enhancer in B-cells 1 (p105) using the numbering given in entry AAF35232.1 (GeneID 4790).

Immunoglobulin concentration was determined by extinction coefficient: absorbance at 280 nm of 1.4 equals 1.0 mg of IgG.

APPLICATIONS Centrifuge tube to remove product from lid. Optimal working dilutions should be determined experimentally by the investigator. Prepare working dilution immediately before use.


Western Blot Not recommended. Use rabbit anti-NF-kappaB1 antibody A301-820A.

Immunoprecipitation 2 – 10 µg/mg lysate

ADDITIONAL INFO <https://www.bethyl.com/product/A301-819A-T>

Use the link above to view SDS, a current list of citations, and other product specific information.
IP-western blot protocol: https://www.bethyl.com/content/protocol_IP_WB

This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc.
Michael Spencer, PhD Date: June 6, 2022

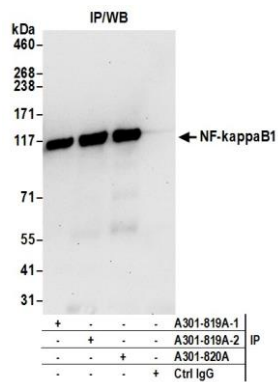

Detection of human NF-kappaB1 by western blot of immunoprecipitates. *Samples:* Whole cell lysate (0.5 or 1.0 mg per IP reaction; 20% of IP loaded) from Jurkat cells prepared using NETN lysis buffer. *Antibodies:* Affinity purified rabbit anti-NF-kappaB1 antibody A301-819A (lot A301-819A-2) used for IP at 6 µg per reaction. NF-kappaB1 was also immunoprecipitated by a previous lot of this antibody (lot A301-819A-1) and rabbit anti-NF-kappaB1 antibody A301-820A. For blotting immunoprecipitated NF-kappaB1, A301-820A was used at 1 µg/ml. *Detection:* Chemiluminescence with an exposure time of 30 seconds.