



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Phospho ZAP70 (Y493) Antibody

Rabbit Polyclonal

Antigen Affinity Purified Protein ID NP_001070.2

Catalog No. A302-008A-T GenelD 7535

Lot No. A302-008A-T-1

APPLICATIONS	IP
SPECIES REACTIVITY	Human
PRESUMED REACTIVITY	Based on 100% sequence identity, this antibody is predicted to react with Mouse
AMOUNT	10 µl
CONCENTRATION	1000 µg/ml
STORAGE/SHELF LIFE	2 – 8°C / 1 year from date of receipt
PHYSICAL STATE	Liquid
BUFFER	Tris-citrate/phosphate buffer, pH 7 to 8 containing 0.09% Sodium Azide
ISOTYPE	IgG
ORIGIN	USA
PRODUCTION PROCEDURES	Antibody was affinity purified using an epitope specific to Phospho ZAP70 (Y493) immobilized on solid support.

The epitope recognized by A302-008A-T maps to a region including tyrosine 493 (when it is phosphorylated) of human zeta-chain (TCR) associated protein kinase 70kDa using the numbering given in entry NP_001070.2 (GenelD 7535).

Immunoglobulin concentration was determined by extinction coefficient: absorbance at 280 nm of 1.4 equals 1.0 mg of IgG.

APPLICATIONS Centrifuge tube to remove product from lid. Optimal working dilutions should be determined experimentally by the investigator. Prepare working dilution immediately before use.

Western Blot Not recommended. A302-008A has not performed satisfactorily when used for WB of phospho ZAP70 in crude preparations (e.g. whole cell lysate). This antibody can be used for WB of enriched (e.g. immunoprecipitated) sources of ZAP70.

Immunoprecipitation 2 – 5 µg/mg lysate

ADDITIONAL INFO <https://www.bethyl.com/product/A302-008A-T>

Use the link above to view SDS, a current list of citations, and other product specific information.


IP-western blot protocol: https://www.bethyl.com/content/protocol_IP_WB

This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc.

Michael Spencer, PhD


Date: June 6, 2022

Immunoprecipitate ZAP70 - Blot phospho ZAP70-Y493

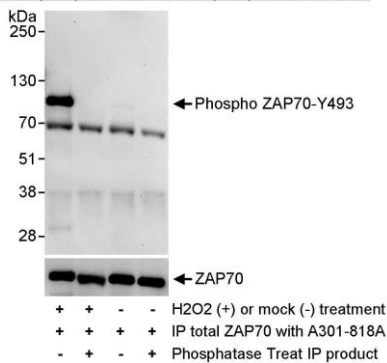

Detection of Phosphorylation of human ZAP70 on Y493 by western blot of immunoprecipitates. *Samples:* Whole cell lysate (1 mg for IP, 20% of IP loaded) from Jurkat cells that had been treated with hydrogen peroxide (+) or mock treated (-). *Antibodies:* Affinity purified rabbit anti-ZAP70 antibody A301-818A was used at 3 µg/mg lysate to immunoprecipitate ZAP70. The immunoprecipitates were treated with phosphatase (+) or mock (-) treated. For blotting immunoprecipitated Phospho ZAP70, anti-Phospho ZAP70 (Y493) antibody A302-008A was used at 1 µg/ml. *Detection:* Chemiluminescence with an exposure time of 30 seconds.

IP Phospho ZAP70 - Blot Zap70

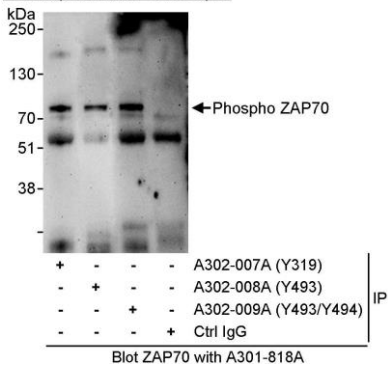

Detection of Phosphorylation of human ZAP70 on Y319, Y493 or Y493/Y494 by western blot of immunoprecipitates. *Samples:* Whole cell lysate (1 mg for IP, 20% of IP loaded) from Jurkat cells that had been treated with hydrogen peroxide. *Antibodies:* Affinity purified rabbit anti-Phospho ZAP70 antibodies A302-007A (for phosphorylation at Y319), A302-008A (for phosphorylation at Y493), and A302-009A (for phosphorylation at Y493/Y494) used for IP at 3 µg/mg lysate. For blotting immunoprecipitated Phospho ZAP70, antibody A301-818A (for total ZAP70) was used at 1 µg/ml. *Detection:* Chemiluminescence with an exposure time of 3 minutes.