



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

**Description**

Cas9 (Streptococcus pyogenes CRISPR associated protein 9) is an endonuclease enzyme that, when recruited to a specific DNA sequence by the sgRNA (single guide RNA), introduces a double stranded break into the DNA. This double stranded break is repaired in mammalian cells either through Non-Homologous End Joining or Homologous Recombination. Non-Homologous End Joining often results in the deletion or insertion of several base pairs at the cut site, which, when resulting in a frameshift, causes the functional inactivation of the targeted gene.

Cas9 Lentivirus can be used to generate Cas9 expressing cells in almost any mammalian cell line. Cells stably expressing Cas9 can then be transduced or electroporated with sgRNA targeting a gene of interest to quickly generate knock-out cell pools or cell lines.

The Cas9 Lentiviruses are replication incompetent, HIV-based VSV-G pseudo-typed lentiviral particles that are ready to be transduced into almost all types of mammalian cells, including primary and non-dividing cells. The particles contain a Cas9 gene driven by an EF1a promoter, along with a Hygromycin selection marker.

**Application**

Transient expression of Cas9 in target cells.

Generation of a stable Cas9 over-expressing cell line following hygromycin selection.

**Formulation**

The lentiviruses were produced from HEK293T cells in medium containing 90% DMEM + 10% FBS.

**Titer**

Two vials (500  $\mu$ l x 2) of lentivirus at a titer  $\geq 5 \times 10^6$  TU/ml. The titer will vary with each lot; the exact value is provided with each shipment.

**Storage**

Lentiviruses are shipped with dry ice. For long term storage, it is recommended to store the lentiviruses at  $-80^{\circ}\text{C}$ . Avoid repeated freeze-thaw cycles. Titers can drop significantly with each freeze-thaw cycle.

**Biosafety**

None of the HIV genes (gag, pol, rev) will be expressed in the transduced cells. Although the pseudotyped lentiviruses are replication-incompetent, they do require the use of a Biosafety Level 2 facility. BPS recommends following all federal, state, local, and institutional regulations and using all appropriate safety precautions.

**Materials Required but Not Supplied**



These materials are not supplied with this lentivirus but are necessary to follow the designed protocol. BPS Bioscience media, reagents, and luciferase assay systems are all validated and optimized for use with this lentivirus and are highly recommended for best results.

Name	Ordering Information
Negative Control Luciferase Lentivirus	BPS Bioscience #79578
ONE-STEP Luciferase Assay System	BPS Bioscience #60690

**License Disclosure**

Visit [bpsbioscience.com/license](https://bpsbioscience.com/license) for the label license and other key information about this product.

**Troubleshooting Guide**

For all questions, please email [support@bpsbioscience.com](mailto:support@bpsbioscience.com).

**Figures and Validation Data**

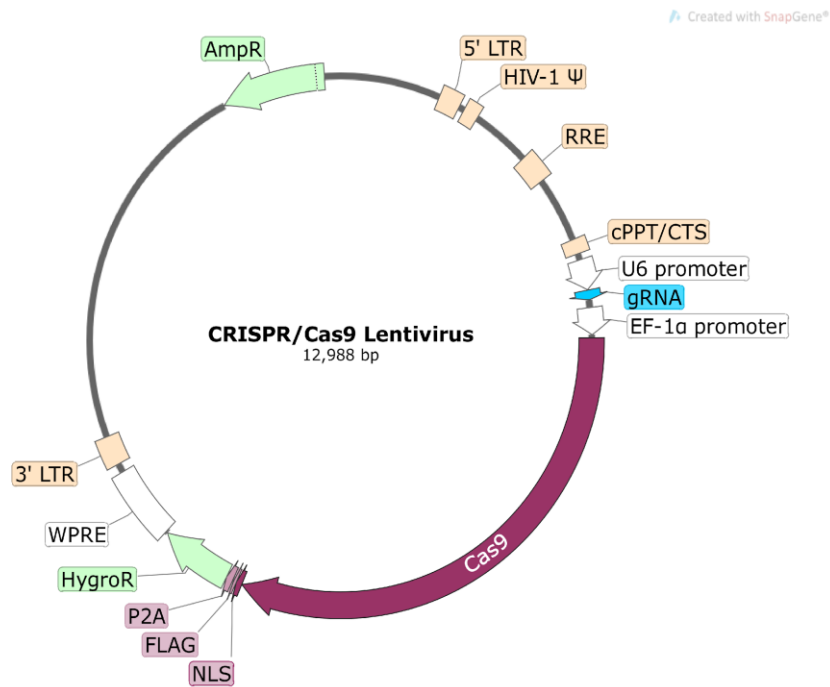


Figure 1. Schematic of the lenti-vector used to generate the firefly luciferase lentivirus

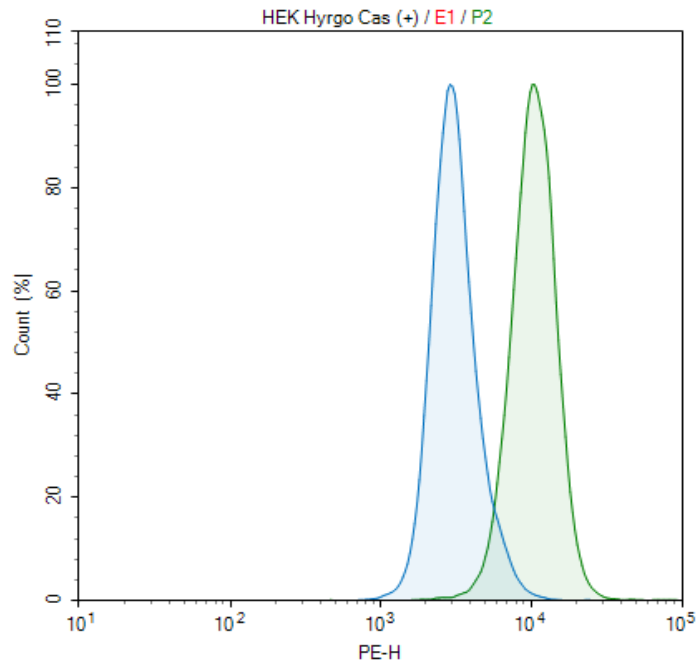


Figure 2. Transduction of Cas9 Lentivirus in HEK293 cells.

HEK293 cells were transduced with 5,000,000 TU/well of Cas9 lentivirus. 72 hours after transduction, cells were stained with PE anti-FLAG antibody (BioLegend, #637309) and analyzed by FACS. Parental HEK293 cells are shown in blue, and the transduced cells are shown in green.

## Related Products

Products	Catalog #	Size
Cas9 Lentivirus (Puromycin)	78066	500 µl x 2
Cas9, His-tag ( <i>S. pyogenes</i> )	100206-1	50 µg
TCR CRISPR/Cas9 Lentivirus (Integrating)	78055	500 µl x 2
TCR CRISPR/Cas9 Lentivirus (Non-Integrating)	78062	500 µl x 2
PD-1 CRISPR/Cas9 Lentivirus (Integrating)	78052	500 µl x 2
PD-1 CRISPR/Cas9 Lentivirus (Non-Integrating)	78059	500 µl x 2
Cas9 Expressing Jurkat cells	78070	2 vials
Cas9 Expressing MDA-MB-231 cells	78069	2 vials
Cas9 Expressing A549 cells	78072	2 vials
Cas9 Expressing HCT116 cells	78073	2 vials
Cas9 Expressing Raji cell	78071	2 vials