



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Description

Adeno-Associated Virus Serotype 5 (AAV5) differs from other parvovirus serotypes according to serological and DNA hybridization data, and AAV5 also uses different inverted terminal repeats (ITRs) compared to other AAV serotypes. AAV5 is the most efficient vector for transducing sensory neurons and is good at mediating gene transfer into human and murine airway epithelia. In addition, AAV5 vectors show a higher tropism for both mouse and human dendritic cells than AAV1, AAV2, AAV7, or AAV8.

These AAV particles constitutively express the firefly (*Photinus pyralis*) luciferase and mCherry genes connected via a T2A linker, under the control of a CMV promoter. The T2A self-cleaving peptide (derived from *Thosea asigna* virus 2A) leads to the efficient cleavage of the transcript, and expression of luciferase and mCherry as two separate proteins.

Application(s)

- Use as a positive control for transduction
- Optimize transduction assays and track protein expression over time

Serotype

Wild-type AAV Serotype 5

Formulation

AAV5 was produced in HEK293-AAV cells and is supplied in PBS-MK (PBS Magnesium-Potassium) buffer containing 0.01% Pluronic F68.

Purification

The purity of the AAV particles was confirmed to be greater than 90% by staining with One-Step Lumitein™ UV Protein Gel Stain (Biotium #21005-1L). Purity will vary with each lot; the exact value will be provided with each shipment.

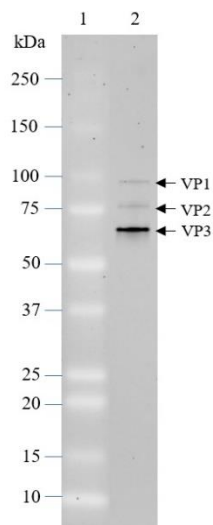


Figure 1. Purified AAV5 Luciferase-mCherry particles.

Staining of a 4-20% SDS-PAGE gel. The protein ladder is in lane 1, and 2×10^9 GC (genome copy number) of AAV5 is shown in lane 2. AAV viral proteins VP1, VP2, and VP3 are labeled.

Titer

Two vials (50 µl x 2) of AAV at a titer $\geq 1 \times 10^{12}$ TU/ml. The titer is determined by qPCR and will vary with each lot; the exact value is provided with each shipment.

Storage



AAV is shipped with dry ice. For long-term storage, it is recommended to store AAV at -80°C. Avoid repeated freeze-thaw cycles. Titters can drop significantly with each freeze-thaw cycle.

Biosafety



Recombinant AAV is inherently replication-deficient and not known to cause any human diseases. Additionally, following transduction, AAV vectors exist episomally and do not integrate into or disrupt the host cell's genome. AAV requires the use of a Biosafety Level 1 facility. BPS Bioscience recommends following all local, federal, state, and institutional regulations and using all appropriate safety precautions.

Validation Data

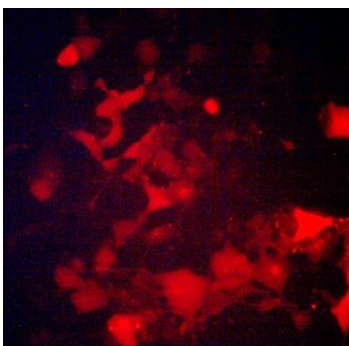


Figure 2. Transduction of HEK293 cells using AAV5 Luciferase-mCherry particles. 1×10^5 cells/well were transduced in a 6-well plate with AAV5 Luciferase-mCherry at an MOI of 2×10^4 . After 72 hours of transduction, mCherry expression in the target cells was observed under a fluorescence microscope. mCherry expression was stable over time and still observed 30 days after transduction.

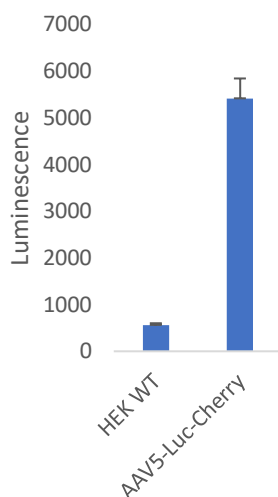


Figure 3. Luciferase activity of HEK293 cells transduced by AAV5 Luciferase-mCherry particles. 1×10^5 cells/well were transduced in a 6-well plate with AAV5 Luciferase-mCherry at an MOI of 2×10^4 . After 72 hours of transduction, transduced cells or parental HEK293 cells were seeded in a 96-well plate at a density of 2×10^4 cells/well, and luciferase activity was measured using the ONE-Step™ Luciferase Assay System (BPS Bioscience #60690).

Troubleshooting Guide

Visit bpsbioscience.com/lentivirus-faq for detailed troubleshooting instructions. For all further questions, please email support@bpsbioscience.com.

Related Products

<i>Products</i>	<i>Catalog #</i>	<i>Size</i>
AAV2 ZsGreen	78444	50 μ l x 2
AAV5 ZsGreen	78447	50 μ l x 2
AAV2 Luciferase-eGFP	78462	50 μ l x 2
AAV3 Luciferase-eGFP	78463	50 μ l x 2
AAV5 Luciferase-eGFP	78465	50 μ l x 2
AAV9 Luciferase-eGFP	78468	50 μ l x 2
AAV6 Luciferase-mCherry	78475	50 μ l x 2
AAV8 Luciferase-mCherry	78476	50 μ l x 2