



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

### Description

Adeno-Associated Virus Serotype 5 (AAV5) differs from other parvovirus serotypes according to serological and DNA hybridization data, and AAV5 also uses different inverted terminal repeats (ITRs) compared to other AAV serotypes. AAV5 is the most efficient vector for transducing sensory neurons and is good at mediating gene transfer into human and murine airway epithelia. In addition, AAV5 vectors show a higher tropism for both mouse and human dendritic cells than AAV1, AAV2, AAV7, or AAV8.

Cas9 is an endonuclease enzyme that introduces a double stranded break into the DNA when recruited to a specific DNA sequence by the sgRNA (single guide RNA). This double stranded break is repaired in mammalian cells either through Non-Homologous End Joining or Homologous Recombination. Non-Homologous End Joining often results in the deletion or insertion of several base pairs at the cut site, which, when resulting in a frameshift, causes the functional inactivation of the gene.

SaCas9 (*Staphylococcus aureus* CRISPR-associated protein 9) has demonstrated high cutting efficiency in mammalian cells, and its smaller size makes it ideal for packaging into AAV. SaCas9 recognizes a longer protospacer adjacent motif (PAM) site, 5'-NNGRRT-3', than the more traditional SpCas9 (*Streptococcus pyogenes* CRISPR-associated protein 9). These AAV particles constitutively express SaCas9 under the control of a CMV promoter.

### Application(s)

Express SaCas9 to enable genome editing

### Serotype

Wild-type AAV Serotype 5

### Formulation

AAV was produced in HEK293-AAV cells and is supplied in PBS-MK (PBS Magnesium-Potassium) buffer containing 0.01% Pluronic F68.

### Purification

The purity of the AAV particles was confirmed to be greater than 90% by staining with One-Step Lumitein™ UV Protein Gel Stain (Biotium #21005-1L). The purity will vary with each lot; the exact value will be provided with each shipment.

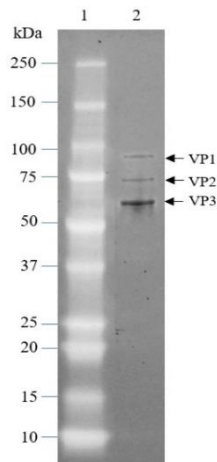


Figure 1. Purified AAV5 SaCas9 particles.

Staining of a 4-20% SDS-PAGE gel. The protein ladder is in lane 1, and  $5 \times 10^9$  GC (genome copy number) of AAV5 is shown in lane 2. AAV viral proteins VP1, VP2, and VP3 are labelled.

**Titer**


Two vials (50  $\mu$ l x 2) of AAV at a titer  $\geq 1 \times 10^{12}$  TU/ml. The titer is determined by qPCR and will vary with each lot; the exact value will be provided with each shipment.

**Storage**

AAV is shipped with dry ice. For long-term storage, it is recommended to store AAV at  $-80^{\circ}\text{C}$ . Avoid repeated freeze-thaw cycles. Titers can drop significantly with each freeze-thaw cycle.

**Biosafety**

Recombinant AAV is inherently replication-deficient and not known to cause any human diseases. Additionally, following transduction, AAV vectors exist episomally and do not integrate into or disrupt the host cell's genome. AAV requires the use of a Biosafety Level 1 facility. BPS Bioscience recommends following all local, federal, state, and institutional regulations and using all appropriate safety precautions.

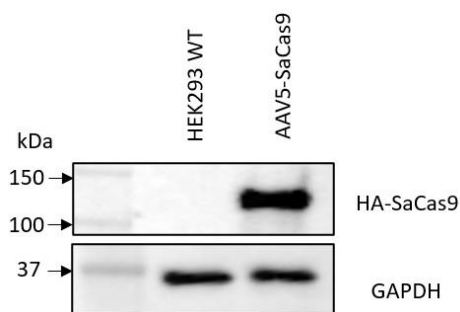**Validation Data**

*Figure 2. Transduction of HEK293 cells using AAV5 SaCas9.*

$1 \times 10^5$  cells were transduced in a 24-well plate with AAV5 SaCas9 particles at an MOI of  $1 \times 10^4$ . After 72 hours of transduction, SaCas9 expression in the target cells was measured by Western blot (Anti-HA.11 Epitope Tag Antibody, Biolegend #901516).

**Troubleshooting Guide**

Visit [bpsbioscience.com/lentivirus-faq](https://bpsbioscience.com/lentivirus-faq) for detailed troubleshooting instructions. For all further questions, please email [support@bpsbioscience.com](mailto:support@bpsbioscience.com).

**Related Products**

Products	Catalog #	Size
AAV3 ZsGreen	78445	50 $\mu$ l x 2
AAV5 ZsGreen	78447	50 $\mu$ l x 2
AAV2 Luciferase	78453	50 $\mu$ l x 2
AAV6 Luciferase-eGFP	78466	50 $\mu$ l x 2
AAV5 Luciferase-eGFP	78465	50 $\mu$ l x 2
AAV9 Luciferase-eGFP	78468	50 $\mu$ l x 2
AAV2 Luciferase-mCherry	78462	50 $\mu$ l x 2
AAV5 Luciferase-mCherry	78474	50 $\mu$ l x 2
AAV8 Luciferase-mCherry	78476	50 $\mu$ l x 2
AAV1 SaCas9	78479	50 $\mu$ l x 2