



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)



Description

The Fibroblast growth factor receptor 4 (FGFR4) Kinase Assay Kit is designed to measure FGFR4 kinase activity for screening and profiling applications using ADP-Glo® as a detection reagent. The assay kit comes in a convenient 96-well format, with enough purified recombinant FGFR4 kinase, kinase substrate, ATP and kinase assay buffer for 100 enzyme reactions.

Background

FGFR4 (also known as CD334) is a transmembrane receptor with tyrosine kinase activity. The extracellular portion of the receptor contains three immunoglobulin-like domains and is responsible for binding fibroblast growth factor 1 (FGF1). It contains a single transmembrane domain followed by the cytosolic kinase domain. FGFR4 may be involved in ovarian cancer as well as adrenocortical carcinoma and is a promising therapeutic target in oncology applications.

Applications

Study enzyme kinetics and screen small molecular inhibitors for drug discovery and high throughput (HTS) applications.

Supplied Materials

Catalog #	Name	Amount	Storage
40213	FGFR4*	20 µg	-80°C
79334	Kinase Assay Buffer 1 (5x)	1.5 ml	-20°C
79686	ATP (500 µM)	50 µl	-20°C
40217	PTK Substrate (Poly-Glu,Tyr 4:1) (10 mg/ml)	50 µl	-20°C
79696	White 96-well plate	1	Room Temperature

*The concentration of the protein is lot-specific and will be indicated on the tube

Materials Required but Not Supplied

Name	Ordering Information
ADP-Glo® Kinase Assay	Promega #V6930
DTT (Dithiothreitol), 1M, optional	
Microplate reader capable of reading luminescence	
Adjustable micropipettor and sterile tips	
30°C incubator	

Storage Conditions



This assay kit will perform optimally for up to 6 months from date of receipt when the materials are stored as directed.

Safety



This product is for research purposes only and not for human or therapeutic use. This product should be considered hazardous and is harmful by inhalation, in contact with skin, eyes, clothing, and if swallowed. If contact occurs, wash thoroughly.

Assay Principle

The **ADP-Glo™ Kinase Assay (Promega #V6930)** quantifies the amount of ADP produced by a kinase upon phosphorylation of a substrate. First, addition of the ADP-Glo™ reagent terminates the reaction and quenches the remaining ATP. Second, addition of the Kinase Detection reagent converts the produced ADP to ATP. The new ATP is quantified by a luciferase reaction. The luminescent signal correlates with the amount of ADP generated by the kinase and is linear to 1 mM ATP.

Contraindications

The final concentration of DMSO in the assay should not exceed 1%.

Assay Protocol

All samples and controls should be tested in duplicate.

1. Thaw **5x Kinase Assay Buffer 1**, (**500 µM**) **ATP**, and **PTK Substrate (Poly-Glu,Tyr 4:1)** (**10 mg/ml**).

*Optional: If desired, add DTT to **5x Kinase Assay Buffer 1** to make a 10 mM DTT concentration (for example, add 10 µl of 1 M DTT to 1 ml of **5x Kinase Assay Buffer 1**).*

2. Prepare 3 ml of **1x Kinase Assay Buffer 1** by mixing 600 µl of **5x Kinase Assay Buffer 1** with 2,400 µl water.

*Note: Three (3 ml) of **1x Kinase Assay Buffer 1** is sufficient for 100 reactions.*

3. Prepare the **Master Mix** (12.5 µl/well): N wells x (6 µl of **5x Kinase Assay Buffer 1** + 0.5 µl of **ATP (500 µM)** + 0.5 µl of **PTK Substrate (Poly-Glu,Tyr 4:1)** (10 mg/ml) + 5.5 µl of distilled water. Add 12.5 µl to every well.
4. Prepare the **Test Inhibitor** (2.5 µl/well): for a titration, prepare serial dilutions at concentrations 10-fold higher than the desired final concentrations. The final volume of the reaction is 25 µl.

If the Test Inhibitor is water-soluble:

- 4.1 Prepare serial dilutions in the **1x Kinase Assay Buffer 1**, 10-fold more concentrated than the desired final concentrations.
- 4.2 For the positive and negative controls, use **1x Kinase Assay Buffer 1** (Diluent Solution).

Or

If the Test inhibitor is soluble in DMSO:

- 4.1 Prepare the test inhibitor at 100-fold the highest desired concentration in DMSO, then dilute the inhibitor 10-fold in **1x Kinase Assay Buffer 1** to prepare the highest concentration of the 10-fold intermediate dilutions. The concentration of DMSO is now 10%.
- 4.2 Prepare serial dilutions of the Test Inhibitor at 10-fold the desired final concentrations using 10% DMSO in **1x Kinase Assay Buffer 1** to keep the concentration of DMSO constant.
- 4.3 For positive and negative controls, prepare 10% DMSO in water (vol/vol) so that all wells contain the same amount of DMSO (Diluent Solution).

Note: The final concentration of DMSO should not exceed 1%.

5. Add 2.5 µl of **Test Inhibitor** to each well labeled "Test Inhibitor." For the "Positive Control" and "Blank," add 2.5 µl of **Diluent Solution** (either kinase assay buffer or 10% DMSO in kinase assay buffer, as described above).

6. To the wells designated as "Blank," add 10 μ l of **1x Kinase Assay Buffer 1**.
7. Thaw **FGFR4 kinase** on ice. Briefly spin the tube to recover its full contents. Dilute the protein kinase to 20 ng/ μ l using **1x Kinase Assay Buffer 1** (you need 10 μ l/well).

Note: The concentration of protein is lot-specific and is indicated on the tube. Verify the initial concentration and dilute accordingly.

Note: This kinase is particularly sensitive to freeze/thaw cycles. Avoid multiple freeze/thaw cycles. Do not re-use the thawed protein and do not re-use the diluted kinase.

8. Initiate the reaction by adding 10 μ l of diluted Kinase to the wells designated "Positive Control" and "Test Inhibitor."

Component	Blank	Positive Control	Test Inhibitor
Master Mix	12.5 μ l	12.5 μ l	12.5 μ l
Test Inhibitor	-	-	2.5 μ l
Diluent Solution	2.5 μ l	2.5 μ l	-
1x Kinase Assay Buffer 1	10 μ l	-	-
FGFR4 (20 ng/ μ l)	-	10 μ l	10 μ l
Total	25 μl	25 μl	25 μl

9. Incubate at 30°C for 45 minutes.
10. During the incubation, thaw the ADP-Glo™ reagent. At the end of the 45-minute reaction, add 25 μ l of ADP-Glo™ reagent to each well. Cover the plate with aluminum foil and incubate at room temperature for 45 minutes.
11. Thaw the Kinase Detection Reagent. At the end of the 45-minute incubation, add 50 μ l of Kinase Detection reagent to each well. Cover the plate with aluminum foil and incubate at room temperature for another 45 minutes.
12. Immediately read in a luminometer or a microplate reader capable of reading luminescence. The "Blank" value is subtracted from all other readings.

Reading Luminescence

Luminescence is the emission of light resulting from a chemical reaction. The detection of luminescence requires no wavelength selection because the method used is emission photometry and not emission spectrophotometry.

To properly read luminescence, make sure the plate reader is set for LUMINESCENCE mode. Typical integration time is 1 second, delay after plate movement is 100 msec. Do not use a filter when measuring light emission. Typical settings for the Synergy 2 BioTek plate reader: use the "hole" position on the filter wheel; Optics position: Top; Read type: endpoint. Sensitivity may be adjusted based on the luminescence of a control assay without enzyme (typically we set this value as 100).

Example Results

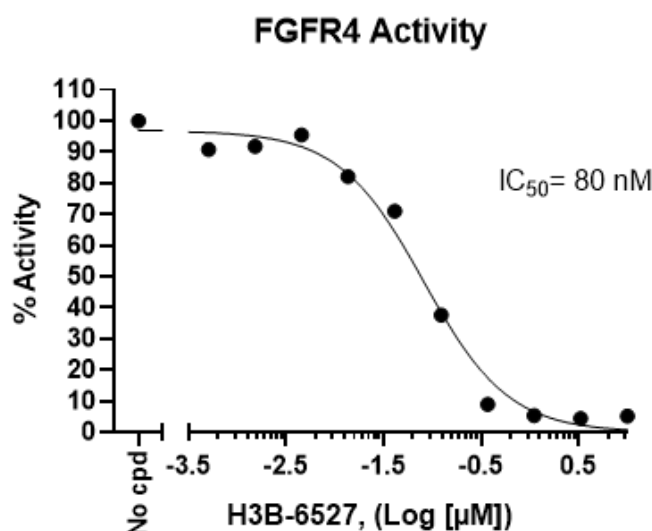


Figure 1: Inhibition of FGFR4 kinase Activity by H3B-6527 (SelleckChem #S8675).

The inhibition of FGFR4 kinase activity was measured in the presence of increasing inhibitor concentrations using the FGFR4 Kinase Assay Kit (BPS Bioscience, #XXXX). The "Blank" value was subtracted from all other values. Results are expressed as the percent of control (kinase activity in the absence of inhibitor, set at 100%).

Data is representative. For lot-specific information, please contact BPS Bioscience, Inc. at support@bpsbioscience.com.

Reference(s)

Tamburello M, *et al.* FGF/FGFR signaling in adrenocortical development and tumorigenesis: novel potential therapeutic targets in adrenocortical carcinoma. *Endocrine*. 2022; 77(3): 411-418.

Troubleshooting Guide

Visit bpsbioscience.com/assay-kits-faq for detailed troubleshooting instructions. For all further questions, please email support@bpsbioscience.com

Related Products

Products	Catalog #	Size
FGFR4, GST-tag Recombinant	40213	10 μg
FGFR3 (V443L) Assay Kit	79819	96 reactions
FGFR3 (L496V) Assay Kit	79818	96 reactions
FGFR2 Assay Kit	79804	96 reactions