



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)



6042 Cornerstone Court W, Ste B
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Data Sheet

CD137/NF- κ B Reporter - HEK293 Recombinant Cell Line

Catalog # 79289

Background

Human CD137 (4-1BB; TNFRS9) is an inducible co-stimulatory molecule that activates T cells. CD137:CD137L-mediated signaling has been shown to be important for proliferation, effector functions and survivals of T cells. CD137 is also expressed in NK and NKT cells. Accumulating evidence shows a role for CD137:CD137L signaling in inflammation, suggesting that inhibition of this pathway may provide a therapeutic avenue to treat autoimmune and inflammatory diseases. Similarly, antibodies targeting CD137 activation in immune cells have demonstrated potent anti-tumor properties in cancer patients.

Description

Recombinant HEK293 cell line expressing a full length human CD137 (NM_001561). The NF- κ B luciferase reporter construct is stably integrated into the genome. The firefly luciferase gene is controlled by 4 copies of NF- κ B response element located upstream of the TATA promoter. Following activation by human CD137 ligand, NF- κ B transcription factors bind to the DNA response elements to induce transcription of the luciferase gene.

Applications

- Screen for activators or inhibitors of CD137 signaling in a cellular context
- Characterize the biological activity of CD137 and its interactions with ligands and antibodies

Format

Each vial contains ~ 2 x 10⁶ cells in 1 ml of 10% DMSO in FBS.

Mycoplasma Testing

This cell line has been screened using the MycoAlert™ Mycoplasma Detection Kit (Lonza, #LT07-118) to confirm the absence of Mycoplasma contamination.

Storage

Store in liquid nitrogen immediately upon receipt.

Culture Medium:

Thaw Medium 1 (BPS Bioscience, #60187): MEM medium (Hyclone #SH30024.01) supplemented with 10% FBS (Thermo Fisher, Cat. #26140079), 1% non-essential amino acids (Hyclone #SH30238.01), 1 mM Na pyruvate (Hyclone #SH30239.01), 1% Penicillin/Streptomycin (Hyclone SV30010.01).

Growth Medium 1G (BPS Bioscience, #79544): Thaw Medium 1, 400 μ g/ml of Geneticin (Thermo Fisher, Cat. #11811031) and 50 μ g/ml of Hygromycin B (Thermo Fisher, Cat. #10687010).

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com



6042 Cornerstone Court W, Ste B
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Recommended Culture Condition

Frozen Cells: Prepare a 50 ml conical tube with 10 ml of pre-warmed Thaw Medium 1 (no G418 or hygromycin). Quickly thaw cells in a 37°C water bath with constant and slow agitation. Clean the outside of the vial with 70% ethanol and immediately transfer the entire content to Thaw Medium 1 (no G418 or hygromycin). Avoid pipetting up and down, and gently rock the conical tube. Spin the cells down at 150 x g for 5 minutes. Discard the medium and re-suspend the cell pellet in fresh Thaw Medium 1 (no G418 or hygromycin). Transfer the entire content to a T25 flask to distribute the cells. Incubate the cells in a humidified 37°C incubator with 5% CO₂. After 24 hours of incubation, change to fresh Thaw Medium 1 (no G418 or hygromycin), being careful to not disturb the attached cells. Begin adding G418 and Hygromycin B to Thaw Medium 1 (Growth Medium 1G) after the first passage.

Subculture: When cells reach 90% confluency, remove the medium and gently wash once with PBS (without Magnesium or Calcium). Treat cells with 0.5 ml of 0.05% trypsin/EDTA and incubate for 2-3 minutes at 37°C. Add 10 ml pre-warmed Growth Medium 1G and gently pipette up and down to dissociate cell clumps. Transfer cells to a 15 ml conical tube and centrifuge at 200 x g for 5 minutes. Remove the medium and re-suspend cells in 10 ml of pre-warmed Growth Medium 1G. Dispense 5 ml of the cell suspension into a new T75 flask containing 10 ml pre-warmed Growth Medium 1G. Incubate cells in a humidified 37°C incubator with 5% CO₂. Freeze cells in freezing medium (10% DMSO in FBS) when cells reach 90% confluency.

Additional reagents required for this assay:

CD137L protein (BPS Bioscience, #71189)
Anti-CD137 Agonist Antibody (BPS Bioscience, #79097-2)
ONE-Step Luciferase Detection Reagent (BPS Bioscience, #60690)
Thaw Medium 1 (BPS Bioscience, #60187)
Growth Medium 1G (BPS Bioscience, #79544)

Assay Protocol

1. Harvest CD137/NF-κB reporter-HEK293 cells from culture in Growth Medium 1G and seed cells at a density of ~35,000 cells per well into a white clear-bottom 96-well microplate in 90 µl of Thaw medium 1. Incubate the plate at 37°C in a CO₂ incubator.
2. 24 hours after seeding, serially dilute the CD137L protein and anti-CD137 agonist antibody in Thaw Medium 1. For an EC₅₀ curve, we recommend a range of approximately 1 ng/ml to 10 µg/ml, final concentration.

Set up each treatment in at least triplicate
Add 10 µl diluted CD137L to the treated wells.
Add 10 µl Thaw Medium 1 to control wells.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**


Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com

- Add 100 μ l Thaw Medium 1 to cell-free control wells (for determining background luminescence)
- Incubate the plate at 37°C in a CO₂ incubator for ~ 6 hours for CD137L and 24 hrs for anti-CD137 agonist antibody.
 - Perform luciferase assay by using the ONE-Step luciferase assay system: Briefly, add 100 μ l of One-Step Luciferase reagent per well and rock at room temperature for ~30 minutes. Measure luminescence using a luminometer.
 - Data Analysis: Subtract the average background luminescence (cell-free control wells) from the luminescence reading of all wells. The fold induction of NF- κ B luciferase reporter expression = background-subtracted luminescence of stimulated well / average background-subtracted luminescence of unstimulated control wells.

Figure 1. Dose response of CD137/ NF- κ B-reporter HEK293 cells


CD137 ligand (BPS Bioscience, #71189) and anti-CD137 agonist (BPS Bioscience, #79097) were diluted and added to the cells (BPS Bioscience, #79289), then incubated at 37°C cell culture incubator for 6 and 24 hrs. respectively. After the treatment, perform Luciferase assay using One-Step Luciferase assay system (BPS Bioscience, #60690)

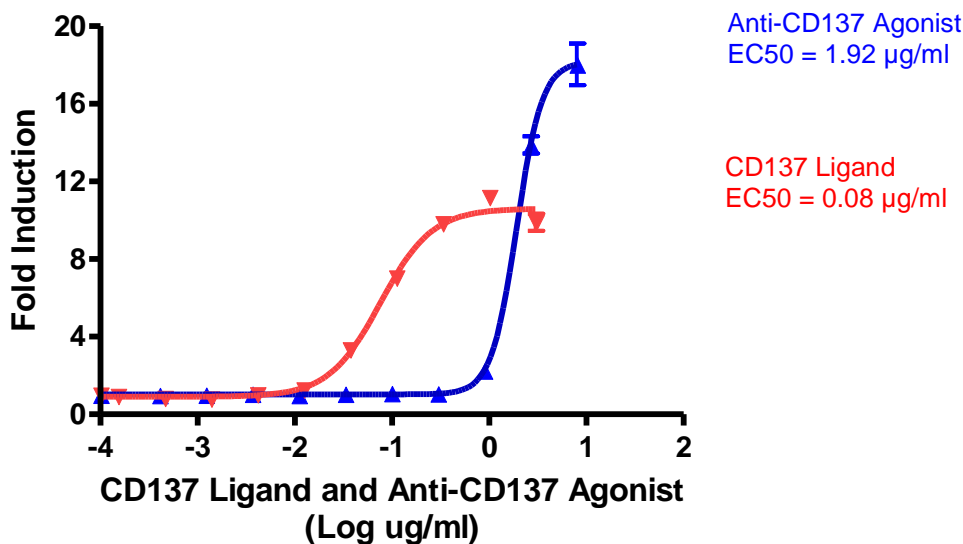

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.
 To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**
 Or you can Email us at: info@bpsbioscience.com
 Please visit our website at: www.bpsbioscience.com

Quality Assurance

Figure 2. Expression of CD137 (4-1BB) protein validated by flow cytometry. Flow cytometry showed PE-conjugated anti-human 4-1BB antibody (Clone 4B4-1, Biolegend, #309803) detects 4-1BB expression cells (green), using wild-type HEK293 cells as a negative control (blue).

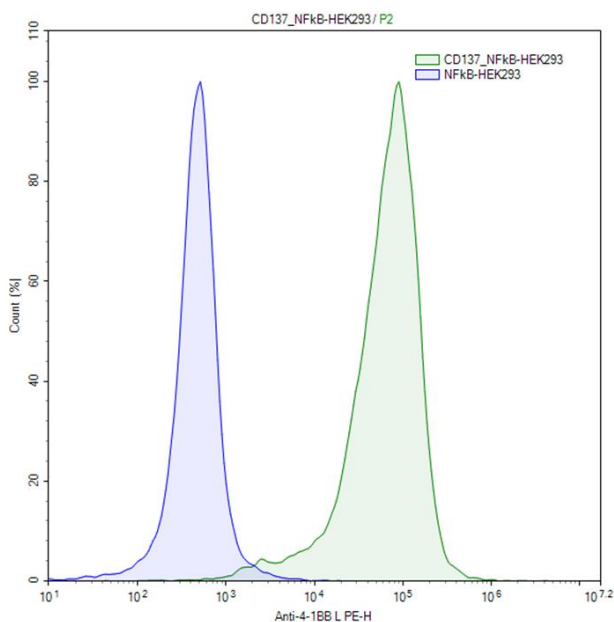

Application References

1. McNamara II JO *et.al.* (2008) Multivalent 4-1BB binding aptamers co-stimulate CD8+ T cells and inhibit tumor growth in mice. *J. Clin. Invest.* **118**: 376-386
2. Ma BY *et.al.* (2005) The expression and regulatory role of OX40 and 4-1BB heterodimers in activated human T cells. *Blood.* **106**: 2002-2010.

Vector and Sequence

Human CD137 (4-1BB), GenBank Accession #NM_001561, was cloned into pIRESHyg.

```

MGNSCYNIVATLLLLVLNFERTRSLQDPCSNCPAGTFCDNNRNQICSPCPPNSFSSAGGQRTCD
ICRQCKGVFRTRKECSSTSNAECDCTPGFHCLGAGCSMCEQDCKQGQELTKKGCKDCCFGT
FNDQKRGICRPWTNCSLDGKSVLVNGTKERDVVCGPSPADLSPGASSVTPAPAREPGHSPQ
IISFFLALTSTALLFLFLLTLRFSVVKRGRKLLYIFKQPFMRPVQTTQEEDGCSCRFPEEEEEGG
CEL
  
```

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com



6042 Cornerstone Court W, Ste B
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

License Disclosure: Purchase of this cell line grants you with a 10-year license to use this cell line in your immediate laboratory, for research use only. This license does not permit you to share, distribute, sell, sublicense, or otherwise make the cell line available for use to other laboratories, departments, research institutions, hospitals, universities, or biotech companies. The license does not permit use of this cell line in humans or for therapeutic or drug use. The license does not permit modification of the cell line in any way. Inappropriate use or distribution of this cell line will result in revocation of the license and result in an immediate cease of sales and distribution of BPS products to your laboratory. BPS does not warrant the suitability of the cell line for any particular use, and does not accept any liability in connection with the handling or use of the cell line. Modifications of this cell line, transfer to another facility, or commercial use of the cells may require a separate license and additional fees; contact sales@bpsbioscience.com for details. Publications using this cell line should reference BPS Bioscience, Inc., San Diego.

Related Products

	<u>Cat. #</u>	<u>Size</u>
Anti-CD137 Agonist Antibody	79097-2	100 µg
CD137L (4-1BBL) CHO-K1 Recombinant Cell Line	60523	2 vials
CD137 (4-1BB) HEK293 Recombinant Cell Line	60691	2 vials
NF-κB Reporter (Luc)-HEK293 cell line	60650	2 vials
CD137L, His-tag	71189	100 µg
CD137, Fc fusion (mIgG2a), Avi-tag (Mouse) HiP™	71254	100 µg
CD137, Fc fusion (mIgG2a), Biotin-labeled (Mouse) HiP™	71255	50 µg
CD137, Fc fusion (hIgG1) (Mouse)	71250	100 µg
CD137 (4-1BB), Fc fusion (Human) HiP™	71170	100 µg
CD137, Fc fusion (mIgG2a), Biotin-labeled (Mouse) HiP™	71171	50 µg
CD137L, His-Tag (Mouse)	72517	100 µg
ONE-Step™ Luciferase Assay System	60690-1	10 ml
ONE-Step™ Luciferase Assay System	60690-2	100 ml
CD137[Biotinylated]:CD137L Inhibitor Screening Kit	72025	96 rxns
Thaw Medium 1	60187	100 ml

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com