



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.202.1401  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

## Data Sheet

### **TrkA / SRE Reporter – HEK293 Recombinant Cell Line**

### **Catalog #: 79798**

#### **Product Description**

Recombinant HEK-293 cells expressing firefly luciferase gene under the control of Serum Response Elements (SRE) with constitutive expression of human TrkA (Tropomyosin receptor kinase A; NTRK1; ref. seq. NM\_002529.3).

#### **Background**

A chronic pain mediator, NGF (Nerve Growth Factor) binds to TrkA resulting in activation of the downstream signaling pathway linked to pro-nociception; therefore, targeting the kinase activity of TrkA and the interaction between NGF/TrkA have been attractive tools in pain management research. In addition to control of chronic pain, a recent clinical success has proven that inhibition of TrkA kinase would be a promising anti-cancer strategy.

#### **Application**

- Screen for inhibitors of NGF/TrkA signaling in a cellular context
- Screen for inhibitors of TrkA kinase activity in a cellular context

#### **Format**

Each vial contains  $2 \times 10^6$  cells in 1 ml of 10% DMSO in FBS

#### **Storage**

Immediately upon receipt, store in liquid nitrogen.

#### **Mycoplasma Testing**

The cell line has been screened using the MycoAlert™ Mycoplasma Detection kit (Lonza, #LT07-218) to confirm the absence of *Mycoplasma* species.

#### **General Culture Conditions**

**Thaw Medium 1 (BPS Bioscience, #60187):** MEM medium (Hyclone, #SH30024.01) supplemented with 10% FBS (Invitrogen, #26140-079), 1% non-essential amino acids (Hyclone, #SH30238.01), 1 mM Na pyruvate (Hyclone, #SH30239.01), 1% Penicillin/Streptomycin (Hyclone, SV30010.01)

**Growth Medium 1A (BPS Bioscience, #79528):** Thaw Medium 1 (BPS Bioscience, #60187) plus 400 µg/ml of Geneticin (Life Technologies #11811031) and 100 µg/ml of Hygromycin B (Life Technologies, #10687-010)

Cells should be grown at 37°C with 5% CO<sub>2</sub> using Growth Medium 1A.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.202.1401  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

**Assay Medium 1A (BPS Bioscience #79805):** Opti-MEM (Thermo Fisher, #31985070) plus 1% Penicillin/Streptomycin (Hyclone, #SV30010.01).

**To thaw the cells,** it is recommended to quickly thaw the frozen cells from liquid nitrogen in a 37°C water-bath, then transfer the entire contents of the vial to a tube containing 10 ml of Thaw Medium 1 (**no Geneticin and Hygromycin B**). Spin down the cells, remove supernatant and resuspend cells in 5 ml of pre-warmed Thaw Medium 1 (**no Geneticin and Hygromycin B**). Transfer the resuspended cells to a T25 flask and incubate at 37°C in a 5% CO<sub>2</sub> incubator. After 24 hours of culture, add an additional ~3 ml of Thaw Medium 1 (**no Geneticin and Hygromycin B**), and continue growing culture in a CO<sub>2</sub> incubator at 37°C until the cells are ready to be split. Cells should be split before they reach ~2.5 x 10<sup>6</sup> cells/ml. At first passage, switch to Growth Medium 1A (**contains Geneticin and Hygromycin B**).

**To passage the cells,** remove the medium, rinse cells with phosphate buffered saline (PBS), and detach cells from culture vessel with 0.05% Trypsin/EDTA. After detachment, add Growth Medium 1B (contains Geneticin) and transfer to a tube, spin down cells, resuspend cells in Growth Medium 1B (contains Geneticin) and seed appropriate aliquots of cell suspension into new culture vessels. Sub cultivation ration: 1:5 to 1:10 weekly or twice a week.

*Note: Just after thawing, the cells may grow at a slower rate. It is recommended to split the cells with ~ 1:4 ratio at the beginning of culturing. After several passages, the cell growth rate increases and the cells can be split with higher ratio.*

**To freeze down the cells,** remove the medium, rinse cells with phosphate buffered saline (PBS), and detach cells from culture vessel with 0.05% Trypsin/EDTA. After detachment, add Thaw Medium 1 (no Geneticin) and count the cells, then transfer to a tube, spin down cells, and resuspend in 4°C Freezing Medium (10% DMSO + 90% FBS) to ~2 x 10<sup>6</sup> cells/ml. Dispense 1 ml of cell aliquots into cryogenic vials. Place vials in an insulated container for slow cooling and store at -80°C overnight. Transfer to liquid nitrogen the next day for storage.

It is recommended to expand the cells and freeze down more than 10 vials of cells at an early passage for future use.


OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)

## Functional Validation and Assay Performance

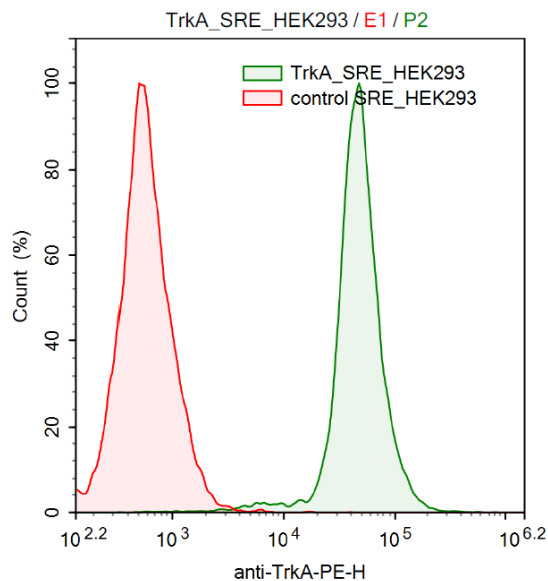

TrkA/SRE-HEK293 cells (green) or control SRE-HEK293 cells (red) were stained with PE-labeled Anti-TrkA Antibody (R&D systems, #FAB1751P) and analyzed by FACS. Y-axis is the % cell number. X-axis is the intensity of PE.

The following assays are designed for 96-well format. To perform the assay in a different format, the cell number and reagent volume should be scaled appropriately.

### Materials Required but Not Supplied

- Human NGF (R&D Systems, #256-GF-100)
- Media: Thaw Medium 1 (BPS Bioscience, #60187), Growth Medium 1A (BPS Bioscience, #79528), Assay Medium 1A (BPS Bioscience, #79805)
- 96-well tissue culture treated white clear-bottom assay plate (Corning, #3610)
- ONE-Step™ Luciferase Assay System (BPS Bioscience, #60690)
- Luminometer

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**


Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)

### A. Dose response of TrkA/SRE Reporter-HEK293 cells to human NGF

1. Harvest TrkA/SRE Reporter-HEK293 cells from culture in the Growth Medium 1A and seed cells at a density of ~50,000 cells per well into a white clear-bottom 96-well microplate in 100  $\mu$ l of Assay Medium 1A. Leave a couple of wells empty for use as a cell-free control.
2. Incubate cells at 37°C in a CO<sub>2</sub> incubator for ~ 16 hours.
3. Add 10  $\mu$ l of two-fold serial dilution of human NGF protein (R&D Systems, #256-GF-100) in the Assay Medium 1A to stimulated wells. Add 10  $\mu$ l of Assay Medium 1A to the unstimulated control wells. Add 100  $\mu$ l of Assay Medium 1A to cell-free control wells (for determining background luminescence). Set up each treatment in at least triplicate. Incubate the plate at 37°C in a CO<sub>2</sub> incubator for 6 hours.
4. Perform luciferase assay using ONE-Step Luciferase Assay buffer, according to the recommended instructions: Add 100  $\mu$ l of the final ONE-Step™ Luciferase reagent per well and rock at room temperature for ~15 to 30 minutes. Measure luminescence using a luminometer.
5. Data Analysis: Subtract average background luminescence (cell-free control wells) from luminescence reading of all wells.

### Example: Dose response of TrkA/SRE-HEK293 cells to human NGF

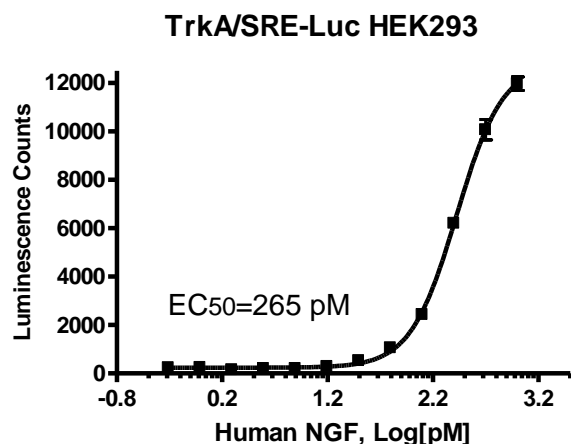

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**


Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)

## B. Inhibition of kinase activity of TrkA in TrkA/SRE Reporter-HEK293 cells

1. Harvest TrkA/SRE Reporter-HEK293 cells from culture in the Growth Medium 1A and seed cells at a density of ~50,000 cells per well into a white clear-bottom 96-well microplate in 80  $\mu$ l of Assay Medium 1A. Leave a couple of wells empty for use as a cell-free control.
2. Incubate cells at 37°C in a CO<sub>2</sub> incubator for ~ 16 hours.
3. Add 10  $\mu$ l TrkA inhibitor or vehicle control in 10  $\mu$ l of Assay Medium 1A to the wells. (Do not exceed a final DMSO concentration above 0.1%)
4. Incubate cells at 37°C in a CO<sub>2</sub> incubator for ~ 1 hours. Add 10  $\mu$ l NGF protein at final concentration of EC<sub>50</sub> in the Assay Medium 1A to the wells.
5. Add 100  $\mu$ l of Assay Medium 1A to cell-free control wells (for determining background luminescence). Set up each treatment in at least triplicate.
6. Incubate the plate at 37°C in a CO<sub>2</sub> incubator for 5 hours.
7. Perform luciferase assay using ONE-Step™ Luciferase Assay kit, according to the recommended instructions: Add 100  $\mu$ l of the final ONE-Step™ Luciferase reagent per well and rock at room temperature for ~15 to 30 minutes. Measure luminescence using a luminometer.
8. Data Analysis: Subtract average background luminescence (cell-free control wells) from luminescence reading of all wells.

### Example: Inhibition of TrkA by Larotrectinib (LOXO101) in TrkA/SRE-HEK293 cells

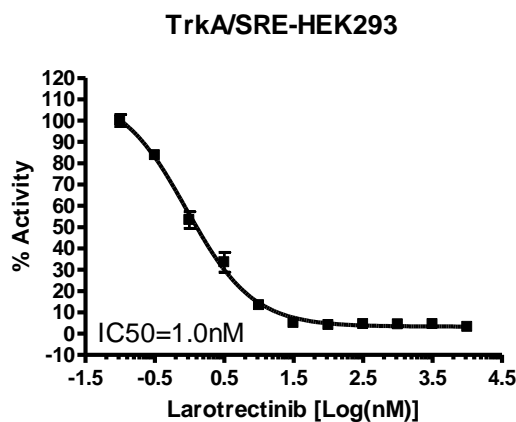

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**


Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)

### C. Inhibition of NGF/TrkA signaling in TrkA/SRE Reporter-HEK293 cells

1. Harvest TrkA/SRE Reporter-HEK293 cells from culture in the Growth Medium 1A and seed cells at a density of ~50,000 cells per well into a white clear-bottom 96-well microplate in 90  $\mu$ l of Assay Medium 1A. Leave a couple of wells empty for use as a cell-free control.
2. Incubate cells at 37°C in a CO<sub>2</sub> incubator for ~ 16 hours.
1. Prepare 100  $\mu$ l serially diluted anti-NGF antibody plus NGF (at 10X EC<sub>50</sub> concentration) solution in the Assay Medium 1A and incubate for 45 minutes at room temperature.
2. Add 10  $\mu$ l the anti-NGF antibody/NGF mixture to the cells and incubate the plate at 37°C in a CO<sub>2</sub> incubator for 5 hours.
3. Perform luciferase assay using ONE-Step™ Luciferase Assay kit, according to the recommended instructions: Add 100  $\mu$ l of the final ONE-Step™ Luciferase reagent per well and rock at room temperature for ~15 to 30 minutes. Measure luminescence using a luminometer.
4. Data Analysis: Subtract average background luminescence (cell-free control wells) from luminescence reading of all wells.

#### Example: Inhibition of NGF/TrkA signaling by anti-NGF antibody in TrkA/SRE-HEK293 cells

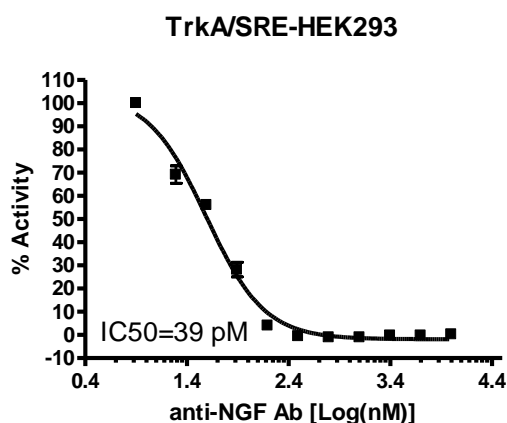

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.202.1401  
**Fax:** 1.858.481.8694  
**Email:** info@bpsbioscience.com

### Sequence

Human TrkA sequence (accession number NM\_002529.3)

MLRGRRGQLGWHWSAAGPGSLLAWLILASAGAAPCPDACC PHGSSGLRCTR DGALDSLHH  
LPGAENLT ELYIENQQHLQHLELRDLRGLGELRNLTIVKSGLR FVAPDAFHFTPRLSRLNLSFNA  
LESLSWKTVQGLSLQELVLSGNPLHCSCALRWLQRWEE EGLGGVPEQKLQCHGQGPLAHMP  
NASC GVP TLKVQVPNASVDVGDDVLLRCQVEGRGLEQAGWILTELEQSATVMKSGGLPSLGL  
TLANVTS DLNRKNVTCWAENDVGRAEVS VQVNV SFPASVQLHTAVEMHHWCIPFSVDGQPAP  
SLRWLFNGSVLNETSFIFTEFLEPAANETVRHGCLRLNQPTHVNNGNYTLAANPFGQASASIM  
AAFMDNPF EFN PEDPIPVSFSPVD TNSTSGDPVEKKDET PFGVSVAVGLAVFACLFLSTLLLVL  
NKCGRRNKFGINRPAVLAPEDGLAMSLHFMTLGGSSLSPTEGKSGSLQGHIIENPQYFSDACV  
HHIKRRDIVLKWELGEGAFGKVFLAECHNLLPEQDKMLVAVKALKEASE SARQDFQREAE LLT  
MLQHQHIVRFFGVCTEGRPLL MVFEYMRHGD LNRFLRSHGPDAKLLAGGEDVAPGPLGLGQL  
LAVASQVAAGMVYLAGLHFVHRDLATRNCLVGGQLVVKIGDFGMSRDIYSTDYR VGGRTMLP  
IRWMPPE SILYRKFTTESDVWSFGVVLWEIFTY GKQPWYQLSNTEAIDCITQGRELERPRACPP  
EYAIMRGCWQREPPQQRHSIKDVHARLQALAQAPPVYLDVLG

### Related Products

<u>Product</u>	<u>Cat. #</u>	<u>Size</u>
TRKA, GST-tag	40280	10 µg
TRKB, GST-tag	40281	10 µg
TRKC, GST-tag	40282	10 µg
TRKC (G623R) Mutant, Active)	40203	10 µg
TRKC (G623E) Mutant, Active)	40204	10 µg
TRKC (G623R L686M), Mutant, Active	40215	10 µg
TRKC (L686M), Mutant, Active	40215	10 µg
Thaw Medium 1	60187	100 ml
Growth Medium 1A	79528	500 ml

### Notes

License Disclosure: Purchase of this cell line grants you with a 10-year license to use this cell line in your immediate laboratory, for research use only. This license does not permit you to share, distribute, sell, sublicense, or otherwise make the cell line available for use to other laboratories, departments, research institutions, hospitals, universities, or biotech companies. The license does not permit the use of this cell line in humans or for therapeutic or drug use. The license does not permit modification of the cell line in any way. Inappropriate use or distribution of this cell line will result in revocation of the license and result in an immediate cease of sales and distribution of BPS products to your laboratory. BPS does not warrant the suitability of the cell line for any particular use, and does not accept any liability in connection with the handling or use of the cell line. Modifications of this cell line, transfer to another facility, or commercial use of the cells may require a separate license and additional fees; contact sales@bpsbioscience.com for details. Publications using this cell line should reference BPS Bioscience, Inc., San Diego.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)