



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.202.1401  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

## **Data Sheet**

### **BAFF-R – CHO K1 Recombinant Cell Line**

#### **Medium Expression**

#### **Catalog # 79921-M**

#### **Background**

BAFF Receptor (B-cell activating factor receptor, BAFF-R) is an essential factor for B cell maturation and survival. In humans it is encoded by the TNFRSF13C gene. BAFF-R function is impressively documented in humans by a homozygous deletion, which leads to an almost complete block of B cell development at the stage of immature/transitional B cells. The resulting immunodeficiency is characterized by B-lymphopenia, agammaglobulinemia, and impaired humoral immune responses.

#### **Description**

Recombinant human BAFF-R-CHO K1 cell line stably expressing full length human BAFF-R/TNFRSF13C receptor (GenBank Accession: NM\_052945.3). Surface expression of BAFF-R was confirmed by flow cytometry. Each stable clonal cell line was selected for different levels of BAFF-R expression (High, Medium, Low) to mimic different stages of cancer target cells with various BAFF-R expression levels.

#### **Application**

- Screen for activators or inhibitors of antibody-mediated signaling for immunotherapy research and drug discovery.
- Characterize BAFF-R antibodies and ligands for binding assays.

#### **Format**

Each vial contains ~ 2 x 10<sup>6</sup> cells in 1 ml of 10% DMSO in FBS.

#### **Storage**

Store in liquid nitrogen immediately upon receipt.

#### **Culture Medium**

**Thaw Medium 3 (BPS Bioscience, #60186):** Ham's F-12 medium (Hyclone, #SH30526.01) supplemented with 10% FBS (Life technologies, #26140-079), 1% Penicillin/Streptomycin (Hyclone, #SV30010.01).

**Growth Medium 3B (BPS Bioscience, #79529):** Thaw Medium 3 (BPS Bioscience, #60186) plus 500 ug/ml Hygromycin B (Invitrogen, #10687010)

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)

### Recommended Culture conditions


**Frozen Cells:** Prepare T-75 culture flask with 20 ml of pre-warmed Thaw Medium 3. Quickly thaw cells in a 37°C water bath with constant and slow agitation. After cleaning the outside of the vial with 70% ethanol, immediately transfer the entire content to Thaw Medium 3 (no Hygromycin). Avoid pipetting up and down, and gently rock the flask to distribute the cells. Incubate the cells in a humidified 37°C incubator with 5% CO<sub>2</sub>. The next day, change to fresh Growth Medium 3B (with Hygromycin) without disturbing the attached cells. Continue to incubate until cells reach desired confluency. If slow cell growth occurs during resuscitation, increase FBS to 15% for the first week of culture.

**Subculture:** When cells reach 90% confluency, remove the medium and wash twice with PBS (without Magnesium or Calcium). Treat cells with 2 ml of 0.25% trypsin/EDTA and incubate for 3 minutes at 37°C. After confirming cell detachment by light microscopy, add 10 ml of prewarmed medium and gently pipette up and down to dissociate cell clumps. Transfer cells to a 15 mL conical tube and centrifuge at 200 x g for 5 minutes. Remove the medium and resuspend cells in 10 ml pre-warmed growth medium. Dispense 1 mL of the cell suspension into a new T75 flask containing pre-warmed 19 ml complete medium (a subcultivation ratio of 1:10 to 1:20 is recommended). Incubate cells in a humidified 37°C incubator with 5% CO<sub>2</sub>.

To freeze cells, resuspend cell pellet in freezing medium (10% DMSO in FBS). BPS Bioscience recommends freezing down at least 10 vials of cells at an early passage number for future use.

### Validation

Cell surface expression of human BAFF-R in CHO K1 cells was confirmed by flow cytometry.

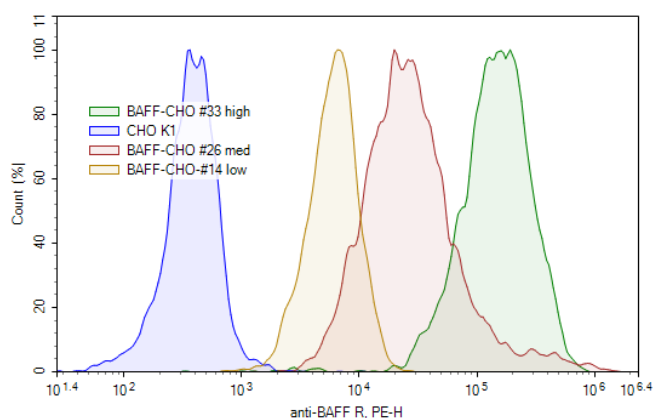

**Figure 1.** Flow cytometry analysis of cell surface expression of BAFF-R in CHO K1 cells. BAFF-R-CHO K1 cells with high (green), medium, (red), or low (orange) expression levels or control CHO K1 cells (blue) were stained with PE-labeled anti-human BAFF-R antibody (Biolegend, #316905) and analyzed by FACS. Y-axis is the % cell number. X-axis is the intensity of PE.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.202.1401  
**Fax:** 1.858.481.8694  
**Email:** info@bpsbioscience.com

### **Mycoplasma Testing**

This cell line has been screened using the MycoAlert™ Mycoplasma Detection Kit (Lonza, #LT07-118) to confirm the absence of Mycoplasma contamination. MycoAlert Assay Control Set (Lonza, #LT07-518) was used as a positive control.

### **Vector and Sequence**

Human BAFF-R/TNFRSF13C (accession number: NM\_052945.3) was cloned into pCMV3.

MRRGPRSLRGRDAPAPTPCVPAECFDLLVRHCVACGLLRTPRKPAGASSPAPRTALQPQES  
VGAGAGEAALPLPGLLFGAPALLGLALVLALVLVGLVSWRRRQRRLRGASSAEAPDGDKDAPE  
PLDKVIILSPGISDATAPAWPPPGEDPGTTTPPGHSPVPATELGSTELVTTKTAGPEQQ

**License Disclosure:** Purchase of this cell line grants you with a 10-year license to use this cell line in your immediate laboratory, for research use only. This license does not permit you to share, distribute, sell, sublicense, or otherwise make the cell line available for use to other laboratories, departments, research institutions, hospitals, universities, or biotech companies. The license does not permit use of this cell line in humans or for therapeutic or drug use. The license does not permit modification of the cell line in any way. Inappropriate use or distribution of this cell line will result in revocation of the license and result in an immediate cease of sales and distribution of BPS products to your laboratory. BPS does not warrant the suitability of the cell line for any particular use and does not accept any liability in connection with the handling or use of the cell line. Modifications of this cell line, transfer to another facility, or commercial use of the cells may require a separate license and additional fees; contact sales@bpsbioscience.com for details. Publications using this cell line should reference BPS Bioscience, Inc., San Diego.

### **Related Product**

	<b><u>Cat. #</u></b>	<b><u>Size</u></b>
Thaw Medium 3	60186	100ml
Growth Medium 3B	79529	500ml
BAFF-R-CHO Recombinant Cell Line (Medium Expression)	79921-M	2 vials
BAFF-R-CHO Recombinant Cell Line (Low Expression)	79921-L	2 vials
BAFF-R-CHO Recombinant Cell Line (High Expression)	79921-H	2 vials
BAFF-R(CD268)	90103-A	10 µg
BAFF, His-Avi-Tag HiP™	100194	50 µg
BAFF-R, Fc-fusion (IgG1), Avi-Tag	100286	100 µg
BAFF-R, Fc-fusion (IgG1), Avi-Tag, Biotinylated	100287-1	25 µg

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)