



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

**Description**

AAV-DJ SP-B Luciferase particles transduce firefly (*Photinus pyralis*) luciferase under the control of a SP-B (Surfactant protein B) promoter that drives luciferase reporter expression mainly in epithelial lung cells.

**Background:**

Adeno-Associated Virus-DJ (AAV-DJ) is a synthetic serotype made from eight different wild-type AAV serotypes (AAV2, 4, 5, 8, 9, avian, bovine, and goat AAV) using DNA shuffling. These modifications allow the AAV-DJ serotype to exhibit improved transduction efficiency *in vitro* and *in vivo* and infect a broader range of cell types compared to the wild-type serotypes.

Surfactant protein B (SP-B) is a developmentally and hormonally regulated lung protein that is required for normal surfactant function. SP-B is selectively expressed in bronchiolar and alveolar epithelial cells of the lung and it is found to be associated with lipids. Surfactant is a mixture of proteins and lipids that allows proper inflation of the lung and breathing. Deficiency in SP-B results in lethality at birth. The human SP-B promoter fragment directs epithelial lung-cell specific transgene expression.

**Application(s)**

- Positive control in the transduction of lung epithelial cells.
- Optimization of transduction assays and tracking of transgene expression over time.

**Serotype**


AAV-DJ

**Formulation**

AAV was produced in HEK293-AAV cells and is supplied in PBS-MK (PBS Magnesium-Potassium) buffer containing 0.01% Pluronic F68. Virus particles can be packaged in custom formulations by special request, for an additional fee.

**Purification**

The purity of the AAV particles was confirmed to be greater than 90% by staining with One-Step Lumitein™ UV Protein Gel Stain (Biotium, 21005-1L). The purity varies with each lot; the exact value will be provided with each shipment.

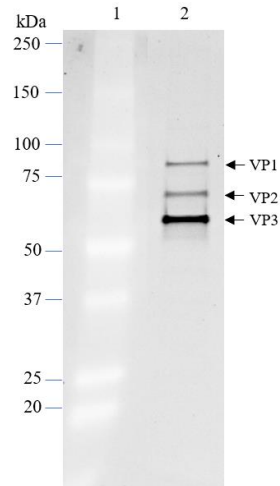

*Figure 1. Purified AAV-DJ SP-B- Luciferase particles.*

Staining of a 4-20% SDS-PAGE gel. The protein marker was loaded in lane 1, and  $2 \times 10^9$  VG (vector genomes) of AAV was loaded in lane 2. AAV viral proteins VP1, VP2, and VP3 are indicated by arrows.

#### Size and Titer

Two vials ( $50 \mu\text{l} \times 2$ ) of AAV at a titer  $\geq 1 \times 10^{12}$  vector genomes/ml. The titer is determined by qPCR and varies with each lot; the exact value will be provided with each shipment.

#### Storage




AAV is shipped with dry ice. For long-term storage, it is recommended to store AAV at  $-80^\circ\text{C}$ . Avoid repeated freeze-thaw cycles. Titers can drop significantly with each freeze-thaw cycle.

#### Biosafety



Recombinant AAV is inherently replication-deficient and not known to cause any human diseases. Additionally, following transduction, AAV vectors exist episomally and do not integrate into or disrupt the host cell's genome. AAV requires the use of a Biosafety Level 1 facility. BPS Bioscience recommends following all local, federal, state, and institutional regulations and using all appropriate safety precautions.

## Validation Data

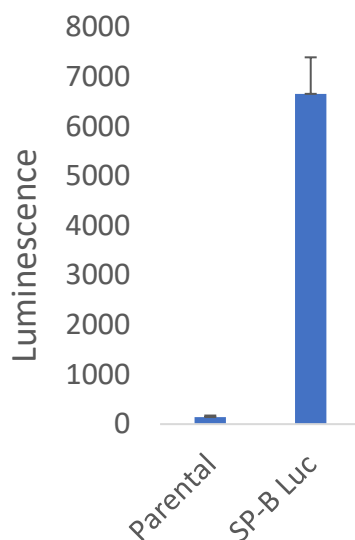

*Figure 2. Luciferase activity in HEK293 cells transduced with AAV-DJ SP-B-Luciferase particles.  $1 \times 10^5$  cells/well were transduced in a 6-well plate with AAV-DJ SP-B-Luciferase at an MOI of  $2 \times 10^4$ . After 72 hours of transduction, transduced cells or parental HEK293 cells were re-seeded in a 96-well plate at a density of  $2 \times 10^4$  cells/well, and luciferase activity was measured with ONE-Step™ Luciferase Assay System (BPS Bioscience #60690).*

*Data shown is representative. For lot-specific information, please contact BPS Bioscience, Inc. at [support@bpsbioscience.com](mailto:support@bpsbioscience.com).*

## Notes

The AAV-DJ viruses are covered under several patents, including U.S. Patent Nos. 7,588,772, 8,067,014, 8,574,583, and 8,906,387, as well as corresponding foreign patents applications and patent rights. AAV-DJ is used under a license agreement.

## Troubleshooting Guide

Visit [bpsbioscience.com/lentivirus-faq](https://bpsbioscience.com/lentivirus-faq) for detailed troubleshooting instructions. For all further questions, please email [support@bpsbioscience.com](mailto:support@bpsbioscience.com).

## Related Products

<i>Products</i>	<i>Catalog #</i>	<i>Size</i>
AAV-DJ SaCas9	78478	50 $\mu$ l x 2
AAV-DJ Luciferase-eGFP	78460	50 $\mu$ l x 2
AAV-DJ Luciferase-mCherry	78469	50 $\mu$ l x 2
AAV-DJ MBP-eGFP	82112	50 $\mu$ l x 2
AAV-DJ SYN1-Luciferase	82134	50 $\mu$ l x 2