



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Coriolis[®] Compact



Portable dry cyclonic air sampler

- Concentrated sample compatible with any downstream analysis
- Light & compact for indoor & outdoor use
- High efficiency with up to 8 hours of autonomy



CORIOLIS® COMPACT

PORTABLE DRY AIR SAMPLER

FOR BIO CONTAMINATION CONTROL


The Coriolis® Compact is a dry cyclonic collector intended for microbial air monitoring. Light and compact, it can be easily transported and positioned on strategic places for indoor and outdoor sampling.

Its innovative dry cyclonic technology aspirates the particles with an airflow of 50L/m & centrifuges them in a cone. The collected particles or microorganisms can then be recovered by rinsing consumable with appropriate buffer

solution, so the sample is compatible with multiple downstream analysis (NGS, qPCR, Culture). The Coriolis® Compact gives access to highly sensitive and specific results in only a few hours.

Due to an improved battery autonomy, it performs up to 8 hours of continuous sampling. Bertin provides single use consumable to avoid cross contamination.

Innovative dry cyclonic technology

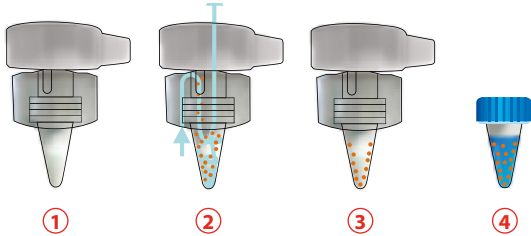

1. The cone is placed on the device and locked
2. Aspirated air runs into the cone & forms a vortex
3. Particles in the air are centrifuged in the cone's walls
4. Collected particles are recovered by rinsing the cone.

Applications



Pollution & Environment



Food / Pharma / Veterinary / Industry



Biomedical & Health



Research & Development

Technical features

Dimensions	255x135x130 mm
Weight	1.2 kg. (battery included)
Collection time	From 1 minute to 8 hours
Temperature range	+5°C to +45°C with a humidity range of 10% to 90%
Flow	50 L./min
Size of particles collected	From 500 nm. to 10 µm
Type of particles collected	Bacteria, viruses, molds...
Autonomy on battery	Max 8 hours in operation at 20°C
Decontamination	H ₂ O ₂
Catalogue reference	P002114-CORC0-A.0

Discover our comprehensive range of solutions

