

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



PRODUCT INFORMATION



Ionomycin

Item No. 10004974

CAS Registry No.: 56092-81-0

Formal Name: (4R,6S,8S,10Z,12R,14R,16E,18R,1

> 9R,20S,21S)-11,19,21-trihydroxy-4,6,8,12,14,18,20-heptamethyl-22-[(2S,2'R,5S,5'S)-octahydro-5'-[(1R)-1hydroxyethyl]-2,5'-dimethyl[2,2'-bifuran]-

5-yl]-9-oxo-10,16-docosadienoic acid


MF: $C_{41}H_{72}O_{9}$ 709.0 FW: **Purity:** ≥95% UV/Vis.:

 λ_{max} : 277 nm Supplied as: A solution in ethanol

-20°C Storage: Stability: ≥2 years

Item Origin: Bacterium/Streptomyces conglobatus

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

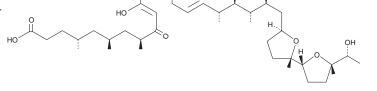

Laboratory Procedures

lonomycin is supplied as a solution in ethanol. To change the solvent, simply evaporate the ethanol under a gentle stream of nitrogen and immediately add the solvent of choice. Solvents such as DMSO and dimethyl formamide purged with an inert gas can be used. The solubility of ionomycin in these solvents is approximately 1.4 and 2.5 mg/ml, respectively.

lonomycin is sparingly soluble in aqueous buffers. For maximum solubility in aqueous buffers, the ethanolic solution of ionomycin should be diluted with the aqueous buffer of choice. Ionomycin has a solubility of approximately 0.1 mg/ml in a 1:10 solution of ethanol:PBS (pH 7.2) using this method. We do not recommend storing the aqueous solution for more than one day.

Description

Ionomycin is a calcium ionophore originally isolated from S. conglobatus. 1 It induces intracellular calcium mobilization in X. laevis oocytes when used at a concentration of 5 μM. Ionomycin (5 μM) induces the formation of neutrophil extracellular traps (NETs) in isolated human neutrophils. 2 Ionomycin in combination with phorbol 12-myristate 13-acetate (PMA; Item No. 10008014) has commonly been used to induce the proliferation and activation of T cells, as well as to stimulate cytokine production in macrophages, in vitro.^{3,4}

References

- 1. Liu, W.-C., Slusarchyk, D.S., Astle, G., et al. lonomycin, a new polyether antibiotic. J. Antibiot. 31(9), 815-819 (1978).
- 2. Locke, M., Francis, R.J., Tsaousi, E., et al. Fibrinogen protects neutrophils from the cytotoxic effects of histones and delays neutrophil extracellular trap formation induced by ionomycin. Sci. Rep. 10(1), 11694 (2020).
- 3. Lehnert, C., Weiswange, M., Jeremias, I., et al. TRAIL-receptor costimulation inhibits proximal TCR signaling and suppresses human T cell activation and proliferation. J. Immunol. 193(8), 4021-4031 (2014).
- Foey, A.D. and Brennan, F.M. Conventional protein kinase C and atypical protein kinase Cζ differentially regulate macrophage production of tumour necrosis factor-a and interleukin-10. Immunology 112(1), 44-53 (2004).

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website

Copyright Cayman Chemical Company, 08/11/2023

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640 CUSTSERV@CAYMANCHEM.COM WWW.**CAYMANCHEM**.COM