



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



N-butyryl-L-Homoserine lactone-d₅

Item No. 10007899

Formal Name: N-[(3S)-tetrahydro-2-oxo-3-furanyl]-butanamide-3,3,4,4,4-d₅
Synonyms: BHL-d₅, C₄-HSL-d₅, PAI-2-d₅
MF: C₈H₈D₅NO₃
FW: 176.2
Chemical Purity: ≥98% (N-butyryl-L-Homoserine lactone)
Deuterium Incorporation: ≥99% deuterated forms (d₁-d₅); ≤1% d₀
Supplied as: A crystalline solid
Storage: -20°C
Stability: ≥1 year

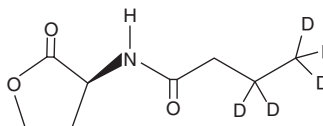

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

N-butyryl-L-Homoserine lactone-d₅ is intended for use as an internal standard for the quantification of N-butyryl-L-homoserine lactone (Item No. 10007898) by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

N-butyryl-L-Homoserine lactone-d₅ is supplied as a crystalline solid. A stock solution may be made by dissolving the N-butyryl-L-homoserine lactone-d₅ in the solvent of choice, which should be purged with an inert gas. N-butyryl-L-Homoserine lactone-d₅ is soluble in organic solvents such as ethanol, DMSO, and dimethyl formamide. The solubility of N-butyryl-L-homoserine lactone-d₅ in these solvents is approximately 30 mg/ml.

Description

N-butyryl-L-Homoserine lactone is a quorum-sensing signaling molecule produced by *P. aeruginosa*.¹ It induces expression of the virulence genes *lasB* and *rhlA* in *P. aeruginosa* when used at a concentration of 10 μM.² N-butyryl-L-Homoserine lactone (50 μM) induces rhamnolipid accumulation in *P. aeruginosa* growth media.³

References

1. Pearson, J.P., Passador, L., Iglewski, B.H., et al. A second N-acylhomoserine lactone signal produced by *Pseudomonas aeruginosa*. *Proc. Natl. Acad. Sci. USA*. **92**(5), 1490-1494 (1995).
2. Ikeda, T., Kajiyama, K., Kita, T., et al. The synthesis of optically pure enantiomers of N-acyl-homoserine lactone autoinducers and their analogues. *Chem. Lett.* **30**(4), 314-315 (2001).
3. Schooling, S.R., Charaf, U.K., Allison, D.G., et al. A role for rhamnolipid in biofilm dispersion. *Biofilm* **1**(2), 91-99 (2004).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 09/16/2021

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM