

Produktinformation



Forschungsprodukte & Biochemikalien
Zellkultur & Verbrauchsmaterial
Diagnostik & molekulare Diagnostik
Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien T. +43(0)1 489 3961-0 F. +43(0)1 489 3961-7 <u>mail@szabo-scandic.com</u> www.szabo-scandic.com


PRODUCT INFORMATION



CI-Amidine (hydrochloride)

Item No. 10599

Sold under license from the University of South Carolina under U.S. Patent No. 7,964,636.

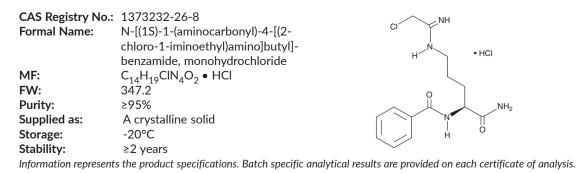

Laboratory Procedures

CI-Amidine (hydrochloride) is supplied as a crystalline solid. A stock solution may be made by dissolving the Cl-amidine (hydrochloride) in the solvent of choice. Cl-Amidine (hydrochloride) is soluble in organic solvents such as ethanol, DMSO, and dimethyl formamide, which should be purged with an inert gas. The solubility of Cl-amidine (hydrochloride) in these solvents is approximately 20, 50, and 14 mg/ml, respectively.

Further dilutions of the stock solution into aqueous buffers or isotonic saline should be made prior to performing biological experiments. Ensure that the residual amount of organic solvent is insignificant, since organic solvents may have physiological effects at low concentrations. Organic solvent-free aqueous solutions of CI-amidine (hydrochloride) can be prepared by directly dissolving the crystalline solid in aqueous buffers. The solubility of Cl-amidine (hydrochloride) in PBS, pH 7.2, is approximately 3 mg/ml. We do not recommend storing the aqueous solution for more than one day.

Description

Cl-Amidine is an irreversible inhibitor of protein arginine deiminases (PADs; $IC_{50}s = 0.8$, 6.2, and 5.9 μ M for PAD1, PAD3, and PAD4, respectively).¹⁻³ It is cytotoxic to HL-60, MCF-7, and HT-29 cancer cells $(IC_{50}s = 0.25, 0.05, and 1 \mu M, respectively)$.⁴ Cl-amidine (50 mg/kg) reduces *ex vivo* extracellular neutrophil extracellular trap (NET) formation and increases survival in a mouse model of sepsis induced by cecal ligation and puncture (CLP).⁶ It also decreases the citrulline content in serum and joints and reduces the development of IgG autoantibodies in a mouse model of collagen-induced arthritis in a dose-dependent manner.⁵

References

- 1. Luo, Y., Arita, K., Bhatia, M., et al. Biochemistry 45, 11727-11736 (2006).
- 2. Knuckley, B., Causey, C.P., Jones, J.E., et al. Biochem. 49(23), 1-28 (2010).
- 3. Muth, A., Subramanian, V., Beaumont, E., et al. J. Med. Chem. 60(7), 3198-3211 (2017).
- 4. Slack, J.L., Causey, C.P., and Thompson, P.R. Cell. Mol. Life Sci. 68(4), 709-720 (2011).
- 5. Willis, V.C., Gizinski, A.M., Banda, N.K., et al. J. Immunol. 186, 4396-4404 (2011).
- 6. Biron, B.M., Chung, C.-S., O'Brien, X.M., et al. J. Innate Immun. 9(1), 22-32 (2017).

WARNING THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFFTY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

subject to Cavman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information Buyer agrees to purchase the material can be found on our website.

Copyright Cayman Chemical Company, 01/18/2021

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA PHONE: [800] 364-9897 [734] 971-3335 FAX: [734] 971-3640 CUSTSERV@CAYMANCHEM.COM WWW.CAYMANCHEM.COM