



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



Rabbit anti-V5 tag IgG conjugated to R-Phycoerythrin

Product Number D5-1834
Lot Number XXXX
Amount 100 µg
Store at 2-8°C

Form/Shipping & Storage


Supplied as a lyophilized powder. Upon receipt, store at 2-8°C in the dark, do not freeze. Phycobiliproteins are sensitive to freeze-thaw cycles.

Handling

Reconstitute to 1.0 ml with distilled deionized water vortex gently and allow vial to sit on ice for 20 minutes. We recommend that the investigator determine the appropriate working concentration for their specific application. Avoid exposure to heat and light.

Buffer

Upon rehydration with 1.0 ml distilled deionized water; the product is in 100 mM sodium phosphate (pH 7.4), 50 mM sucrose, 100 mM sodium chloride, 0.1% BSA as a stabilizer, and 0.05% sodium azide as a preservative. The concentration of the conjugate is 100 µg/ml

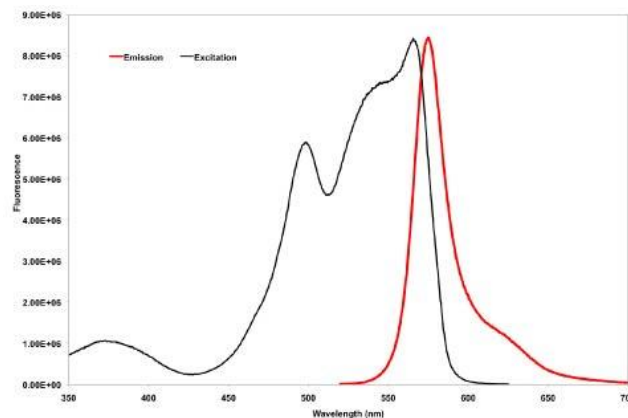

Fluorescence excitation and emission spectra of R-phycoerythrin in 100 mM sodium phosphate (pH 7.2) + 1 mM EDTA and 1 mM sodium azide. Emission scan was taken with excitation at 498 nm. Excitation scan was taken with emission at 575 nm. Scans were normalized to equalize peak heights.

Spectral Characteristics

Visible absorption maxima	565>540>498
Emission maximum	578

Concentration

After reconstitution to 1.0 ml

Fluor:Protein = ~1.2:1

Antigen Info

Antibody was raised against V5-tag with amino acid sequence GKPIPPLLGLDST.

References:

Nygaard NR, Bono C, Brown LR, Gorka J, Giacometto KS, Schaiff WT, Graham MB, McCourt DW, Kabeer M, Braciale VL, et al. Antibody recognition of an immunogenic influenza hemagglutinin-human leukocyte antigen class II complex. J Exp Med. 1991 Jul 1;174(1):243-51.

Rogers BE, Chaudhuri TR, Reynolds PN, Della Manna D, Zinn KR. Non-invasive gamma camera imaging of gene transfer using an adenoviral vector encoding an epitope-tagged receptor as a reporter. Gene Ther. 2003 Jan;10(2):105-14.