



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



RCS-4 N-pentanoic acid metabolite-d₅

Item No. 14365

CAS Registry No.: 2748469-51-2
Formal Name: 5-(3-(4-methoxybenzoyl)-1H-indol-1-yl)-2,4,5,6,7-d₅-pentanoic acid
MF: C₂₁H₁₆D₅NO₄
FW: 356.4
Chemical Purity: ≥98% (RCS-4 N-pentanoic acid metabolite)
Deuterium
Incorporation: ≥99% deuterated forms (d₁-d₅); ≤1% d₀
UV/Vis.: λ_{max}: 213, 264, 313 nm
Supplied as: A solution in acetonitrile
Storage: -20°C
Stability: ≥4 years

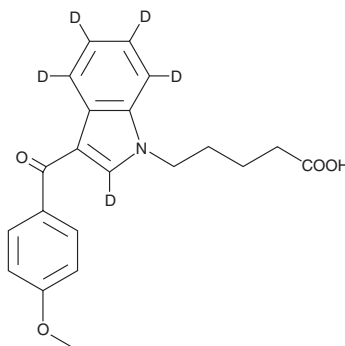

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Description

RCS-4 N-pentanoic acid metabolite-d₅ (Item No. 14365) is intended for use as an internal standard for the quantification of RCS-4 N-pentanoic acid metabolite by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

RCS-4 N-pentanoic acid metabolite is an expected phase I metabolite of RCS-4 (Item No. 10645), a synthetic cannabinoid that has been detected in herbal mixtures.¹ This product is intended for research and forensic applications

Reference

1. Kikura-Hanajiri, R., Uchiyama, N., and Goda, Y. Survey of current trends in the abuse of psychotropic substances and plants in Japan. *Leg. Med. (Tokyo)* **13**(3), 109-15 (2011).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 06/16/2023

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM