



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

# PRODUCT INFORMATION

## 5-(N-ethyl-N-isopropyl)-Amiloride

Item No. 14406

CAS Registry No.: 1154-25-2

Formal Name: 3-amino-N-(aminoiminomethyl)-6-chloro-5-[ethyl(1-methylethyl)amino]-2-pyrazinecarboxamide

Synonyms: EIPA, L-593,754, MH 12-43

MF: C<sub>11</sub>H<sub>18</sub>ClN<sub>7</sub>O

FW: 299.8


Purity: ≥98%

UV/Vis.: λ<sub>max</sub>: 230, 291, 380 nm

Supplied as: A crystalline solid

Storage: -20°C

Stability: ≥2 years

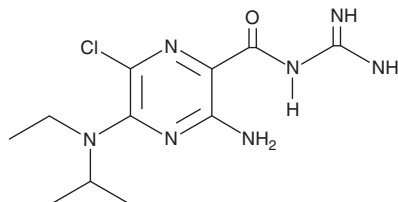

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

### Laboratory Procedures

5-(N-ethyl-N-isopropyl)-Amiloride (EIPA) is supplied as a crystalline solid. A stock solution may be made by dissolving the EIPA in the solvent of choice. EIPA is soluble in organic solvents such as ethanol, DMSO, and dimethyl formamide (DMF), which should be purged with an inert gas. The solubility of EIPA in these solvents is approximately 2, 14, and 20 mg/ml, respectively.

EIPA is sparingly soluble in aqueous buffers. For maximum solubility in aqueous buffers, EIPA should first be dissolved in DMF and then diluted with the aqueous buffer of choice. EIPA has a solubility of approximately 0.5 mg/ml in a 1:1 solution of DMF:PBS (pH 7.2) using this method. We do not recommend storing the aqueous solution for more than one day.

### Description

Sodium-hydrogen exchangers (NHE) are involved in maintaining sodium and pH balance in a variety of tissues. They are also known as sodium-hydrogen antiporters and solute carrier family 9 members. EIPA is a potent inhibitor of several NHE isoforms, inhibiting NHE1, NHE2, NHE3, and NHE5 with K<sub>i</sub> values of 0.02, 0.5, 2.4, and 0.42 μM, respectively.<sup>1,2</sup> It less effectively inhibits NHE4 (IC<sub>50</sub> ≥10 μM).<sup>1</sup> EIPA is commonly used at a concentration of 5-10 μM to inhibit cellular HNE activity.<sup>3,4</sup>

### References

1. Masereel, B., Pochet, L., and Laeckmann, D. An overview of inhibitors of Na<sup>+</sup>/H<sup>+</sup> exchanger. *Eur. J. Med. Chem.* **38**(6), 547-554 (2003).
2. Attaphitaya, S., Park, K., and Melvin, J.E. Molecular cloning and functional expression of a rat Na<sup>+</sup>/H<sup>+</sup> exchanger (NHE5) highly expressed in brain. *J. Biol. Chem.* **274**(7), 4383-4388 (1999).
3. Hwang, S.-M., Koo, N.-Y., Jin, M., et al. Intracellular acidification is associated with changes in free cytosolic calcium and inhibition of action potentials in rat trigeminal ganglion. *J. Biol. Chem.* **286**(3), 1719-1729 (2011).
4. Xue, L., Aihara, E., Wang, T.C., et al. Trefoil factor 2 requires Na/H exchanger 2 activity to enhance mouse gastric epithelial repair. *J. Biol. Chem.* **286**(44), 38375-38382 (2011).

#### WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

#### SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

#### WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 12/12/2018

#### CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD  
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM

WWW.CAYMANCHEM.COM