

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



PRODUCT INFORMATION



Compound E

Item No. 15579

CAS Registry No.: 209986-17-4

Formal Name: N-[(1S)-2-[[(3S)-2,3-dihydro-1-

methyl-2-oxo-5-phenyl-1H-1,4benzodiazepin-3-yl]amino]-1methyl-2-oxoethyl]-3,5-difluoro-

benzeneacetamide


Synonym: γ-Secretase Inhibitor XXI

MF: $C_{27}H_{24}F_2N_4O_3$

FW: 490.5 **Purity:** ≥98% UV/Vis.: λ_{max} : 228 nm Supplied as: A crystalline solid

-20°C Storage: Stability: ≥2 years

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

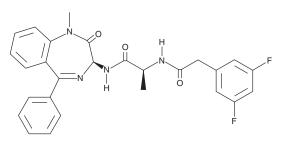

Laboratory Procedures

Compound E is supplied as a crystalline solid. A stock solution may be made by dissolving the compound E in the solvent of choice, which should be purged with an inert gas. Compound E is soluble in organic solvents such as ethanol, DMSO, and dimethyl formamide (DMF). The solubility of compound E in these solvents is approximately 2, 20, and 25 mg/ml, respectively.

Compound E is sparingly soluble in aqueous buffers. For maximum solubility in aqueous buffers, compound E should first be dissolved in DMF and then diluted with the aqueous buffer of choice. Compound E has a solubility of approximately 0.03 mg/ml in a 1:30 solution of DMF:PBS (pH 7.2) using this method. We do not recommend storing the aqueous solution for more than one day.

Description

y-Secretase is a multimeric aspartyl protease that regulates signaling pathways by proteolytically cleaving substrates, abrogating or releasing signaling molecules. Two well-known substrates are the carboxylterminal fragments (CTFs) of the receptor Notch, which has key roles in development, and that of amyloid precursor protein (APP), which is important in Alzheimer's disease. Compound E is a potent, cell-permeable, and selective inhibitor of γ -secretase, blocking the cleavage of both APP and Notch CTFs with IC₅₀ values of ~0.3 nM.²⁻⁴ Compound E induces neuronal differentiation, impairs ovarian folliculogenesis, and suppresses thymocyte development by preventing Notch activation by y-secretase. 5-7

References

- 1. Jurisch-Yaksi, N., Sannerud, R., and Annaert, W. Biochim Biophys. Acta. 1828(12), 2815-2827 (2013).
- 2. Seiffert, D., Bradley, J.D., Rominger, C.M., et al. J. Biol. Chem. 275(44), 34086-34091 (2000).
- 3. Beher, D., Wrigley, J.D.J., Nadin, A., et al. J. Biol. Chem. 276(48), 45394-45402 (2001).
- 4. Zhao, G., Mao, G., Tan, J., et al. J. Biol. Chem. 279(49), 50647-50650 (2004).
- 5. Ferrari-Toninelli, G., Bonini, S.A., Uberti, D., et al. Neuro. Oncol. 12(12), 1231-1243 (2010).
- 6. Jovanovic, V.P., Sauer, C.M., Shawber, C.J., et al. Reprod. Biol. Endocrinol. 11, 43 (2013).
- 7. Doerfler, P., Shearman, M.S., and Perlmutter, R.M. Proc. Natl. Acad. Sci. USA 98(16), 9312-9317 (2001).

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information Buyer agrees to purchase the mater can be found on our website.

Copyright Cayman Chemical Company, 02/24/2023

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640 CUSTSERV@CAYMANCHEM.COM WWW.**CAYMANCHEM**.COM