



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com


[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

PRODUCT INFORMATION



ON-123300
Item No. 19902

CAS Registry No.: 1357470-29-1
Formal Name: 8-cyclopentyl-7,8-dihydro-2-[[4-(4-methyl-1-piperazinyl)phenyl]amino]-7-oxo-pyrido[2,3-d]pyrimidine-6-carbonitrile
MF: C₂₄H₂₇N₇O
FW: 429.5
Purity: ≥98%
UV/Vis.: λ_{max}: 224, 269, 315, 392 nm
Supplied as: A crystalline solid
Storage: -20°C
Stability: ≥2 years

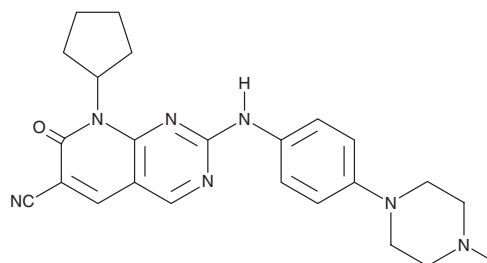

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

ON-123300 is supplied as a crystalline solid. A stock solution may be made by dissolving the ON-123300 in the solvent of choice. ON-123300 is soluble in organic solvents such as ethanol, DMSO, and dimethyl formamide (DMF), which should be purged with an inert gas. The solubility of ON-123300 in these solvents is approximately 0.25, 10, and 16 mg/ml, respectively.

ON-123300 is sparingly soluble in aqueous buffers. For maximum solubility in aqueous buffers, ON-123300 should first be dissolved in DMF and then diluted with the aqueous buffer of choice. ON-123300 has a solubility of approximately 0.3 mg/ml in a 1:2 solution of DMF:PBS (pH 7.2) using this method. We do not recommend storing the aqueous solution for more than one day.

Description

ON-123300 is a multi-kinase inhibitor.^{1,2} It inhibits ARK5, Cdk4, and Cdk6 (IC₅₀s = 3.87, 9.82, and 4.95 nM, respectively) as well as FYN, PDGFRβ, Abl, and PI3Kδ (IC₅₀s = 11.09, 26, 53.32, and 144 nM, respectively) among others.² It inhibits proliferation in a panel of 38 cancer cell lines (GI₅₀s = 0.05-5 μM).² ON-123300 also inhibits proliferation of U87 glioma cells (IC₅₀ = 3.4 μM), decreases phosphorylation of Akt and activates ERK, as well as halts the cell cycle in the G₂/M phase and induces apoptosis.¹ When used in combination with the EGFR inhibitor gefitinib (Item No. 13166), it synergistically induces cytotoxicity of U87, U87 VIII, and U87 PTEN cells. ON-123300 reduces tumor growth in an MDA-MB-231 mouse xenograft model when administered at a dose of 50 mg/kg every other day.²

References

1. Zhang, X., Lv, H., Zhou, Q., *et al.* Preclinical pharmacological evaluation of a novel multiple kinase inhibitor, ON123300, in brain tumor models. *Mol. Cancer Ther.* **13**(5), 1105-1116 (2014).
2. Reddy, M.V.R., Akula, B., Conseza, S.C., *et al.* Discovery of 8-cyclopentyl-2-[[4-(4-methyl-piperazin-1-yl)-phenylamino]-7-oxo-7,8-dihydro-pyrido[2,3-d]pyrimidine-6-carbonitrile (7x) as a potent inhibitor of cyclin-dependent kinase 4 (CDK4) and AMPK-related kinase 5 (ARK5). *J. Med. Chem.* **57**(3), 578-599 (2014).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 11/09/2018

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM