



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



MMP-2/MMP-9 Inhibitor I

Item No. 20315

CAS Registry No.: 193807-58-8
Formal Name: N-([1,1'-biphenyl]-4-ylsulfonyl)-D-phenylalanine
Synonym: Matrix Metalloproteinase-2/Matrix Metalloproteinase-9 Inhibitor I
MF: C₂₁H₁₉NO₄S
FW: 381.4
Purity: ≥99%
Supplied as: A solid
Storage: -20°C
Stability: ≥2 years

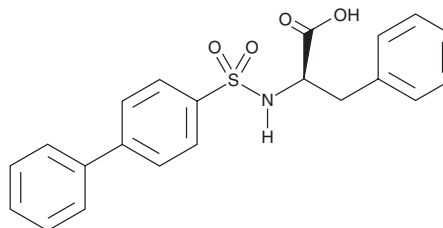

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

MMP-2/MMP-9 inhibitor I is supplied as a solid. A stock solution may be made by dissolving the MMP-2/MMP-9 inhibitor I in the solvent of choice. MMP-2/MMP-9 inhibitor I is soluble in DMSO to 25 mM.

Description

MMP-2/MMP-9 inhibitor I is a potent inhibitor of matrix metalloproteinase-2 (MMP-2) and MMP-9 (IC₅₀s = 310 and 240 nM, respectively).¹ It acts by binding zinc at the active site of these MMPs. MMP-2/MMP-9 inhibitor I is active *in vivo* as well as *in vitro* and blocks MMP-2/MMP-9-dependent invasion, tumor growth, and metastasis in both cell culture and mouse tumor models.¹ This compound has been used to elucidate the roles of MMP-2 and MMP-9 in diverse systems, including mammary epithelial cell transformation, neuronal dysfunction, lymphocyte recruitment, and progressive hereditary kidney disease.²⁻⁵

References

1. Tamura, Y., Watamane, F., Nakatani, T., *et al.* Highly selective and orally active inhibitors of type IV collagenase (MMP-9 and MMP-2): N-sulfonylamino acid derivatives. *J. Med. Chem.* **41**(4), 640-649 (1998).
2. Dechow, T.B., Pedranzini, L., Leitch, A., *et al.* Requirement of matrix metalloproteinase-9 for the transformation of human mammary epithelial cells by Stat3-C. *Proc. Natl. Acad. Sci. USA* **101**(29), 10602-10607 (2004).
3. Lim, K., Hyun, Y.-M., Lambert-Emo, K., *et al.* Neutrophil trails guide influenza-specific CD8⁺ T cells in the airways. *Science* **349**(6252), aaa4352 (2015).
4. Mizoguchi, H., Yamada, K., and Nabeshima, T. Matrix metalloproteinases contribute to neuronal dysfunction in animal models of drug dependence, Alzheimer's disease, and Epilepsy. *Biochem. Res. Int.* **2011**:681385, (2011).
5. Zeisberg, M., Khurana, M., Rao, V.H., *et al.* Stage-specific action of matrix metalloproteinases influences progressive hereditary kidney disease. *PLoS Med.* **3**(4), 535-546 (2006).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 03/07/2018

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM