



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

PRODUCT INFORMATION



Ritonavir-d₆ Item No. 25462

CAS Registry No.: 1217720-20-1
Formal Name: thiazol-5-ylmethyl ((2S,3S,5S)-3-hydroxy-5-((S)-3-methyl-2-(3-methyl-3-((2-(propan-2-yl-1,1,1,3,3,3-d₆)thiazol-4-yl)methyl)ureido)butanamido)-1,6-diphenylhexan-2-yl)carbamate

MF: C₃₇H₄₂D₆N₆O₅S₂
FW: 727.0

Chemical Purity: ≥98% (Ritonavir)


Deuterium

Incorporation: ≥99% deuterated forms (d₁-d₆); ≤1% d₀

Supplied as: A solid

Storage: -20°C

Stability: ≥4 years

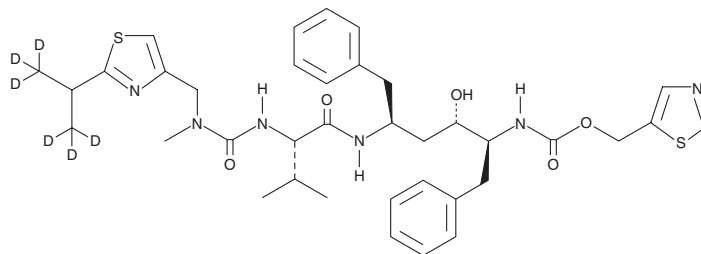

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

Ritonavir-d₆ is intended for use as an internal standard for the quantification of ritonavir (Item No. 13872) by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated *versus* unlabeled).

Ritonavir-d₆ is supplied as a solid. A stock solution may be made by dissolving the ritonavir-d₆ in the solvent of choice, which should be purged with an inert gas. Ritonavir-d₆ is slightly soluble in chloroform, DMSO, and methanol.

Description

Ritonavir is an HIV protease inhibitor.¹ It inhibits recombinant HIV-1 protease by 79% when used at a concentration of 0.5 nM. It inhibits HIV-1_{3B}-induced cell death in MT-4 human T cell leukemia cells (EC₅₀ = 25 nM) as well as cell death induced by HIV-1_{LAI}, HIV-2_{ROD}, and HIV-2_{EHO} in human MT-2 cells (IC₅₀s = 0.045, 0.13, and 0.24 μM, respectively).^{1,2} Ritonavir also inhibits the cytochrome P450 (CYP) isoform CYP3A (IC₅₀ = 0.14 μM).³ It inhibits CYP-mediated oxidative metabolism of the HIV protease inhibitors saquinavir (Item No. 9001893), indinavir (Item No. 15150), nelfinavir (Item No. 15144), and amprenavir (Item No. 15369) in rat and human liver microsomes in a concentration-dependent manner.⁴ Ritonavir (10 mg/kg) also prevents decreases in plasma levels of these four compounds in rats. Formulations containing ritonavir have been used in the treatment of HIV-1 infection.

References

1. Kempf, D.J., Shan, H.L., Marsh, K.C., *et al.* *J. Med. Chem.* **41**(4), 602-617 (1998).
2. Koh, Y., Nakata, H., Maeda, K., *et al.* *Antimicrob. Agents Chemother.* **47**(10), 3123-3129 (2003).
3. Kumar, G.N., Dykstra, J., Roberts, E.M., *et al.* *Drug Metab. Dispos.* **27**(8), 902-908 (1999).
4. Kempf, D.J., Marsh, K.C., Kumar, G., *et al.* *Antimicrob. Agents Chemother.* **41**(3), 654-660 (1997).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 12/14/2022

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM