

# Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



## Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

## Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

## SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



# PRODUCT INFORMATION



## Difloxacin-d<sub>3</sub> Item No. 27617

CAS Registry No.: 1173147-93-7

Formal Name: 6-fluoro-1-(4-fluorophenyl)-7-(4-

> (methyl-d<sub>2</sub>)piperazin-1-yl)-4-oxo-1,4dihydroquinoline-3-carboxylic acid

MF:  $C_{21}H_{16}D_3F_2N_3O_3$ 402.4 FW:


**Chemical Purity:** ≥98% (Difloxacin)

Deuterium

Incorporation:  $\geq$ 99% deuterated forms (d<sub>1</sub>-d<sub>3</sub>);  $\leq$ 1% d<sub>0</sub>

Supplied as: A solid -20°C Storage: Stability: ≥2 years

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

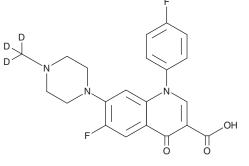

#### **Laboratory Procedures**

Difloxacin-d3 is intended for use as an internal standard for the quantification of difloxacin by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

Difloxacin-d<sub>3</sub> is supplied as a solid. A stock solution may be made by dissolving the difloxacin-d<sub>3</sub> in the solvent of choice, which should be purged with an inert gas. Difloxacin-d3 is soluble in ethanol and chloroform.

#### Description

Difloxacin is a fluoroquinolone antibiotic. 1 It is active against isolates of anaerobic bacteria, including B. fragilis, as well as Fusobacterium and Actinomyces species with MIC values ranging from ≤0.125 to 8 µg/ml. It eliminates E. coli and B. fragilis infection in a rat intra-abdominal abscess model when administered at a dose of 40 mg/kg three times per day.<sup>2</sup>

### Reference

- 1. Bansal, M.B. and Thadepalli, H. Activity of difloxacin (A-56619) and A-56620 against clinical anaerobic bacteria in vitro. Antimicrob. Agents Chemother. 31(4), 619-621 (1987).
- Thadepalli, H., Gollapudi, S.V., and Chuah, S.K. Therapeutic evaluation of difloxacin (A-56619) and A-56620 for experimentally induced Bacteroides fragilis-associated intra-abdominal abscess. Antimicrob. Agents Chemother. 30(4), 574-576 (1986).

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

## WARRANTY AND LIMITATION OF REMEDY

uyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website

Copyright Cayman Chemical Company, 04/17/2019

#### **CAYMAN CHEMICAL**

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640 CUSTSERV@CAYMANCHEM.COM WWW.**CAYMANCHEM**.COM