



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



3 α ,5 α -THDOC

Item No. 28364

CAS Registry No.: 567-02-2
Formal Name: (3 α ,5 α)-3,21-dihydroxy-pregnan-20-one
Synonyms: 3 α ,5 α -THDOC,
3 α ,5 α -Tetrahydrodeoxycorticosterone,
5 α ,3 α -THDOC, NSC 113927
MF: C₂₁H₃₄O₃
FW: 334.5
Purity: \geq 95%
Supplied as: A solid
Storage: -20°C
Stability: \geq 2 years

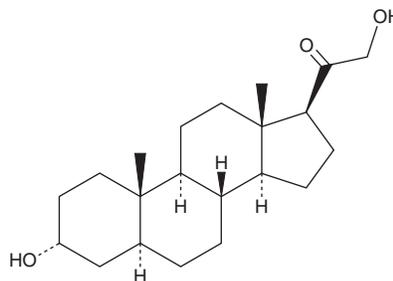

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

3 α ,5 α -THDOC is supplied as a solid. A stock solution may be made by dissolving the 3 α ,5 α -THDOC in the solvent of choice, which should be purged with an inert gas. 3 α ,5 α -THDOC is soluble in organic solvents such as ethanol and DMSO. The solubility of 3 α ,5 α -THDOC in these solvents is approximately 100 mM.

Description

3 α ,5 α -THDOC is a neurosteroid and positive allosteric modulator of GABA_A receptors.¹ It inhibits binding of the convulsant *t*-butylbicyclophosphorothionate (TBPS) and increases binding of the benzodiazepine flunitrazepam to rat synaptosomal membrane preparations in a concentration-dependent manner.² 3 α ,5 α -THDOC potentiates GABA-induced chloride currents in cultured rat hippocampal and spinal cord neurons. It inhibits seizures induced by pilocarpine or pentylenetetrazol (PTZ; Item No. 18682) in mice (ED₅₀s = 7.3 and 15 mg/kg, respectively).³ 3 α ,5 α -THDOC (20 mg/kg) increases the number of entries into and percentage of time spent in the open arms of the elevated plus maze in mice, indicating anxiolytic activity.⁴

References

1. Usami, N., Yamamoto, T., Shintani, S., *et al.* Substrate specificity of human 3(20) α -hydroxysteroid dehydrogenase for neurosteroids and its inhibition by benzodiazepines. *Biol. Pharm. Bull.* **25**(4), 441-445 (2002).
2. Majewska, M.D., Harrison, N.L., Schwartz, R.D., *et al.* Steroid hormone metabolites are barbiturate-like modulators of the GABA receptor. *Science* **232**(4753), 1004-1007 (1986).
3. Kokate, T.G., Cohen, A.L., Karp, E., *et al.* Neuroactive steroids protect against pilocarpine- and kainic acid-induced limbic seizures and status epilepticus in mice. *Neuropharmacology* **35**(8), 1049-1056 (1996).
4. Rodgers, R.J. and Johnson, N.J. Behaviorally selective effects of neuroactive steroids on plus-maze anxiety in mice. *Pharmacol. Biochem. Behav.* **59**(1), 221-232 (1998).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 11/23/2021

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM