



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)




PRODUCT INFORMATION



L-Glutamine-d₅

Item No. 34839

CAS Registry No.: 14341-78-7
Formal Name: L-Glutamine-2,3,3,4,4-d₅
Synonyms: (S)-2,5-Diamino-5-Oxopentanoic Acid-d₅,
(+)-Glutamine-d₅
MF: C₅H₅D₅N₂O₃
FW: 151.2
Chemical Purity: ≥98% (L-Glutamine)
Deuterium
Incorporation: ≥99% deuterated forms (d₁-d₅); ≤1% d₀
Supplied as: A solid
Storage: -20°C
Stability: ≥2 years

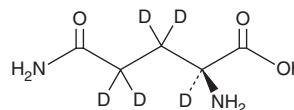

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

L-Glutamine-d₅ is intended for use as an internal standard for the quantification of L-glutamine (Item No. 23716) by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

Description

L-Glutamine is a non-essential amino acid involved in many biochemical processes. It is synthesized *in vivo* by glutamate and ammonia.² It serves as a substrate for glutamine synthetase in neurons for the biosynthesis of the major excitatory and inhibitory neurotransmitters glutamate and GABA.¹ L-Glutamine decreases adhesion of sickle red blood cells (RBCs) to human umbilical vein endothelial cells (HUVECs) when incubated *ex vivo* with patient-derived autologous plasma either alone or with LPS.³ It has commonly been used in cell culture media. Formulations containing L-glutamine have been used in the treatment of sickle cell disease.

References

1. Bak, L.K., Schousboe, A., and Waagepetersen, H.S. The glutamate/GABA-glutamine cycle: Aspects of transport, neurotransmitter homeostasis and ammonia transfer. *J. Neurochem.* **98**(3), 641-653 (2006).
2. Eisenberg, D., Gill, H.S., Pfluegl, G.M.U., *et al.* Structure-function relationships of glutamine synthetases. *Biochim Biophys. Acta.* **1477**(1-2), 122-145 (2000).
3. Niihara, Y., Matsui, N.M., Shen, Y.M., *et al.* L-glutamine therapy reduces endothelial adhesion of sickle red blood cells to human umbilical vein endothelial cells. *BMC Blood Disord.* **5**(4), (2005).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 08/23/2021

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM