



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



Ophthalmic Acid

Item No. 35343

CAS Registry No.: 495-27-2
Formal Name: L-γ-glutamyl-(2S)-2-aminobutanoyl-glycine
MF: C₁₁H₁₉N₃O₆
FW: 289.3
Purity: ≥98%
Supplied as: A solution in ethanol
Storage: -20°C
Stability: ≥2 years

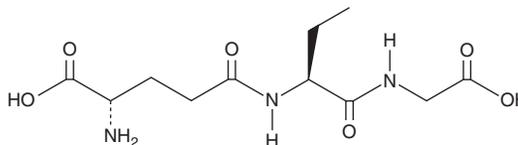

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

Ophthalmic acid is supplied as a solution in ethanol. To change the solvent, simply evaporate the ethanol under a gentle stream of nitrogen and immediately add the solvent of choice. Solvents such as ethanol, DMSO, and dimethyl formamide (DMF) purged with an inert gas can be used. The solubility of ophthalmic acid in ethanol is approximately 12 mg/ml and approximately 5 mg/ml in DMSO and DMF.

Ophthalmic acid is sparingly soluble in aqueous buffers. For maximum solubility in aqueous buffers, the ethanolic solution of ophthalmic acid should be diluted with the aqueous buffer of choice. The solubility of ophthalmic acid in PBS (pH 7.2) is approximately 3 mg/ml. We do not recommend storing the aqueous solution for more than one day.

Description

Ophthalmic acid is a glutathione (GSH) analog and biomarker of hepatic glutathione depletion.¹ It is formed via the GSH synthesis pathway, substituting 2-aminobutyrate for cysteine. Plasma and hepatic levels of ophthalmic acid are increased in fasted mice and in a mouse model of acetaminophen-induced hepatotoxicity.^{1,2}

References

1. Kobayashi, S., Lee, J., Takao, T., *et al.* Increased ophthalmic acid production is supported by amino acid catabolism under fasting conditions in mice. *Biochem. Biophys. Res. Commun.* **491(3)**, 649-655 (2017).
2. Soga, T., Baran, R., Suematsu, M., *et al.* Differential metabolomics reveals ophthalmic acid as an oxidative stress biomarker indicating hepatic glutathione consumption. *J. Biol. Chem.* **281(24)**, 16768-16776 (2006).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 11/30/2022

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM