



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

EX302 Instant Exosomes™ from LnCAP cell line (Human prostate adenocarcinoma)

Description

Instant Exosomes™ isolated from LnCAP cell line (human prostate adenocarcinoma). Exosomes are isolated following a combination of precipitation and size exclusion chromatography (SEC). The exosome samples are characterized for overall protein content, using a BCA assay, expression of commonly expressed exosomal markers (CD9, CD63 and CD81) using the ExoLISA™ exosome detection assay and, particle concentration and size distribution by nanoparticle tracking analysis (NTA).

Properties


Protein content per vial	25 µg
Method of isolation	Precipitation and size exclusion chromatography (SEC)
Characterized by	ExoLISA™ exosome detection assay, nanoparticle tracking analysis (NTA) and BCA assay
Removal of bovine Evs	Serum-free media replaces FBS-supplemented media for the final 72 hours of the cell culture process, significantly depleting the number of residual bovine EVs present in Instant Exosomes™
Nanoparticles/ml (average)	1 x 10 ⁹
State	Freeze-dried
Stability and Storage	12 months from the date of receipt when stored at -20°C as supplied. Up to 1 month when stored at -20°C, or 6 months when stored at -80°C, after reconstituting as directed. Avoid repeated freeze-thaw cycles.

Reconstitution Procedure


Reconstitute each vial with 250 µL of deionised or other ultrapure water for a final concentration of 100 µg/mL. Any further dilution once reconstituted may be done in 1X PBS or cell media to maintain osmolality. Resuspend the exosomes by pipetting the solution up and down, whilst avoiding bubbles. Vortex the reconstituted sample for 60 seconds. Briefly centrifuge the sample to ensure that the solution is collected at the bottom of the tube.

Data

TRF Counts

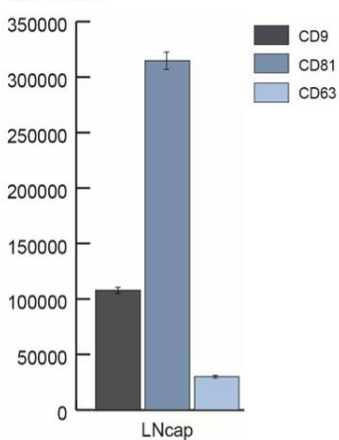

Commonly expressed exosomal markers assessed using the ExoLISA™ exosome detection assay are shown to be differentially expressed from Instant Exosomes™ isolated from LnCAP cell line.

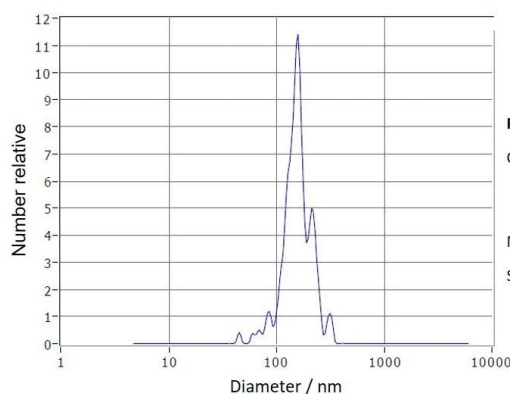

Results:

Concentration: 9.3E+9 Particles / mL

	Number	Concentration	Volume
Median (X50)	153.7	153.7	200.1
Span	48.9	48.7	59.2

Particle size distribution of Instant Exosomes™.