



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



Goat vitamin B12 (VB12) ELISA kit

Product Code	CSB-EQ028114GO																
Abbreviation	VB12																
Target Name	vitamin B12 (VB12)																
Product Type	ELISA Kit																
Immunogen Species	Capra hircus (Goat)																
Sample Types	serum, plasma, tissue homogenates																
Detection Range	1.56 pg/mL-100 pg/mL																
Sensitivity	0.39 pg/mL																
Assay Time	1-5h																
Sample Volume	50-100ul																
Detection Wavelength	450 nm																
Lead Time	3-5 working days after you place the order, and it takes another 3-5 days for delivery via DHL or FedEx.																
Research Area	Metabolism																
Tag Info	quantitative																
Protein Length	Competitive																
Description	This Goat VB12 ELISA Kit was designed for the quantitative measurement of Goat VB12 protein in serum, plasma, tissue homogenates. It is a Competitive ELISA kit, its detection range is 1.56 pg/mL-100 pg/mL and the sensitivity is 0.39 pg/mL.																
Product Precision	Intra-assay Precision (Precision within an assay): CV%<8% Three samples of known concentration were tested twenty times on one plate to assess. Inter-assay Precision (Precision between assays): CV%<10% Three samples of known concentration were tested in twenty assays to assess.																
Linearity	To assess the linearity of the assay, samples were spiked with high concentrations of goat VB12 in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay. <table> <tr> <td>?</td><td>Sample</td><td>Serum(n=4)</td></tr> <tr> <td>1:10</td><td>Average %</td><td>85</td></tr> <tr> <td></td><td>Range %</td><td>80-92</td></tr> <tr> <td>1:20</td><td>Average %</td><td>97</td></tr> <tr> <td></td><td>Range %</td><td>90-105</td></tr> </table>		?	Sample	Serum(n=4)	1:10	Average %	85		Range %	80-92	1:20	Average %	97		Range %	90-105
?	Sample	Serum(n=4)															
1:10	Average %	85															
	Range %	80-92															
1:20	Average %	97															
	Range %	90-105															



1:40	Average %	97
	Range %	92-104
1:80	Average %	92
	Range %	86-98


Recovery

The recovery of goat VB12 spiked to levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

Sample Type	Average % Recovery	Range
Serum (n=5)	95	89-100
EDTA plasma (n=4)	97	90-100

Typical

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

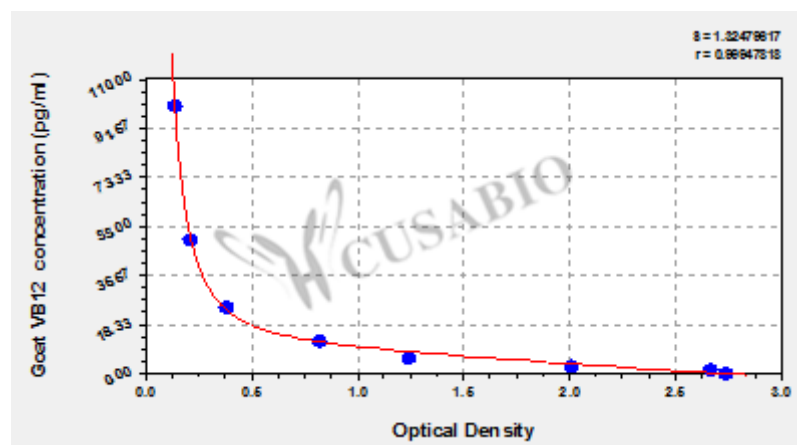

pg/ml	OD1	OD2	Average	?
100	0.145	0.148	0.147	?
50	0.221	0.227	0.224	?
25	0.396	0.394	0.395	?
12.5	0.837	0.815	0.826	?
6.25	1.252	1.243	1.248	?
3.12	2.018	2.002	2.010	?
1.56	2.679	2.641	2.660	?
0	2.754	2.716	2.735	?