



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## Prostaglandin E Receptor 4 (EP4R) ACTOne™ Stable Cell Line


CATALOG NUMBER: CL-01-EP4R

### Introduction

EP4R is a member of the G-protein coupled receptor family. This protein is one of four receptors identified for prostaglandin E2 (PGE2). This receptor can activate T-cell factor signaling. It has been shown to mediate PGE2 induced expression of early growth response 1 (EGR1), regulate the level and stability of cyclooxygenase-2 mRNA, and lead to the phosphorylation of glycogen synthase kinase-3. Knockout studies in mice suggest that this receptor may be involved in the neonatal adaptation of circulatory system, osteoporosis, as well as initiation of skin immune responses.

### Description

Human EP4R ACTOne™ is a HEK-293 CNG cell line that expresses recombinant human EP4R. HEK-293 CNG cells express a modified CNG (Cyclic Nucleotide Gated) channel that opens in response to elevated intracellular cAMP levels and consequently result in ion flux (often detectable by calcium-responsive dye, Cat# CA-C155) and cell membrane depolarization which can be easily measured with fluorescent Membrane Potential Dye (Cat# CA-M165). The assay allows both end-point and kinetic measurement of intracellular cAMP changes with a FDSS, FLIPR, or a fluorescence microplate reader.

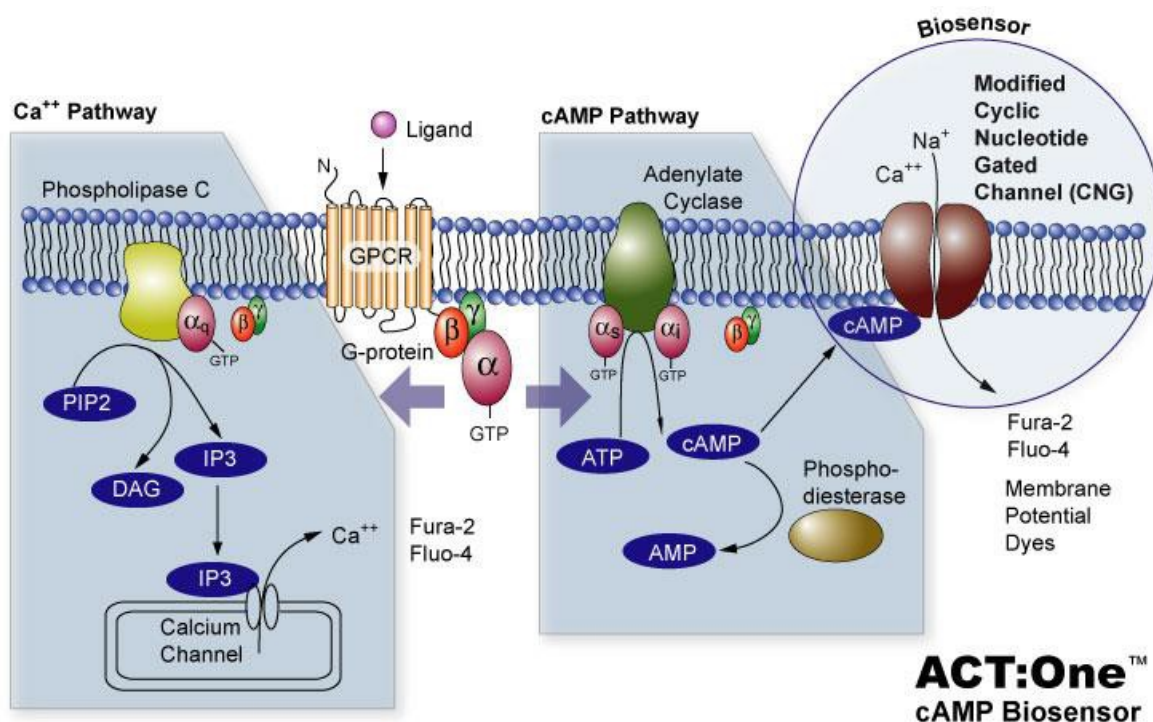

### Parental Cells

HEK-293 CNG cells (originally developed by BD Biosciences by introducing CNG in HEK-293 cells) (Cat# CL-03-PC20)

### Gene/Enzyme Introduced

EP4R (Genbank Accession No. NP\_000949)



### Applications

- cAMP dependent human EP4R receptor cell based assay
- cell based high-throughput screening of human EP4R receptor agonists/antagonists

### Functional Test

- this cell line has been tested positive for EP4R receptor specific response
- surviving rate: More than 2.5 million/vial on the second day after thawing
- the receptor specific activity is stable for 10 weeks continuous passage

### Mycoplasma Contamination Test

This lot of cells has been tested and found to be free of mycoplasma contamination.

### Content

- Stable EP4R receptor cells: 1 mL ( $1 \times 10^6$  cells/mL in 70% DMEM, 20% FBS, 10% DMSO)

### Growth Properties

Adherent

### Cell Culture Medium

- Growth medium: DMEM-10% FBS supplemented with 250  $\mu$ g/ml G418, 1  $\mu$ g/ml Puromycin
- Freezing medium: 10% DMSO, 90% complete cell culture medium

### Subculturing Procedure


1. Thaw the frozen cryovial of cells within 1-2 min by gentle agitation in a 37°C water bath. Decontaminate the cryovial by wiping the surface of the vial with 70% ethanol and transfer into a 75 cm<sup>2</sup> flask with 20 ml of complete DMEM growth medium.
2. Remove and discard culture medium next day, and then add fresh DMEM complete medium.
3. Monitor cell density daily. Cells should be passaged (1:3) when the culture reaches 90% confluence. Expected cell yield is between  $1.5 \times 10^5$  and  $2 \times 10^5$  viable cells/cm<sup>2</sup>.
4. Add 2.0 to 3.0 mL of 0.25% (w/v) trypsin-0.53 mM EDTA solution to the flask and observe cells under an inverted microscope until the cell layer is dispersed (usually within 15 to 20 minutes).  
**Note:** To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Place at 37°C to facilitate dispersal.
5. Transfer cell suspension to a 15mL centrifuge tube and spin at approximately 250 x g for 5 to 10 minutes.
6. Discard supernatant and resuspend cells in fresh growth medium. Add appropriate aliquots of the cell suspension to new culture vessels. An inoculum of 4 to 6 x 10<sup>4</sup> viable cells/cm<sup>2</sup> is recommended.
7. Incubate cultures at 37°C (5% CO<sub>2</sub>).

### Storage

Remove the frozen cells from the dry ice packaging and immediately place the cells at a temperature below -130°C, preferably in liquid nitrogen vapor, until ready for use.



## Data Example

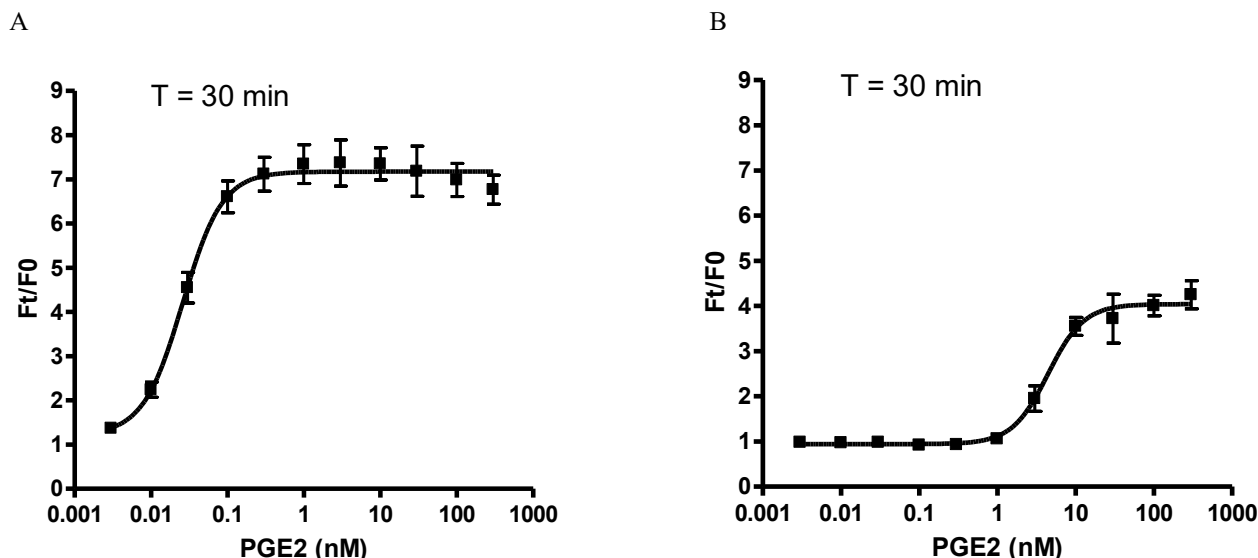

**Figure 1. Response of ACTOne™ EP4R cell line & parental cell line to PGE2.**

ACTOne™ EP4R cells and parental cells (Cat# CL-03-PC20) were plated overnight in 20  $\mu$ l culture medium on a 384 well Biocoat plate. The next day, cells were dye-loaded with 20  $\mu$ l/well of 1x Dye-loading solution (membrane potential dye kit, Cat# CA-M165). After 2 hours of incubation at room temperature, two readings were obtained prior to and 30 min after the addition of glucagon. Ratios of the two readings (F/F0) are plotted in the figure.

- A. Dose response curve of PGE2 in ACTOne™ EP4R cell line.** EC50 = 26.2 nM in the presence of PDE inhibitor Ro 20-1724, and EC50 = 121 pM in the absence of Ro20-1724 (data shown).
- B. Parental cells do not respond to PGE2 when the concentrations are less than 1 nM (in presence of PDE inhibitor). In the absence of PDE inhibitor, parental cells do not respond to PGE2 when the concentrations are less than 10 nM (data not shown).**

## Notice to Purchaser

1. This cell line is to be used for research purposes only. It may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products or to provide a service to third parties without written approval of eEnzyme LLC.
2. Refer to the license agreement for details on the usage restrictions.



Please consider the environment before printing