



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## Metabotropic Glutamate Receptor 8 (GRM8) ACTOne™ Stable Cell Line


CATALOG NUMBER: CL-11-GRM8

### Introduction

The metabotropic glutamate receptors are a family of G protein-coupled receptors, that have been divided into 3 groups on the basis of sequence homology, putative signal transduction mechanisms, and pharmacologic properties. Group I includes GRM1 and GRM5 and these receptors have been shown to activate phospholipase C. Group II includes GRM2 and GRM3 while Group III includes GRM4, GRM6, GRM7 and GRM8. Group II and III receptors are linked to the inhibition of the cyclic AMP cascade but differ in their agonist selectivities.

### Description

Human GRM8 ACTOne™ is a HEK-293 CNG-Slca3 cell line that expresses recombinant human GRM8. HEK-293 CNG-Slca3 cells express a modified CNG (Cyclic Nucleotide Gated) channel that opens in response to elevated intracellular cAMP levels and consequently result in ion flux (often detectable by calcium-responsive dye, Cat# CA-C155) and cell membrane depolarization which can be easily measured with fluorescent Membrane Potential Dye (Cat# CA-M165). The assay allows both end-point and kinetic measurement of intracellular cAMP changes with a FDSS, FLIPR, or a fluorescence microplate reader.

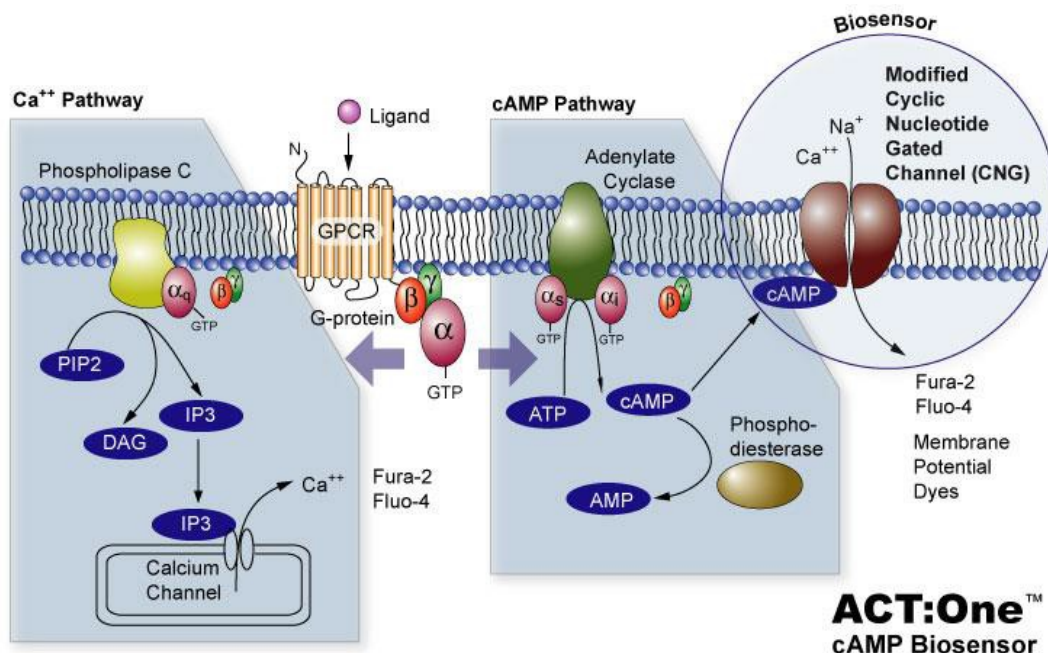

### Parental Cells

HEK-293 CNG-Slca3 cells (Cat# CL-03-PC40)

### Gene/Enzyme Introduced

GRM8 (Genbank Accession No. AAB51764.1)



Please consider the environment before printing

**Applications**

- cAMP dependent human GRM8 receptor cell based assay
- cell based high-throughput screening of human GRM8 receptor agonists/antagonists

**Functional Test**

- this cell line has been tested positive for GRM8 receptor specific response
- surviving rate: More than 2.5 million/vial on the second day after thawing
- the receptor specific activity is stable for 10 weeks continuous passage

**Mycoplasma Contamination Test**

This lot of cells has been tested and found to be free of mycoplasma contamination.

**Content**

- Stable GRM8 receptor cells: 1 mL (1 x 10<sup>6</sup> cells/mL in 70% DMEM, 20% FBS, 10% DMSO)

**Growth Properties**

Adherent

**Cell Culture Medium**

- Growth medium: DMEM-10% FBS supplemented with 250 µg/ml G418, 1 µg/ml Puromycin
- Freezing medium: 10% DMSO, 90% complete cell culture medium

**Subculturing Procedure**

1. Thaw the frozen cryovial of cells within 1-2 min by gentle agitation in a 37°C water bath. Decontaminate the cryovial by wiping the surface of the vial with 70% ethanol and transfer into a 75 cm<sup>2</sup> flask with 20 ml of complete DMEM growth medium.
2. Remove and discard culture medium next day, and then add fresh DMEM complete medium.
3. Monitor cell density daily. Cells should be passaged (1:3) when the culture reaches 90% confluence. Expected cell yield is between 1.5 x 10<sup>5</sup> and 2x 10<sup>5</sup> viable cells/cm<sup>2</sup>.
4. Add 2.0 to 3.0 mL of 0.25% (w/v) trypsin-0.53 mM EDTA solution to the flask and observe cells under an inverted microscope until the cell layer is dispersed (usually within 15 to 20 minutes).

**Note:** To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Place at 37°C to facilitate dispersal.


5. Transfer cell suspension to a 15mL centrifuge tube and spin at approximately 250 x g for 5 to 10 minutes.
6. Discard supernatant and resuspend cells in fresh growth medium. Add appropriate aliquots of the cell suspension to new culture vessels. An inoculum of 4 to 6 x 10<sup>4</sup> viable cells/cm<sup>2</sup> is recommended.
7. Incubate cultures at 37°C (5% CO<sub>2</sub>).

**Storage**

Remove the frozen cells from the dry ice packaging and immediately place the cells at a temperature below -130°C, preferably in liquid nitrogen vapor, until ready for use.



Data Example

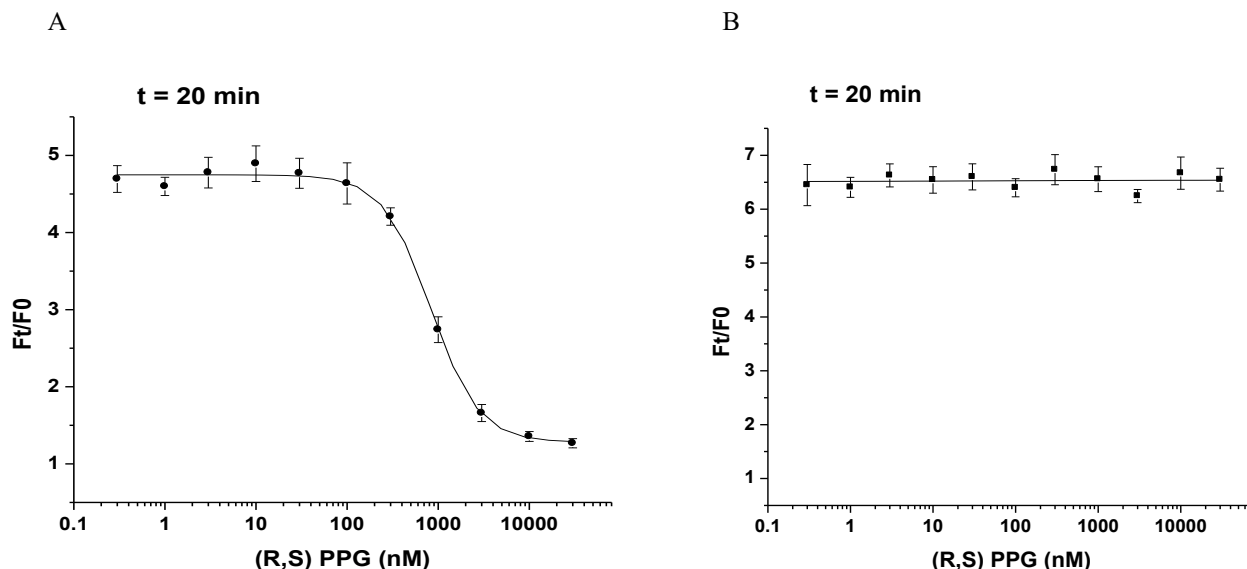

**Figure 1. Response of ACTOne™ GRM8 cell line & parental cell line to (R,S)-PHOSPHONOPHENYLGLYCINE.**

ACTOne™ GRM8 receptor cells and parental cells (Cat# CL-03-PC20) were plated overnight in 20  $\mu$ l culture medium on a 384 well Biocoat plate. The next day, cells were dye-loaded with 20  $\mu$ l/well of 1x Dye-loading solution (membrane potential dye kit, Cat# CA-M165). After 2 hours of incubation at room temperature, two readings were obtained prior to and 30 min after the addition of (R,S)-PHOSPHONOPHENYLGLYCINE. Ratios of the two readings (F/F<sub>0</sub>) are plotted in the figure.

- A. Dose response curve of (R,S)-PHOSPHONOPHENYLGLYCINE in ACTOne™ GRM8 cell line. EC<sub>50</sub> = 828 nM in the presence of PDE inhibitor Ro 20-1724 and  $\beta$ -adrenoceptor agonist isoproterenol.**
- B. Parental cells do not respond to (R,S)-PHOSPHONOPHENYLGLYCINE.**

**Notice to Purchaser**

1. This cell line is to be used for research purposes only. It may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products or to provide a service to third parties without written approval of eEnzyme LLC.
2. Refer to the license agreement for details on the usage restrictions.