



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Neuropeptide Y Receptor Y5 (NPY5R) ACTOne™ Stable Cell Line


CATALOG NUMBER: CL-11-NPY5R

Introduction

Neuropeptide Y (NPY) receptors are a family of Gi/o-protein-coupled receptors that are currently divided into four subtypes: Y1, Y2, Y4 and Y5. NPY receptors mediate a diverse range of biological actions including stimulation of food intake, anxiolysis, modulation of circadian rhythm, pain transmission and control of pituitary hormone release. All four NPY receptor subtypes are found at high concentrations in the brain and colon and are also found in a wide variety of peripheral tissues at lower levels.

Description

Human NPY5R ACTOne™ is a HEK-293 CNG cell line that expresses recombinant human NPY5R. HEK-293 CNG cells express a modified CNG (Cyclic Nucleotide Gated) channel that opens in response to elevated intracellular cAMP levels and consequently result in ion flux (often detectable by calcium-responsive dye, Cat# CA-C155) and cell membrane depolarization which can be easily measured with fluorescent Membrane Potential Dye (Cat# CA-M165). The assay allows both end-point and kinetic measurement of intracellular cAMP changes with a FDSS, FLIPR, or a fluorescence microplate reader.

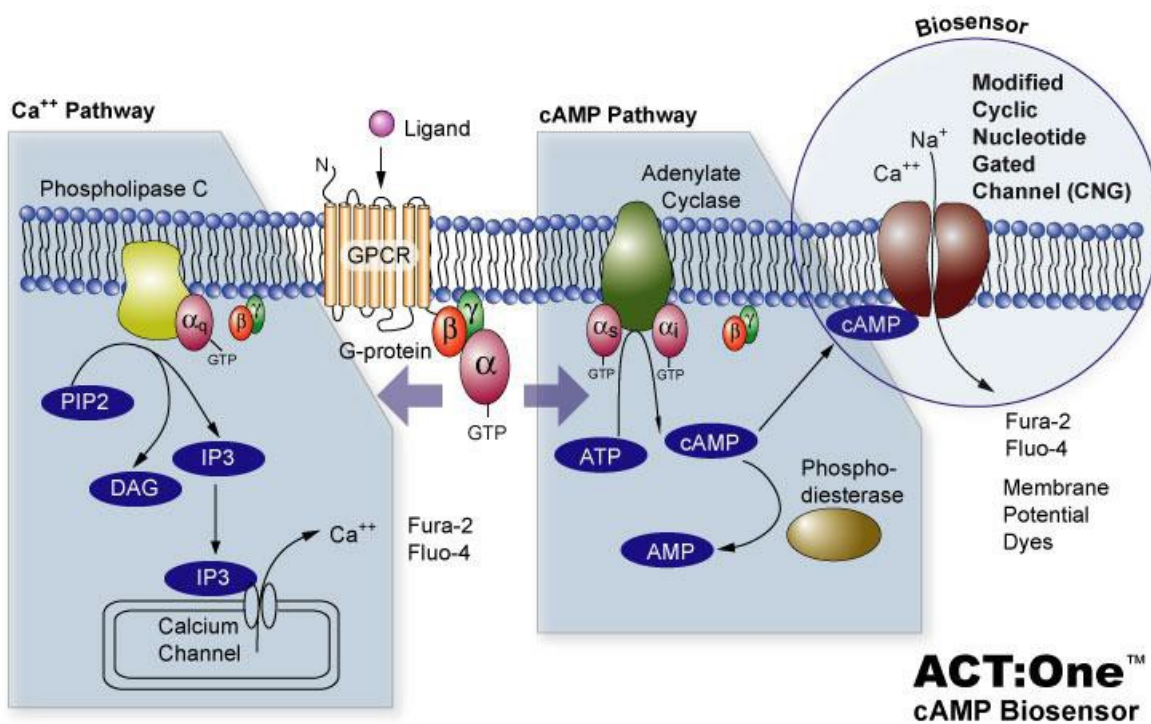

Parental Cells

HEK-293 CNG cells (originally developed by BD Biosciences by introducing CNG in HEK-293 cells) (Cat# CL-03-PC20)

Gene/Enzyme Introduced

NPY5R (Genbank Accession No. EAX-04840.1)



Applications

- cAMP dependent human NPY5R receptor cell based assay
- cell based high-throughput screening of human NPY5R receptor agonists/antagonists

Functional Test

- this cell line has been tested positive for NPY5R receptor specific response
- surviving rate: More than 2.5 million/vial on the second day after thawing
- the receptor specific activity is stable for 10 weeks continuous passage

Mycoplasma Contamination Test

This lot of cells has been tested and found to be free of mycoplasma contamination.

Content

- Stable cells: 1 mL (1×10^6 cells/mL in 70% DMEM, 20% FBS, 10% DMSO)

Growth Properties

Adherent

Cell Culture Medium

- Growth medium: DMEM-10% FBS supplemented with 250 μ g/ml G418, 1 μ g/ml Puromycin
- Freezing medium: 10% DMSO, 90% complete cell culture medium

Subculturing Procedure

1. Thaw the frozen cryovial of cells within 1-2 min by gentle agitation in a 37°C water bath. Decontaminate the cryovial by wiping the surface of the vial with 70% ethanol and transfer into a 75 cm² flask with 20 ml of complete DMEM growth medium.

Note: Pay special attention to the difference of DMEM growth medium used for stable cells and parental cells

2. Remove and discard culture medium next day, and then add fresh DMEM complete medium.
3. Monitor cell density daily. Cells should be passaged (1:3) when the culture reaches 90% confluence. Expected cell yield is between 1.5×10^5 and 2×10^5 viable cells/cm².
4. Add 2.0 to 3.0 mL of 0.25% (w/v) trypsin-0.53 mM EDTA solution to the flask and observe cells under an inverted microscope until the cell layer is dispersed (usually within 15 to 20 minutes).

Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Place at 37°C to facilitate dispersal.

5. Transfer cell suspension to a 15mL centrifuge tube and spin at approximately 250 x g for 5 to 10 minutes.
6. Discard supernatant and resuspend cells in fresh growth medium. Add appropriate aliquots of the cell suspension to new culture vessels. An inoculum of 4 to 6×10^4 viable cells/cm² is recommended.
7. Incubate cultures at 37°C (5% CO₂).

Storage

Remove the frozen cells from the dry ice packaging and immediately place the cells at a temperature below -130°C, preferably in liquid nitrogen vapor, until ready for use.



Data Example

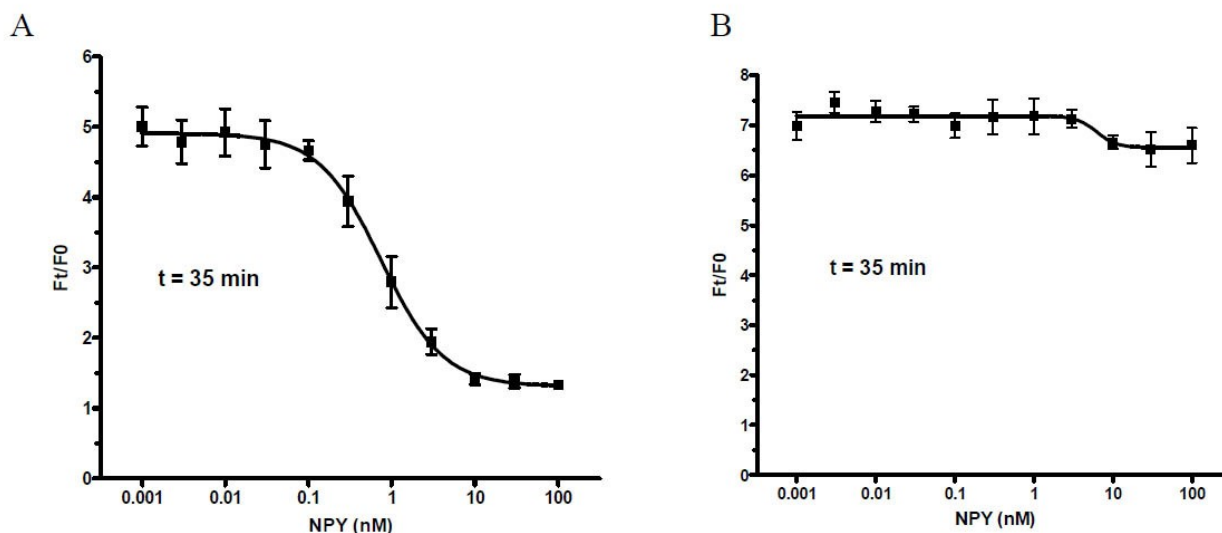


Figure 1. Response of ACTOne™ NPY5R cell line & parental cell line to NPY.

ACTOne™ NPY5R receptor cells and parental cells (Cat# CL-03-PC20) were plated overnight in 20 μ l culture medium on a 384 well BD Biocoat plate. The next day, cells were dye-loaded with 20 μ l/well of 1x Dye-loading solution (membrane potential dye kit, Cat# CA-M165). After 2 hours of incubation at room temperature, two readings were obtained prior to and 50 min after the addition of NPY. Ratios of the two readings (F/F0) are plotted in the figure.

- A. Dose response curve of NPY in ACTOne™ NPY5R cell line. EC50 = 0.74 pM in the presence of PDE inhibitor Ro 20-1724 and β -adrenoceptor agonist isoproterenol.**
- B. Parental cells do not respond to NPY.**

Notice to Purchaser

1. This cell line is to be used for research purposes only. It may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products or to provide a service to third parties without written approval of eEnzyme LLC.
2. Refer to the license agreement for details on the usage restrictions.