

# Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



# Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

# Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

## SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in







A Reliable Research Partner in Life Science and Medicine

# Elab Fluor® Red 780 Anti-Mouse CD49b/pan-NK cells Antibody[DX5]

Catalog No.E-AB-F1116USReactivityMouseStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

### **Antigen Information**

Alternate Names Integrin alpha-2,CD49 antigen-like family member B,Collagen receptor,Platelet membrane

glycoprotein Ia, GPIa, VLA-2 subunit alpha, CD49b

**Uniprot ID** Q62469 **Gene ID** 16398

**Background** DX5 antigen has been recently characterized as CD49b. It is a 150 kD integrin  $\alpha$  chain also

known as  $\alpha 2$  integrin, VLA-2  $\alpha$  chain, and integrin  $\alpha 2$  chain. CD49b non-covalently associates with CD29 ( $\beta 1$  integrin) to form the CD49b/CD29 complex known as VLA-2, a receptor for collagen and laminin. CD49b is expressed on platelets, the majority of NK cells, NKT cells, and a small subset of CD8+ T cells (this population can be significantly increased following viral infection). DX5 is used for the identification and isolation of NK cells, and is especially useful

for identifying NK cells in mice lacking the NK1.1 antigen.

#### **Product Details**

 $\begin{tabular}{lll} Form & Liquid \\ Concentration & 0.5 mg/mL \\ Size & 25 \mu g/100 \mu g \\ Clone No. & DX5 \\ Host & Rat \\ Isotype & Rat IgM, <math>\kappa$  \\ \end{tabular}

Reactivity Mouse
Application FCM

Isotype Control [Product E-AB-F09773S]

**Storage Buffer** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping** Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086
 Fax: 1-832-243-6017

 Web: <a href="https://www.elabscience.com">www.elabscience.com</a>
 Email: <a href="mailto:techsupport@elabscience.com">techsupport@elabscience.com</a>

Fax: 1-832-243-6017

A Reliable Research Partner in Life Science and Medicine

## **Fluorophore**


Conjugation: Elab Fluor® Red 780

Elab Fluor<sup>®</sup> Red 780 is designed to be excited by the Red (627-640 nm) laser and detected using an optical filter centered near 770 nm (e.g., a 780/60 nm bandpass filter).

### Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use. We suggest each investigator should titrate the reagent to obtain optimal results [The recommended concentration is 0.1-1  $\mu g/10^6$  cells in  $100~\mu L$  volume].

#### **Product data**

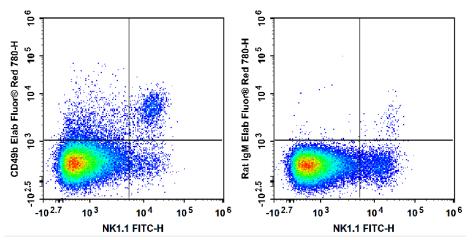

C57BL/6 murine splenocytes are stained with FITC Anti-Mouse NK1.1 Antibody and Elab Fluor<sup>®</sup> Red 780 Anti-Mouse CD49b Antibody (Left). Splenocytes are stained with FITC Anti-Mouse NK1.1 Antibody and Elab Fluor<sup>®</sup> Red 780 Rat IgM,  $\kappa$  Isotype Control (Right).

### **Related Information**

- 1. Sample Preparation for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5594.html">https://www.elabscience.com/List-detail-5594.html</a>
- 2. Staining Cell Surface Targets for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5568.html">https://www.elabscience.com/List-detail-5568.html</a>
- 3. Flow Cytometry Troubleshooting Tips <a href="https://www.elabscience.com/List-detail-5593.html">https://www.elabscience.com/List-detail-5593.html</a>
- 4. How to select the appropriate detection channel through the spectrogram? <a href="https://www.elabscience.com/List-detail-459742.html">https://www.elabscience.com/List-detail-459742.html</a>

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: <a href="www.elabscience.com">www.elabscience.com</a>
 Email: <a href="techsupport@elabscience.com">techsupport@elabscience.com</a>