

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



Elabscience Bionovation Inc.



A Reliable Research Partner in Life Science and Medicine

Elab Fluor® 647 Anti-Mouse CD150/SLAM Antibody[TC15-12F12.2]

Catalog No. E-AB-F1177M Reactivity Mouse **Storage** Store at 2~8°C, Avoid freeze / thaw cycles **Applications** FCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Slam, Slamf 1, CD 150, Signaling lymphocytic activation molecule

Uniprot ID Q9QUM4 Gene ID 27218

Background CD150 is a 75-95 kD member of the immunoglobulin superfamily, also known as SLAM

(signaling lymphocyte activation molecule) or IPO-3. CD150, a single chain type I

transmembrane molecule, is expressed on thymocytes, T cell subsets, B cells, dendritic cells, and endothelial cells. The expression is upregulated upon activation. CD150 expression has been shown to be maintained on Th1 but not Th2 clones. T regulatory cells express a relatively high level of CD150. Antibodies against CD150 have been shown to augment IFN- γ production by Th1 cells, especially when co-stimulated through the TCR. CD150 associates with the src homology 2-domain-containing protein tyrosine phosphatase SHP-2, and this association is thought to be involved in signal transduction. In combination with CD48, CD150 is a useful

marker for hematopoietic stem cell studies.

Product Details

Form Liquid

Size 50Tests/100Tests/100Tests×2

Clone No. TC15-12F12.2

Host Rat

 $\begin{array}{lll} \textbf{Isotype} & \text{Rat IgG2a, } \lambda \\ \textbf{Reactivity} & \text{Mouse} \\ \textbf{Application} & \text{FCM} \\ \end{array}$

Isotype Control Elab Fluor[®] 647 Rat IgG2a, λ Isotype Control[B39-4] [Product E-AB-F09742M]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at $2\sim8^{\circ}\text{C}$ and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Fax: 1-832-243-6017

Web: www.elabscience.com Email: techsupport@elabscience.com

Elabscience Bionovation Inc.



A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: Elab Fluor® 647

Elab Fluor $^{\otimes}$ 647 is designed to be excited by the Red laser (627-640 nm) and detected using an optical filter centered near 670 nm (e.g., a 660/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

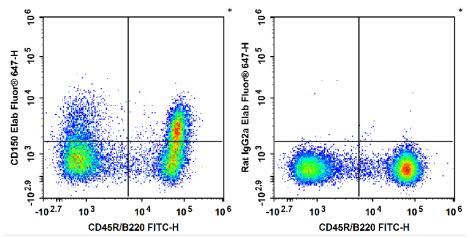

C57BL/6 murine splenocytes are stained with FITC Anti-Mouse CD45R/B220 Antibody and Elab Fluor[®] 647 Anti-Mouse CD150 Antibody (Left). Splenocytes are stained with FITC Anti-Mouse CD45R/B220 Antibody and Elab Fluor[®] 647 Rat IgG2a, λ Isotype Control (Right).

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com

Fax: 1-832-243-6017