



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Staphylococcus aureus Enterotoxin A + Enterotoxin B antibody [sa-01]

Cat No. GTX00864

Host	Mouse
Clonality	Monoclonal
Isotype	IgG2a
Application	WB, ELISA, Lateral Flow
Reactivity	Campylobacter jejuni, E. coli, Staphylococcus aureus, Salmonella Enteritidis

Package
100 µg

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Dilution
WB	1:500-1:1000
ELISA	Assay dependent
Lateral Flow	Assay dependent

Not tested in other applications.

Specificity/Sensitivity It reacts with enterotoxin A and B of Staphylococcus aureus.

PROPERTIES

Form	Liquid
Buffer	Filter-sterilized PBS, 50% glycerol
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	1 mg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Culture supernatant of Staphylococcus aureus
Purification	Protein A/G purified
Conjugation	Unconjugated
Note	For laboratory use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.



For full product information, images and publications, please visit our [website](#).


Date 2020 / 07 / 27 Page 1 of 2

DATA IMAGES

	ELISA	WB
<i>Staphylococcus aureus</i> (NBRC15306)	+	28K
Other 5 isolated strains	+	—
<i>Campylobacter jejuni</i>	—	—
<i>Salmonella enteritidis</i>	—	—
<i>Escherichia coli</i>	—	—
<i>E. coli</i> O157:H7	—	—
<i>Clostridium perfringens</i>	—	—
<i>Bacillus cereus</i>	—	—
<i>S. aureus</i> enterotoxin A (SEA)	+	28K
<i>S. aureus</i> enterotoxin B (SEB)	+	28K
<i>S. aureus</i> enterotoxin (Set A, B, C, D, E)	+	—
Set A, B, C, D, E positive control	+	—


GTX00864 Image

Reactivity of GTX00864 *Staphylococcus aureus* Enterotoxin A + Enterotoxin B antibody [sa-01] with various food poisoning bacteria.

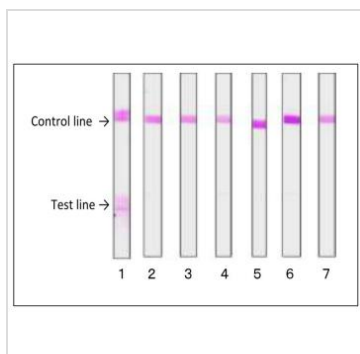

GTX00864 Lateral Flow Image

Reactivity of GTX00864 *Staphylococcus aureus* Enterotoxin A + Enterotoxin B antibody [sa-01] with various food poisoning bacteria in immunochromatographic strip test. A mouse mono antibody was coated onto a specific area (test line) of a nitrocellulose membrane, while goat anti-mouse IgG was coated onto another specific area (control line) on the same membrane. Extract of each strain of food poisoning bacteria was mixed with the MAb (sa-01) conjugated with colloidal gold. The strip was soaked and reacted with the mixture fluid. (1) *S. aureus*, (2) *Escherichia coli* O157:H7, (3) *E. coli* K12, (4) *Salmonella* Enteritidis, (5) *Campylobacter jejuni*, (6) *Vibrio parahaemolyticus*, (7) PBS. MAb (sa-01) reacted with *S. aureus*, but not with other food poisoning bacteria.

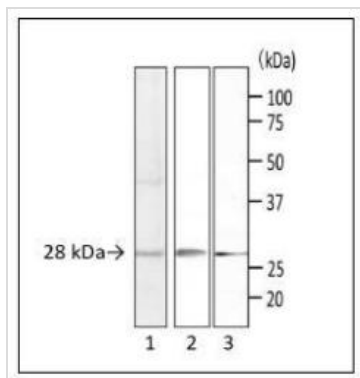


GTX00864 WB Image

WB analysis of various samples using GTX00864 *Staphylococcus aureus* Enterotoxin A + Enterotoxin B antibody [sa-01]. The molecular masses of enterotoxin A and B are 30 kDa and 29 kDa, respectively.

Lane 1 : Culture medium of *S. aureus*

Lane 2 : Purified *S. aureus* enterotoxin A

Lane 3 : Purified *S. aureus* enterotoxin B

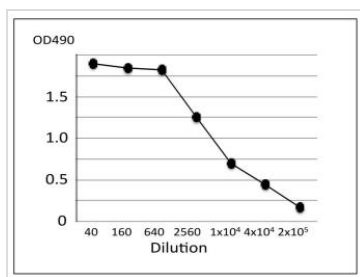

GTX00864 ELISA Image

ELISA analysis of crude extract of *S. aureus* using GTX00864 *Staphylococcus aureus* Enterotoxin A + Enterotoxin B antibody [sa-01].

Protein coating amount : 100μl (1 μg/ml) per well

Antibody amount : 100μl at indicated dilution



For full product information, images and publications, please visit our [website](http://www.genetex.com).