



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

MTAP antibody [MSVA-741R] HistoMAX™

Cat. No. GTX04964

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Application	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

Summary This antibody was validated on 76 different Normal Tissues by IHC-P.

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:50 - 1:100

Note : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Tris-EDTA-based Target Retrieval Solution buffer.

Not tested in other applications.

PROPERTIES

Form	Liquid
Buffer	10mM PBS, BSA (Please contact us for sodium azide-free or sodium azide and BSA-free format)
Preservative	sodium azide
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	22.3 µg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Recombinant human MTAP protein fragment (95-197 aa)
Purification	Protein A/G purified
Conjugation	Unconjugated

Note


For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.



For full product information, images and publications, please visit our [website](#).

DATA IMAGES

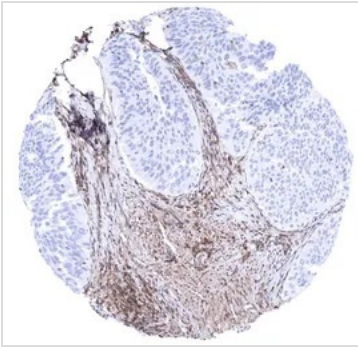

GTx04964 IHC-P Image

IHC-P analysis of human non-invasive papillary urothelial carcinoma (NIPUC) tissue using GTx04964 MTAP antibody [MSVA-741R] HistoMAX™.

Non invasive urothelial carcinoma pTa with complete absence of MTAP staining in all tumor cells. MTAP positive stroma cells serve as an internal control.



GTx04964 IHC-P Image

IHC-P analysis of human ovarian stroma tissue using GTx04964 MTAP antibody [MSVA-741R] HistoMAX™.

Moderate predominantly cytoplasmic MTAP staining of ovarian stroma cells.



GTx04964 IHC-P Image

IHC-P analysis of human urinary bladder tissue using GTx04964 MTAP antibody [MSVA-741R] HistoMAX™.

Urothelium with intense cytoplasmic and nuclear MTAP staining of all epithelial cells.



For full product information, images and publications, please visit our [website](https://www.genetex.com).