



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

EasyBlot anti-Rabbit IgG Kit (Optimized for Prot A/G)

Cat No. GTX225856-01

Application	WB
Reactivity	Rabbit

Reference (8)
 Package
 250 µl, 50 µl

PRODUCT

	EasyBlot Kit Package	Content	Package	Storage Temp
Content	GTX225856-01 (250µl)	GTX221666-01 EasyBlot anti Rabbit IgG (HRP)	250µl	-20 °C
		GTX425858 EasyBlocker	60g	4 °C
	GTX225856-01 (50µl)	GTX221666-01 EasyBlot anti Rabbit IgG (HRP)	50µl	-20 °C
		GTX425858 EasyBlocker	12g	4 °C

EasyBlot anti-Rabbit IgG is an HRP-conjugated second step reagent that specifically reacts with the native, non-reduced form of rabbit IgG and does not bind to the reduced, denatured form. Using Easyblot decreases the interference caused by the heavy (~50 kDa) and light chains (~25 kDa) of the IgG used for immunoprecipitation (IP).

Summary

The protein A/G beads used in the IP procedure sometimes bind to IgG and produce contaminating signals at approximately 50 kDa and 25 kDa. This kit also includes an EasyBlocker (GTX425858) that can minimize the background caused by Protein A/G contamination. GeneTex's EasyBlot kit is simple to use and detects only the target band in your IP experiments, resulting in cleaner western blots.

APPLICATION

Application Note

Suggested starting dilutions for WB is 1:1000. Use 2-5 µg of the IP antibody; do not exceed 5 µg. Completely reduce all IP samples by 100mM DTT (fresh prepared). For further detail, please refer to the protocol. Not yet tested in other applications. Optimal working dilutions should be determined experimentally by the end user.

PROPERTIES

Buffer	0.05M Tris (pH7.4), 0.15M NaCl, 1%BSA, 40% Glycerol
Preservative	0.005% Thimerosal
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.318 mg/ml (Please refer to the vial label for the specific concentration.)
Purification	Affinity purified by Protein G.




For full product information, images and publications, please visit our [website](#).

Date 2022 / 09 / 14 Page 1 of 2

Conjugation Horseradish peroxidase(HRP)

Note For laboratory use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

DATA IMAGES

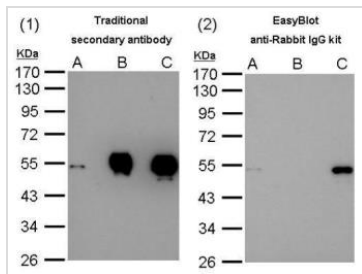

GTX225856-01 WB Image

Comparison of western blot analysis after immunoprecipitation with traditional secondary antibody (1) and EasyBlot anti-Rabbit IgG kit (GTX225856-01)(2). (A) Input, 20ug of HCT116 lysate. (B) and (C): HCT116 lysates were immunoprecipitated using 2.5 ug of normal rabbit IgG (B) or anti-p53 antibody (GTX102965) (C).

The precipitated protein was detected by GTX102965 diluted at 1:10000.



For full product information, images and publications, please visit our [website](#).